



THE UNIVERSITY *of* EDINBURGH

## Edinburgh Research Explorer

### Generation of twenty four induced pluripotent stem cell lines from twenty four members of the Lothian 4 Birth Cohort 1936

**Citation for published version:**

Toombs, J, Panther, L, Ornelas, L, Liu, C, Gomez, E, Martin-Ibanez, R, Cox, S, Ritchie, SJ, Harris, S, Taylor, A, Redmond, P, Russ, T, Murphy, L, Cooper, JD, Burr, K, Thangaraj Selvaraj, B, Browne, C, Svendsen, CN, Cowley, S, Deary, I, Chandran, S, Spires-Jones, T & Sareen, D 2020, 'Generation of twenty four induced pluripotent stem cell lines from twenty four members of the Lothian 4 Birth Cohort 1936', *Stem Cell Research*. <https://doi.org/10.1016/j.scr.2020.101851>

**Digital Object Identifier (DOI):**

[10.1016/j.scr.2020.101851](https://doi.org/10.1016/j.scr.2020.101851)

**Link:**

[Link to publication record in Edinburgh Research Explorer](#)

**Document Version:**

Peer reviewed version

**Published In:**

Stem Cell Research

**General rights**

Copyright for the publications made accessible via the Edinburgh Research Explorer is retained by the author(s) and / or other copyright owners and it is a condition of accessing these publications that users recognise and abide by the legal requirements associated with these rights.

**Take down policy**

The University of Edinburgh has made every reasonable effort to ensure that Edinburgh Research Explorer content complies with UK legislation. If you believe that the public display of this file breaches copyright please contact [openaccess@ed.ac.uk](mailto:openaccess@ed.ac.uk) providing details, and we will remove access to the work immediately and investigate your claim.



# Lab Resource: Multiple Stem Cell Lines

**Title:** Generation of twenty four induced pluripotent stem cell lines from twenty four members of the Lothian Birth Cohort 1936

**Authors:** Jamie Toombs<sup>a</sup>, Lindsay Panther<sup>b,c</sup>, Loren Ornelas<sup>b,c</sup>, Chunyan Liu<sup>b,c</sup>, Emilda Gomez<sup>b,c</sup>, Raquel Martín-Ibáñez<sup>c</sup>, Simon R. Cox<sup>d</sup>, Stuart J. Ritchie<sup>d</sup>, Sarah E. Harris<sup>d</sup>, Adele Taylor<sup>d</sup>, Paul Redmond<sup>d</sup>, Tom C. Russ<sup>d</sup>, Lee Murphy<sup>e</sup>, James D. Cooper<sup>f,k</sup>, Karen Burr<sup>f,k</sup>, Bhuvaneish T. Selvaraj<sup>f,k</sup>, Cathy Browne<sup>i</sup>, Clive N. Svendsen<sup>g,h</sup>, Sally A. Cowley<sup>i,j</sup>, Ian J. Deary<sup>d</sup>, Siddharthan Chandran<sup>f,k</sup>, Tara Spires-Jones<sup>a,†</sup>, Dhruv Sareen<sup>b,c,g,h,†</sup>

## Affiliations:

- a) Centre for Discovery Brain Sciences, UK Dementia Research Institute, The University of Edinburgh, UK.
  - b) iPSC Core, The David Janet Polak Foundation Stem Cell Core Laboratory, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
  - c) Cedars-Sinai Biomanufacturing Center, West Hollywood, CA, 90069, USA.
  - d) Lothian Birth Cohorts, Department of Psychology, University of Edinburgh, Edinburgh, UK.
  - e) Edinburgh Clinical Research Facility, University of Edinburgh, Edinburgh, UK.
  - f) Dementia Research Institute at the University of Edinburgh, UK
  - g) Board of Governors-Regenerative Medicine Institute, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
  - h) Department of Biomedical Sciences, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
  - i) James Martin Stem Cell Facility, Sir William Dunn School of Pathology, University of Oxford, Oxford, UK.
  - j) Oxford Parkinson's Disease Centre, Oxford, UK.
  - k) Centre for Clinical Brain Sciences, University of Edinburgh, Edinburgh, Scotland, UK.
- † Corresponding authors.

## Abstract:

Cognitive decline is among the most feared aspects of ageing. We have generated induced pluripotent stem cells (iPSCs) from 24 people from the Lothian Birth Cohort 1936, whose cognitive ability was tested in childhood and in older age. Peripheral blood mononuclear cells (PBMCs) were reprogrammed using non-integrating oriP/EBNA1 backbone plasmids expressing six iPSC reprogramming factors (OCT3/4 (POU5F1), SOX2, KLF4, L-Myc, shp53, Lin28, SV40LT). All lines demonstrated STR matched karyotype and pluripotency was validated by multiple methods. These iPSC lines are a valuable resource to study molecular mechanisms underlying individual differences in cognitive ageing and resilience to age-related neurodegenerative diseases.

1 **Resource Table:**

Unique stem cell lines identifier	EDi021-A EDi022-A EDi023-A EDi025-A EDi026-A EDi027-A EDi028-A EDi029-A EDi030-A EDi031-A EDi032-A EDi033-A EDi034-A EDi035-A EDi036-A EDi037-A EDi038-A EDi039-A EDi040-A EDi041-A EDi042-A EDi043-A EDi044-A EDi045-A
Alternative names of stem cell lines	N/A
Institution	Cedars-Sinai Medical Center, Los Angeles, USA
Contact information of distributor	USA distributor: Dhruv Sareen - dhruv.sareen@cshs.org  UK distributor: Karen Burr – Karen.burr@ed.ac.uk  Clinical data distributor: Paul Redmond – paul.redmond@ed.ac.uk
Type of cell lines	iPSC

Origin	Human
Cell Source	Peripheral Blood Mononuclear Cell
Clonality	Clonal
Method of reprogramming	Non-integrating episomal plasmids
Multiline rationale	24 cell lines from a shared birth year/region cohort
Gene modification	NO
Type of modification	N/A
Associated disease	N/A
Gene/locus	N/A
Method of modification	N/A
Name of transgene or resistance	N/A
Inducible/constitutive system	N/A
Date archived/stock date	EDi021-A: 14/07/2017 EDi022-A: 26/04/2017 EDi023-A: 29/03/2017 EDi025-A: 23/02/2018 EDi026-A: 30/06/2017 EDi027-A: 03/05/2017 EDi028-A: 14/06/2017 EDi029-A: 28/07/2017 EDi030-A: 19/05/2017 EDi031-A: 21/03/2018 EDi032-A: 18/01/2017 EDi033-A: 31/08/2016 EDi034-A: 16/12/2016 EDi035-A: 22/03/2017 EDi036-A: 13/01/2017 EDi037-A: 24/02/2017 EDi038-A: 23/06/2017 EDi039-A: 06/06/2018 EDi040-A: 21/06/2017 EDi041-A: 18/08/2017 EDi042-A: 14/06/2017 EDi043-A: 03/02/2017 EDi044-A: 03/05/2017



	EDi045-A: 21/02/2018
Cell line repository/bank	<p>The following lines have been added to the Cedars-Sinai iPSC Core Repository which can be viewed by the public online at <a href="https://biomanufacturing.cedars-sinai.org">https://biomanufacturing.cedars-sinai.org</a>. Direct links to each database record are included below.</p> <p>Edi021-A (<a href="#">Link</a>)</p> <p>Edi022-A (<a href="#">Link</a>)</p> <p>Edi023-A (<a href="#">Link</a>)</p> <p>Edi025-A (<a href="#">Link</a>)</p> <p>Edi026-A (<a href="#">Link</a>)</p> <p>Edi027-A (<a href="#">Link</a>)</p> <p>Edi028-A (<a href="#">Link</a>)</p> <p>Edi029-A (<a href="#">Link</a>)</p> <p>Edi030-A (<a href="#">Link</a>)</p> <p>Edi031-A (<a href="#">Link</a>)</p> <p>Edi032-A (<a href="#">Link</a>)</p> <p>Edi033-A (<a href="#">Link</a>)</p> <p>Edi034-A (<a href="#">Link</a>)</p> <p>Edi035-A (<a href="#">Link</a>)</p> <p>Edi036-A (<a href="#">Link</a>)</p> <p>Edi037-A (<a href="#">Link</a>)</p> <p>Edi038-A (<a href="#">Link</a>)</p> <p>Edi040-A (<a href="#">Link</a>)</p> <p>Edi041-A (<a href="#">Link</a>)</p> <p>Edi042-A (<a href="#">Link</a>)</p> <p>Edi043-A (<a href="#">Link</a>)</p> <p>Edi044-A (<a href="#">Link</a>)</p> <p>Edi045-A (<a href="#">Link</a>)</p>
Ethical approval	<p>NHS Lothian Research Ethics Committee: 10/S1103/10.</p> <p>CSMC Induced Pluripotent Stem Cell (iPSC) Core Facility</p> <p>Repository and Stem Cell program IRB Protocol: Pro00032834</p>

1

## 2 **Resource utility**

3 The neurobiology of cognitive ability and its decline during ageing are poorly understood. Human iPSC lines  
4 from the Lothian Birth Cohort 1936 comprise individuals with rich life-course cognitive performance data

(Taylor et al., 2018; Wardlaw et al., 2011), affording a rare model to investigate molecular mechanisms relevant to differences in brain development, cellular resilience, and vulnerability to pathology.

## Resource Details

Human peripheral blood mononuclear cells (PBMCs) were obtained from 24 unrelated members of the Lothian Birth Cohort 1936. Demographic parameters are 50% female (n = 12), 100% white Scottish (Table 1). Line donors can be grouped into 'successful', 'typical', and 'poor' cognitive ageing categories (sFig.1). Exclusion criteria were: self-reported dementia, Parkinson's disease or stroke, Mini Mental State Examination (MMSE) score <24, as well as standardised childhood IQ scores (<65, Moray House Test No. 12 at age 11), and standardised adult IQ scores (<85, average of Moray House Test No. 12 at age 70 and 76).

PBMCs were reprogrammed to generate induced pluripotent stem cells (iPSCs) using episomal plasmids encoding human OCT3/4 (POU5F1), SOX2, KLF4, L-Myc, shp53, Lin28, SV40LT. All lines were reprogrammed and stored within 22 months of each other. EBNA-related gene analysis demonstrated that iPSCs were EBNA transgene-free (and therefore exogenous reprogramming factors were no longer present) by passage 17-21 (depending on line). Qualitative tests for parental cell type by TCR- $\alpha\beta$  and TCR- $\gamma\delta$  T-cell clonality assay revealed that 83% (n = 20) of lines were non-T cell-derived, 17% (n = 4) were T-cell derived. T-cell derived lines are: EDi021-A, EDi025-A, EDi026-A, and EDi035-A. All lines have been confirmed mycoplasma negative (sFig.27).

All lines demonstrated stem cell-like morphology (Fig1F, sFig2-24F) and expressed six pluripotency markers (OCT3/4, NANOG, SOX2, TRA-1-60, TRA-1-81, SSEA4) evaluated by immunocytochemistry (Fig.1B, sFig.2-24B). Additionally, all lines demonstrated positive alkaline phosphatase AP staining (Fig.1A, sFig.2-24A) and self-renewal in undifferentiated iPSCs as assessed by PluriTest (Fig.1C, sFig.2-24C) and TaqMan<sup>®</sup>hPSC Scorecard<sup>™</sup> Panel (Fig.1D, sFig.2-24E). However, whilst EDi035-A had a positive PluriTest and Scorecard<sup>™</sup> pluripotency result, the PluriTest novelty score was borderline (1.688) (sFig.14C,E). Furthermore, EDi027-A also had a borderline positive ectoderm score as assessed by Scorecard<sup>™</sup> (sFig.6E). At 14 days of embryoid body differentiation, all lines demonstrated tri-lineage potential except EDi022-A (negative endoderm, borderline mesoderm score, sFig.2E), EDi035-A (negative mesoderm, borderline endoderm score, sFig.14E), and EDi042-A (negative endoderm score, sFig.21E), as assessed by Scorecard<sup>™</sup>.

All lines showed a normal karyotype (Fig.1D, sFig. 2-24D) between passages 6-22, with one exception. All five clones of EDi-038-A (a male) karyotyped as monosomy (45,X) (sFig.18D), and thus very likely stems from the source PBMCs. Mosaicism is a relatively common and probably harmless finding in blood cultures from normal

females and, though rarer, also in males (Bukvic et al., 2001). No differences were detected between the original PBMC samples and the corresponding iPSC lines.

All lines were confirmed to be of human origin and iPSCs matched the profile of parent PBMCs by Short Tandem Repeat (STR) analysis. Parent line data was not available for EDi026-A and EDi028-A. Genetic profiles for these lines were compared to the cell line genetic profiles available in the DSMZ STR database and did not match any other reported profiles in the DSMZ database. These profiles were found to be unique and did not match to any previously submitted profiles from the iPSC Core. The genetic profiles established here can be used for future comparisons for these cell lines. Whole genome sequence data for all 24 lines has been deposited at the European Genome-phenome Archive (EGA), which is hosted by the EBI and the CRG, under accession number EGAS00001003819.

An overview of iPSC line characterisation can be found in Table 2. Figure 1 presents example characterisation data from EDi021-A. Data for all other lines can be found in Supplementary Figures 2-27.

## Materials and Methods

### *PBMC isolation*

Blood samples were collected with NHS Lothian Research Ethics Committee Approval (10/S1103/10). Blood samples were collected in Sodium Citrate BD Vacutainer CPT tubes (BD, Cat. 362761) (three tubes per participant). For samples EDi021-A, EDi025-A, EDi028-A, EDi030-A, EDi031-A, EDi032-A, EDi033-A, EDi034-A, and EDi035-A PBMC isolation was performed by Roslin Cells. For all other lines, PBMC isolation was performed by the Edinburgh Clinical Research Facility (ECRF).

### *Generation of human iPSCs*

Generation of human iPSCs lines from PBMCs was performed using nucleofection of episomal plasmids containing POU5F1, SOX2, KLF4, LIN28, L-MYC, TP53shRNA, and SV40LT.

Briefly,  $\sim 5 \times 10^6$  cells per nucleofection of PBMCs were nucleofected with the Amaxa Human T-cell Nucleofector® Kit (Lonza, Cat. VVPA-1002) and a 5p plasmid mixture using program V-024 on a Amaxa Nucleofector 2D Device (Lonza, Cat. AAB-1001). Each transfection contained the following seven factors: OCT4, SOX2, KLF4, LMYC, LIN28, SV40LT and p53 shRNA. These were delivered on the following plasmids from Addgene, together with an EBNA1 plasmid for episomal plasmid maintenance: pEP4 E02S ET2K (Cat. 20927), pCXLE-hOCT3/4-shp53-F (Cat. 27077), pCXLEhUL (Cat. 27080), pCXLE-hSK (Cat. 27078), and pCXWB-EBNA1 (Cat. 37624). Each transfection used 0.5µg of plasmid pCXWB-EBNA1 and 0.83µg of each of the remaining four plasmids. After nucleofection, cells were immediately plated in either αβ T-cell medium (X-vivo10 [Lonza, Cat.

04-380Q] supplemented with 30U/ml IL-2 [ThermoFisher Scientific, Cat. PHC0026] and 5µl/well Dynabeads Human T-activator CD3/CD28 [Life Technologies, Cat. 11161D]) or non T-cell medium (αMEM [Life Technologies, Cat. 12561056] supplemented with 10% Heat Inactivated-FBS [Life Technologies, Cat. 10437028], 10ng/ml IL-3 [StemCell Technologies, Cat. 78040.1], 10ng/ml IL-6 [StemCell Technologies, Cat. 78050.1], 10ng/ml G-CSF [StemCell Technologies, Cat. 78012.1] and 10ng/ml GM-CSF [StemCell Technologies, Cat. 78015.1]) onto mitomycin treated mouse embryonic fibroblasts (MEF) and placed in a 37°C incubator with 20% O<sub>2</sub> and 5% CO<sub>2</sub>.

Two days after nucleofection, 2mL/well of Primate ESC medium (ReproCell, Cat. RCHEMD001) containing 5ng/ml bFGF (for MEF condition) was added to the wells without aspirating the previous medium. Beginning on day four, the medium was gently aspirated from each well and 2mL of the appropriate fresh reprogramming media was added to each well. Medium was replaced every other day. At approximately day 18 post nucleofection, individual colonies were observed in all wells of each condition. Individual PBMC-iPSC colonies with ES/iPSC-like morphology appeared between day 25-32 and those with best morphology were mechanically isolated, transferred onto 12-well plates with fresh Matrigel™ Matrix (Corning/BD Biosciences, Cat. 354230), and maintained in mTeSR®1 medium (StemCell Technologies, Cat. 85850). The iPSC clones were further expanded and scaled up for further analysis. All cultures were maintained at 37°C, 20% O<sub>2</sub>, and 5% CO<sub>2</sub> throughout the reprogramming process.

#### ***iPSC maintenance and storage***

Human iPSCs were cultured in mTeSR®1 medium (StemCell Technologies, Cat. 85850) on growth factor-reduced Matrigel™ Matrix (Corning, Cat. 354230) -coated plates at 37°C in a 20% O<sub>2</sub>, 5% CO<sub>2</sub> incubator. Briefly, 70–90% confluent human iPSC colonies were passaged every 7 days chemically (Versene, Life Technologies, Cat. 15040-066 or ReLeSR, StemCell Technologies, Cat. 05872) or mechanically by StemPro® EZPassage™ Disposable Stem Cell Passaging Tool (Life Technologies, Cat. 23181-010) and re-plated at a 1:6 or 1:9 ratio depending on the cell line. The iPSCs were passaged every 5-7 days. The iPSCs were expanded for 6-22 passages during which period various characterization assays were performed. The iPSCs were cryopreserved using CryoStor CS10 (StemCell Technologies, Cat. 07930) and an isopropanol freezing vessel at -80°C overnight. The cryopreserved vials were subsequently stored in liquid nitrogen tanks for long-term storage. Working Cell Banks (WCB) of iPSCs were cryopreserved at passage 9-14 and then Distribution Cells Banks (DCB) were created between passages 18-22.

#### ***Mycoplasma testing***

The absence of mycoplasma contamination in the iPSC lines were confirmed monthly using the MycoAlert Detection Kit, a selective biochemical test (LONZA, Cat. LT07-1188).

### ***EBNA-related gene analysis***

250ng of genomic DNA was extracted using the KingFisher™ DUO Prime purification system (Thermo Fisher Scientific) and the MagMAX™ DNA Multi-Sample Ultra 2.0 Kit (Applied Biosystems, A36570). An embryonic stem cell line (H9) was included alongside LBC lines as a negative control. DNA Amplification was conducted using TaKaRa Ex Taq® DNA Polymerase (TaKaRa Bio, RR001) and a Bio Rad 1000 Touch Thermal Cycler. Primers that recognize EBNA1 along with housekeeping gene Glyceraldehyde 3-phosphate dehydrogenase (GAPDH), which was used as a housekeeping gene, were included in this study (Table 2). PCR was run for 35 cycles at 95°C for 30 seconds, 60°C for 30 seconds, and 72°C for 30 seconds.

### ***TCRB and TCRG T-Cell Clonality Assay***

TCRB and TCRG T-Cell Clonality testing was conducted using Gene Rearrangement and Translocation assays from Invivoscribe Technologies, Inc. Genomic DNA was harvested from all iPSC lines using the MagMAX™ DNA Multi-Sample Ultra 2.0 Kit (Cat. A36570) from Applied Biosystems and it was re-suspended to a final concentration of 100-400µg per ml in dilution buffer. Three Clonal Control DNA and one Polyclonal Control DNA provided with the kit were used. PCR was carried out as per the manufacturer's protocol. PCR products were analysed using 6% TBE gel electrophoresis with gel red staining.

### ***Karyotyping***

Human PBMC-iPSCs were incubated in Colcemid (100ng/mL; Life Technologies) for 30 minutes at 37°C and then dissociated using TrypLE for 5 minutes. They were then washed in phosphate buffered saline (PBS) and incubated at 37°C in 5mL hypotonic solution (1g KCl, 1g Na Citrate in 400mL water) for 30 minutes. The cells were centrifuged for 2.5 minutes at 1500RPM and re-suspended in fixative (methanol: acetic acid, 3:1) at room temperature for 5 minutes. This was repeated twice, and finally cells were re-suspended in 500µl of fixative solution and submitted to the Cedars-Sinai Clinical Cytogenetics Core for G-band karyotyping. Karyotyping of each iPSC line was conducted at early and late passage, between passages 6-22. Approximately 20 metaphase spreads were counted per line.

### ***Immunocytochemistry***

iPSCs were plated on Matrigel™ (Corning, Cat. 354230) -coated glass coverslips or optical-bottom 96-well plates (ThermoFisher Scientific, Cat. 165305) and subsequently fixed in 4% paraformaldehyde (10 minutes, room temperature (RT)). The blocking buffer used was 5% goat serum (Millipore, Cat. S26-100ML) and 5% donkey serum (Millipore, S30-100ML) with 0.15% Triton X-100 in PBS, except for SSEA4 and OCT4 staining, for which 5% goat serum with 0.15% Triton X-100 in PBS was used as the block. All cells were blocked for one hour at RT, then incubated with primary antibodies (Table 3) for either 3 hours at RT or overnight at 4°C. Cells were then rinsed and incubated in species-specific AF488 or AF594-conjugated secondary antibodies (1:500, diluted in the same block as the primary antibodies) for one hour at RT, followed by DAPI (0.5-1µg/ml; Sigma)

to counterstain nuclei (10 minutes, RT). Cells were imaged using Nikon/Leica microscopes or Image Express. The iPSCs exhibited an embryonic stem cell like morphology, and expressed a range of pluripotency markers (OCT3/4, NANOG, SOX2, TRA-1-60, TRA-1-81, SSEA4) (Figure 1B, Supplementary Figures 2-24B).

### ***Alkaline phosphatase staining***

Alkaline phosphatase staining was performed using the Alkaline Phosphatase Staining Kit II (Stemgent, Cat. 00-0055) according to the manufacturer's instructions.

### ***PluriTest***

PluriTest was used to assess the pluripotency of undifferentiated iPSCs (Figure 1C, Supplementary Figures 2-24C). Cell pellets were sent to Life Technologies Corporation for the PluriTest Service. Total RNA was isolated using the PureLink™ RNA Mini Kit (Thermo Fisher Scientific) and quantified using NanoDrop™. 100ng total RNA was used to prepare the GeneChip® for the PluriTest™. In this assay, 36,000 transcripts and variants against a >450 sample reference set are assessed for gene expression analysis. A non-iPSC sample was used in this experiment to serve as a control for non-pluripotency. The transcriptome of all samples were analysed and processed in the PluriTest™ algorithm to generate a pluripotency and novelty score. These two scores determine the pluripotency signature of the cell line which is represented in the pluripotency plot. The threshold for pluripotency was >20, and the threshold for novelty was <1.6.

### ***hPSC Scorecard Data Analysis***

Applied Biosystems TaqMan®hPSC Scorecard™ Panel (Thermo Fisher Scientific) was used as an additional technique to assess pluripotency and tri-lineage differentiation potential of iPSC lines using real-time qPCR assays (Figure 1E, Supplementary Figures 2-24E). Total RNA from undifferentiated and EB differentiated iPSC lines was isolated using MagMAX™ mirVana™ Total RNA Isolation Kit (A27828), and 1µg of RNA was used to make cDNA using the High Capacity cDNA Reverse Transcription Kit (4368813), both from Applied Biosystems. TaqMan qRT-PCR was carried out using the hPSC Scorecard 384w Fast plate (Life technologies, A15870) and QuantStudio 12k Flex, following manufacturer protocol. We analysed the gene expression data from the TaqMan®hPSC Scorecard™ Panel using the web-based hPSC Scorecard™ Analysis Software (Thermo Fisher Scientific).

### ***Embryoid Body (EB) Formation***

iPSC lines were allowed to differentiate by EB formation. Briefly, iPSCs were lifted from 3 wells of a 6 well plate using a cell scraper and seeded in a T25 flask treated with poly-HEMA to prevent cell attachment in EB media containing: IMDM basal media (Cat. 12440061), 17% KnockOut Serum Replacement (KOSR; Cat. 10828028), 1% non-essential amino acids (Cat. 11140050), 1% Antibiotic-Antimycotic (Cat. 15240062) and 110µM β-Mercaptoethanol (Cat. 21985023), all from Thermo Fisher. EBs were allowed to form by self-

aggregation, grow and differentiate for 14 days in EB culture media replacing it twice a week. Differentiation to endoderm, mesoderm and ectoderm was assessed by TaqMan® hPSC Scorecard™ Assay (Figure 1E, Supplementary Figures 2-24E).

## **STR Analysis**

Short Tandem Repeat (STR) Analysis is conducted to confirm iPSC genetic identity. For that, a frozen vial of the parent PBMCs and a frozen vial of the reprogrammed iPSC line at late passage (18-21, depending on the cell line) are sent to IDEXX BioResearch. STR profile and interspecies contamination testing is analysed. iPSC line human authentication was conducted at IDEXX BioResearch by Cell Check™. Profiling included using a nine marker STR profile (AMEL, CSF1PO, D13S317, D16S539, D5S818, D7S820, TH01, TPOX and vWA) and interspecies contamination check for human, mouse, rat, African green monkey and Chinese hamster cells. Comparative analysis was conducted between parent PBMCs and reprogrammed iPSC lines.

## **Funding**

This work was funded by MRC Dementias Platform UK Stem Cell Network Capital Equipment MC\_EX\_MR/N50192X/1 and Partnership Award MR/N013255/1, the UK Dementia Research Institute which receives its funding from DRI Ltd, funded by the UK Medical Research Council, Alzheimer's Society, and Alzheimer's Research UK, and the European Research Council (ERC) under the European Union's Horizon 2020 research and innovation programme (Grant agreement No. 681181). Funding for the Lothian Birth Cohort 1936 (LBC1936) has been received from Research Into Ageing programme grant and the Age UK-funded Disconnected Mind project. Additional funding from the UK Medical Research Council (MRC; G0701120, G1001245, MR/M013111/1; MR/R024065/1), National Institutes of Health (R01AG054628) and the University of Edinburgh is gratefully acknowledged. This work was undertaken as part of the Cross Council and University of Edinburgh Centre for Cognitive Ageing and Cognitive Epidemiology (CCACE), funded by the Biotechnology and Biological Sciences Research Council (BBSRC) and the MRC (MR/K026992/1). A portion of the personnel support for the generation and maintenance of iPSCs was supported by Cedars-Sinai Institutional Funds. The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

## **Acknowledgements**

We thank the LBC1936 participants who took part in this study, and gratefully acknowledge the contribution of the late Professor John M. Starr, who was the Lothian Birth Cohort 1936's research medical doctor from its beginning in 2004 until 2018. Additionally, we would like to thank the Edinburgh Clinical Research Facility nursing staff, Roslin Cells, and the LBC research team members for their contributions to sample collection,

sample processing, and data processing respectively. Finally, we thank The David and Janet Polak Foundation for their support of the Cedars-Sinai iPSC Core laboratory.

**Competing interests**

US patent US 10,221,395 B2 has been granted describing some of the methods to reprogram to iPSCs. Apart from this issued patent filing the authors have declared that no other competing financial interests exist.

**References**

Bukvic, N., Gentile, M., Susca, F., Fanelli, M., Serio, G., Buonadonna, L., Capurso, A., Guanti, G., 2001. Sex chromosome loss, micronuclei, sister chromatid exchange and aging: a study including 16 centenarians. *Mutat. Res. Toxicol. Environ. Mutagen.* 498, 159–167. [https://doi.org/10.1016/S1383-5718\(01\)00279-0](https://doi.org/10.1016/S1383-5718(01)00279-0)

Taylor, A.M., Pattie, A., Deary, I.J., 2018. Cohort Profile Update: The Lothian Birth Cohorts of 1921 and 1936. *Int. J. Epidemiol.* 47, 1042-1042r. <https://doi.org/10.1093/ije/dyy022>

Wardlaw, J.M., Bastin, M.E., Valdés Hernández, M.C., Maniega, S.M., Royle, N.A., Morris, Z., Clayden, J.D., Sandeman, E.M., Eadie, E., Murray, C., Starr, J.M., Deary, I.J., 2011. Brain Aging, Cognition in Youth and Old Age and Vascular Disease in the Lothian Birth Cohort 1936: Rationale, Design and Methodology of the Imaging Protocol. *Int. J. Stroke* 6, 547–559. <https://doi.org/10.1111/j.1747-4949.2011.00683.x>



1 **Table 1: Summary of lines**

iPSC line names	Abbreviation in figures	Gender	Age at collection	Ethnicity	Genotype of locus	Disease
EDi021-A		M	78.8	White Scottish	N/A	N/A
EDi022-A		M	79.22	White Scottish	N/A	N/A
EDi023-A		F	79.1	White Scottish	N/A	N/A
EDi025-A		M	78	White Scottish	N/A	N/A
EDi026-A		M	79.45	White Scottish	N/A	N/A
EDi027-A		F	79.65	White Scottish	N/A	N/A
EDi028-A		M	79.1	White Scottish	N/A	N/A
EDi029-A		M	80.13	White Scottish	N/A	N/A
EDi030-A		F	78.98	White Scottish	N/A	N/A
EDi031-A		F	78	White Scottish	N/A	N/A
EDi032-A		F	79.29	White Scottish	N/A	N/A
EDi033-A		F	78.67	White Scottish	N/A	N/A
EDi034-A		F	78.68	White Scottish	N/A	N/A
EDi035-A		F	78.79	White Scottish	N/A	N/A
EDi036-A		F	79.22	White Scottish	N/A	N/A
EDi037-A		M	79.1	White Scottish	N/A	N/A

EDi038-A		M	79.19	White Scottish	N/A	N/A
EDi039-A		M	78	White Scottish	N/A	N/A
EDi040-A		M	79.67	White Scottish	N/A	N/A
EDi041-A		F	80.13	White Scottish	N/A	N/A
EDi042-A		F	79.42	White Scottish	N/A	N/A
EDi043-A		M	80.26	White Scottish	N/A	N/A
EDi044-A		F	79.85	White Scottish	N/A	N/A
EDi045-A		M	80.32	White Scottish	N/A	N/A

1 **Table 2: Characterization and validation**

2

Classification	Test	Result	Data
Morphology	Photography of phase contrast.	Normal. Colonies of small rounded cells with large nuclei.	Figure 1F; Supplementary Figures 2-24F.
Phenotype	Qualitative analysis: Immunofluorescence, Alkaline Phosphatase Staining.	OCT3/4+, NANOG+, SOX2+, TRA-1-60+, TRA-1-81+, SSEA4+, Alkaline Phosphatase+.	Figure 1A,B; Supplementary Figures 2-24A,B.
	Quantitative analysis: Pluritest.	Pluripotency score $\geq 20$ and novelty score $\leq 1.6$ .	Figure 1C; Supplementary Figures 2-24C.
Genotype	Karyotype (G-banding).	Normal XX and XY corresponding to gender (Table 1). Resolution 400 bands.	Figure 1D; Supplementary Figures 2-24D.
Identity	STR analysis.	9 loci tested. 100% match for lines where original PBMCs were available (22/24 lines).	Available with the authors.
		N/A	N/A
Mutation analysis (IF APPLICABLE)	Sequencing.	N/A	N/A
	Southern Blot OR WGS.	N/A	N/A
Microbiology and virology	Mycoplasma.	Negative.	sFig.27
Differentiation potential	TaqMan® hPSC Scorecard™ Assay.	Endoderm, mesoderm, ectoderm negative at day 0, positive at day 14.	Figure 1E; Supplementary Figures 2-24E.

<b>Donor screening (OPTIONAL)</b>	HIV 1 + 2 Hepatitis B, Hepatitis C.	N/A	N/A
<b>Genotype additional info (OPTIONAL)</b>	Blood group genotyping.	N/A	N/A
	HLA tissue typing.	N/A	N/A

1

2

1    **Table 3: Reagents details**

Antibodies used for immunocytochemistry/flow-cytometry			
	Antibody	Dilution	Company Cat # and RRID
Pluripotency Markers	SSEA4 (mIgG3)	1:250	Stemgent (cat. 09-0006, RRID: AB_1512169)
	TRA-1-60 (mIgM, <sub>κ</sub> )	1:250	Stemgent (cat. 09-0010, RRID: AB_1512170)
	TRA-1-81 (mIgM, <sub>κ</sub> )	1:250	Stemgent (cat. 09-0011, RRID: AB_1512171)
	OCT4 (Rabbit, IgG)	1:250	Stemgent (cat. 09-0023, RRID: AB_2167689)
	NANOG (Rabbit, IgG)	1:250	Stemgent (cat. 09-0020, RRID: AB_2298294)
	SOX2 (Rabbit, IgG)	1:250	Stemgent (cat. 09-0024, RRID: AB_2195775)
N/A	N/A	N/A	N/A
Secondary antibodies	Donkey anti-Mouse IgG AF488	1:500	Life Technologies (cat. A-21202)
	Donkey anti-Rabbit IgG AF594		Life Technologies (cat. A-21207)
	Goat anti-Mouse IgG, IgM, IgA AF488		Life Technologies (cat. A-10667)
Primers			
	Target	Forward/Reverse primer (5'-3')	
Episomal Plasmids (qPCR)	Epstein-Barr virus nuclear antigen (EBNA)	GGTCCCGAGAATCCCCATCC/ TTCATGGTCGCTGTCAGACAG	
N/A	N/A	N/A	
House-Keeping Genes (qPCR)	Glyceraldehyde 3-phosphate dehydrogenase (GAPDH)	GTGGACCTGACCTGCCGTCT/ GGAGGAGTGGGTGTCGCTGT	
N/A	N/A	N/A	
N/A	N/A	N/A	

1

2 **Figure 1: Characterization for iPSC line EDi021-A**

3

# Lab Resource: Multiple Stem Cell Lines

**Title:** Generation of twenty four induced pluripotent stem cell lines from twenty four members of the Lothian Birth Cohort 1936

**Authors:** Jamie Toombs<sup>a</sup>, Lindsay Panther<sup>b,c</sup>, Loren Ornelas<sup>b,c</sup>, Chunyan Liu<sup>b,c</sup>, Emilda Gomez<sup>b,c</sup>, Raquel Martín-Ibáñez<sup>c</sup>, Simon R. Cox<sup>d</sup>, Stuart J. Ritchie<sup>d</sup>, Sarah E. Harris<sup>d</sup>, Adele Taylor<sup>d</sup>, Paul Redmond<sup>d</sup>, Tom C. Russ<sup>d</sup>, Lee Murphy<sup>e</sup>, James D. Cooper<sup>f,k</sup>, Karen Burr<sup>f,k</sup>, Bhuvaneish T. Selvaraj<sup>f,k</sup>, Cathy Browne<sup>i</sup>, Clive N. Svendsen<sup>g,h</sup>, Sally A. Cowley<sup>i,j</sup>, Ian J. Deary<sup>d</sup>, Siddharthan Chandran<sup>f,k</sup>, Tara Spires-Jones<sup>a,†</sup>, Dhruv Sareen<sup>b,c,g,h,†</sup>

## Affiliations:

- a) Centre for Discovery Brain Sciences, UK Dementia Research Institute, The University of Edinburgh, UK.
- b) iPSC Core, The David Janet Polak Foundation Stem Cell Core Laboratory, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
- c) Cedars-Sinai Biomanufacturing Center, West Hollywood, CA, 90069, USA.
- d) Lothian Birth Cohorts, Department of Psychology, University of Edinburgh, Edinburgh, UK.
- e) Edinburgh Clinical Research Facility, University of Edinburgh, Edinburgh, UK.
- f) Dementia Research Institute at the University of Edinburgh, UK
- g) Board of Governors-Regenerative Medicine Institute, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
- h) Department of Biomedical Sciences, Cedars-Sinai Medical Center, Los Angeles, CA, 90048, USA.
- i) James Martin Stem Cell Facility, Sir William Dunn School of Pathology, University of Oxford, Oxford, UK.
- j) Oxford Parkinson's Disease Centre, Oxford, UK.
- k) Centre for Clinical Brain Sciences, University of Edinburgh, Edinburgh, Scotland, UK.

† Corresponding authors.

## Abstract:

Cognitive decline is among the most feared aspects of ageing. We have generated induced pluripotent stem cells (iPSCs) from 24 people from the Lothian Birth Cohort 1936, whose cognitive ability was tested in childhood and in older age. Peripheral blood mononuclear cells (PBMCs) were reprogrammed using non-integrating oriP/EBNA1 backbone plasmids expressing six iPSC reprogramming factors (OCT3/4 (POU5F1), SOX2, KLF4, L-Myc, shp53, Lin28, SV40LT). All lines demonstrated STR matched karyotype and pluripotency was validated by multiple methods. These iPSC lines are a valuable resource to study molecular mechanisms underlying individual differences in cognitive ageing and resilience to age-related neurodegenerative diseases.

1 **Resource Table:**

Unique stem cell lines identifier	EDi021-A EDi022-A EDi023-A EDi025-A EDi026-A EDi027-A EDi028-A EDi029-A EDi030-A EDi031-A EDi032-A EDi033-A EDi034-A EDi035-A EDi036-A EDi037-A EDi038-A EDi039-A EDi040-A EDi041-A EDi042-A EDi043-A EDi044-A EDi045-A
Alternative names of stem cell lines	N/A
Institution	Cedars-Sinai Medical Center, Los Angeles, USA
Contact information of distributor	USA distributor: Dhruv Sareen - dhruv.sareen@cshs.org  UK distributor: Karen Burr – Karen.burr@ed.ac.uk  Clinical data distributor: Paul Redmond – paul.redmond@ed.ac.uk
Type of cell lines	iPSC



Origin	Human
Cell Source	Peripheral Blood Mononuclear Cell
Clonality	Clonal
Method of reprogramming	Non-integrating episomal plasmids
Multiline rationale	24 cell lines from a shared birth year/region cohort
Gene modification	NO
Type of modification	N/A
Associated disease	N/A
Gene/locus	N/A
Method of modification	N/A
Name of transgene or resistance	N/A
Inducible/constitutive system	N/A
Date archived/stock date	EDi021-A: 14/07/2017 EDi022-A: 26/04/2017 EDi023-A: 29/03/2017 EDi025-A: 23/02/2018 EDi026-A: 30/06/2017 EDi027-A: 03/05/2017 EDi028-A: 14/06/2017 EDi029-A: 28/07/2017 EDi030-A: 19/05/2017 EDi031-A: 21/03/2018 EDi032-A: 18/01/2017 EDi033-A: 31/08/2016 EDi034-A: 16/12/2016 EDi035-A: 22/03/2017 EDi036-A: 13/01/2017 EDi037-A: 24/02/2017 EDi038-A: 23/06/2017 EDi039-A: 06/06/2018 EDi040-A: 21/06/2017 EDi041-A: 18/08/2017 EDi042-A: 14/06/2017 EDi043-A: 03/02/2017 EDi044-A: 03/05/2017

	EDi045-A: 21/02/2018
Cell line repository/bank	<p>The following lines have been added to the Cedars-Sinai iPSC Core Repository which can be viewed by the public online at <a href="https://biomanufacturing.cedars-sinai.org">https://biomanufacturing.cedars-sinai.org</a>. Direct links to each database record are included below.</p> <p>Edi021-A (<a href="#">Link</a>)</p> <p>Edi022-A (<a href="#">Link</a>)</p> <p>Edi023-A (<a href="#">Link</a>)</p> <p>Edi025-A (<a href="#">Link</a>)</p> <p>Edi026-A (<a href="#">Link</a>)</p> <p>Edi027-A (<a href="#">Link</a>)</p> <p>Edi028-A (<a href="#">Link</a>)</p> <p>Edi029-A (<a href="#">Link</a>)</p> <p>Edi030-A (<a href="#">Link</a>)</p> <p>Edi031-A (<a href="#">Link</a>)</p> <p>Edi032-A (<a href="#">Link</a>)</p> <p>Edi033-A (<a href="#">Link</a>)</p> <p>Edi034-A (<a href="#">Link</a>)</p> <p>Edi035-A (<a href="#">Link</a>)</p> <p>Edi036-A (<a href="#">Link</a>)</p> <p>Edi037-A (<a href="#">Link</a>)</p> <p>Edi038-A (<a href="#">Link</a>)</p> <p>Edi040-A (<a href="#">Link</a>)</p> <p>Edi041-A (<a href="#">Link</a>)</p> <p>Edi042-A (<a href="#">Link</a>)</p> <p>Edi043-A (<a href="#">Link</a>)</p> <p>Edi044-A (<a href="#">Link</a>)</p> <p>Edi045-A (<a href="#">Link</a>)</p>
Ethical approval	<p>NHS Lothian Research Ethics Committee: 10/S1103/10.</p> <p>CSMC Induced Pluripotent Stem Cell (iPSC) Core Facility</p> <p>Repository and Stem Cell program IRB Protocol: Pro00032834</p>

1

## 2 Resource utility

3 The neurobiology of cognitive ability and its decline during ageing are poorly understood. Human iPSC lines  
4 from the Lothian Birth Cohort 1936 comprise individuals with rich life-course cognitive performance data

(Taylor et al., 2018; Wardlaw et al., 2011), affording a rare model to investigate molecular mechanisms relevant to differences in brain development, cellular resilience, and vulnerability to pathology.

## Resource Details

Human peripheral blood mononuclear cells (PBMCs) were obtained from 24 unrelated members of the Lothian Birth Cohort 1936. Demographic parameters are 50% female (n = 12), 100% white Scottish (Table 1). Line donors can be grouped into 'successful', 'typical', and 'poor' cognitive ageing categories (sFig.1). Exclusion criteria were: self-reported dementia, Parkinson's disease or stroke, Mini Mental State Examination (MMSE) score <24, as well as standardised childhood IQ scores (<65, Moray House Test No. 12 at age 11), and standardised adult IQ scores (<85, average of Moray House Test No. 12 at age 70 and 76).

PBMCs were reprogrammed to generate induced pluripotent stem cells (iPSCs) using episomal plasmids encoding human OCT3/4 (POU5F1), SOX2, KLF4, L-Myc, shp53, Lin28, SV40LT. All lines were reprogrammed and stored within 22 months of each other. EBNA-related gene analysis demonstrated that iPSCs were EBNA transgene-free (and therefore exogenous reprogramming factors were no longer present) by passage 17-21 (depending on line). Qualitative tests for parental cell type by TCR- $\alpha\beta$  and TCR- $\gamma\delta$  T-cell clonality assay revealed that 83% (n = 20) of lines were non-T cell-derived, 17% (n = 4) were T-cell derived. T-cell derived lines are: EDi021-A, EDi025-A, EDi026-A, and EDi035-A. All lines have been confirmed mycoplasma negative (sFig.27).

All lines demonstrated stem cell-like morphology (Fig1F, sFig2-24F) and expressed six pluripotency markers (OCT3/4, NANOG, SOX2, TRA-1-60, TRA-1-81, SSEA4) evaluated by immunocytochemistry (Fig.1B, sFig.2-24B). Additionally, all lines demonstrated positive alkaline phosphatase AP staining (Fig.1A, sFig.2-24A) and self-renewal in undifferentiated iPSCs as assessed by PluriTest (Fig.1C, sFig.2-24C) and TaqMan<sup>®</sup>hPSC Scorecard<sup>™</sup> Panel (Fig.1D, sFig.2-24E). However, whilst EDi035-A had a positive PluriTest and Scorecard<sup>™</sup> pluripotency result, the PluriTest novelty score was borderline (1.688) (sFig.14C,E). Furthermore, EDi027-A also had a borderline positive ectoderm score as assessed by Scorecard<sup>™</sup> (sFig.6E). At 14 days of embryoid body differentiation, all lines demonstrated tri-lineage potential except EDi022-A (negative endoderm, borderline mesoderm score, sFig.2E), EDi035-A (negative mesoderm, borderline endoderm score, sFig.14E), and EDi042-A (negative endoderm score, sFig.21E), as assessed by Scorecard<sup>™</sup>.

All lines showed a normal karyotype (Fig.1D, sFig. 2-24D) between passages 6-22, with one exception. All five clones of EDi-038-A (a male) karyotyped as monosomy (45,X) (sFig.18D), and thus very likely stems from the source PBMCs. Mosaicism is a relatively common and probably harmless finding in blood cultures from normal

females and, though rarer, also in males (Bukvic et al., 2001). No differences were detected between the original PBMC samples and the corresponding iPSC lines.

All lines were confirmed to be of human origin and iPSCs matched the profile of parent PBMCs by Short Tandem Repeat (STR) analysis. Parent line data was not available for EDi026-A and EDi028-A. Genetic profiles for these lines were compared to the cell line genetic profiles available in the DSMZ STR database and did not match any other reported profiles in the DSMZ database. These profiles were found to be unique and did not match to any previously submitted profiles from the iPSC Core. The genetic profiles established here can be used for future comparisons for these cell lines. Whole genome sequence data for all 24 lines has been deposited at the European Genome-phenome Archive (EGA), which is hosted by the EBI and the CRG, under accession number EGAS00001003819.

An overview of iPSC line characterisation can be found in Table 2. Figure 1 presents example characterisation data from EDi021-A. Data for all other lines can be found in Supplementary Figures 2-27.

## **Materials and Methods**

### ***PBMC isolation***

Blood samples were collected with NHS Lothian Research Ethics Committee Approval (10/S1103/10). Blood samples were collected in Sodium Citrate BD Vacutainer CPT tubes (BD, Cat. 362761) (three tubes per participant). For samples EDi021-A, EDi025-A, EDi028-A, EDi030-A, EDi031-A, EDi032-A, EDi033-A, EDi034-A, and EDi035-A PBMC isolation was performed by Roslin Cells. For all other lines, PBMC isolation was performed by the Edinburgh Clinical Research Facility (ECRF).

### ***Generation of human iPSCs***

Generation of human iPSCs lines from PBMCs was performed using nucleofection of episomal plasmids containing POU5F1, SOX2, KLF4, LIN28, L-MYC, TP53shRNA, and SV40LT.

Briefly,  $\sim 5 \times 10^6$  cells per nucleofection of PBMCs were nucleofected with the Amaxa Human T-cell Nucleofector<sup>®</sup> Kit (Lonza, Cat. VVPA-1002) and a 5p plasmid mixture using program V-024 on a Amaxa Nucleofector 2D Device (Lonza, Cat. AAB-1001). Each transfection contained the following seven factors: OCT4, SOX2, KLF4, LMYC, LIN28, SV40LT and p53 shRNA. These were delivered on the following plasmids from Addgene, together with an EBNA1 plasmid for episomal plasmid maintenance: pEP4 E02S ET2K (Cat. 20927), pCXLE-hOCT3/4-shp53-F (Cat. 27077), pCXLEhUL (Cat. 27080), pCXLE-hSK (Cat. 27078), and pCXWB-EBNA1 (Cat. 37624). Each transfection used 0.5 $\mu$ g of plasmid pCXWB-EBNA1 and 0.83 $\mu$ g of each of the remaining four plasmids. After nucleofection, cells were immediately plated in either  $\alpha\beta$  T-cell medium (X-vivo10 [Lonza, Cat.

04-380Q] supplemented with 30U/ml IL-2 [ThermoFisher Scientific, Cat. PHC0026] and 5µl/well Dynabeads Human T-activator CD3/CD28 [Life Technologies, Cat. 11161D]) or non T-cell medium (αMEM [Life Technologies, Cat. 12561056] supplemented with 10% Heat Inactivated-FBS [Life Technologies, Cat. 10437028], 10ng/ml IL-3 [StemCell Technologies, Cat. 78040.1], 10ng/ml IL-6 [StemCell Technologies, Cat. 78050.1], 10ng/ml G-CSF [StemCell Technologies, Cat. 78012.1] and 10ng/ml GM-CSF [StemCell Technologies, Cat. 78015.1]) onto mitomycin treated mouse embryonic fibroblasts (MEF) and placed in a 37°C incubator with 20% O<sub>2</sub> and 5% CO<sub>2</sub>.

Two days after nucleofection, 2mL/well of Primate ESC medium (ReproCell, Cat. RCHEMD001) containing 5ng/ml bFGF (for MEF condition) was added to the wells without aspirating the previous medium. Beginning on day four, the medium was gently aspirated from each well and 2mL of the appropriate fresh reprogramming media was added to each well. Medium was replaced every other day. At approximately day 18 post nucleofection, individual colonies were observed in all wells of each condition. Individual PBMC-iPSC colonies with ES/iPSC-like morphology appeared between day 25-32 and those with best morphology were mechanically isolated, transferred onto 12-well plates with fresh Matrigel™ Matrix (Corning/BD Biosciences, Cat. 354230), and maintained in mTeSR®1 medium (StemCell Technologies, Cat. 85850). The iPSC clones were further expanded and scaled up for further analysis. All cultures were maintained at 37°C, 20% O<sub>2</sub>, and 5% CO<sub>2</sub> throughout the reprogramming process.

#### ***iPSC maintenance and storage***

Human iPSCs were cultured in mTeSR®1 medium (StemCell Technologies, Cat. 85850) on growth factor-reduced Matrigel™ Matrix (Corning, Cat. 354230) -coated plates at 37°C in a 20% O<sub>2</sub>, 5% CO<sub>2</sub> incubator. Briefly, 70–90% confluent human iPSC colonies were passaged every 7 days chemically (Versene, Life Technologies, Cat. 15040-066 or ReLeSR, StemCell Technologies, Cat. 05872) or mechanically by StemPro® EZPassage™ Disposable Stem Cell Passaging Tool (Life Technologies, Cat. 23181-010) and re-plated at a 1:6 or 1:9 ratio depending on the cell line. The iPSCs were passaged every 5-7 days. The iPSCs were expanded for 6-22 passages during which period various characterization assays were performed. The iPSCs were cryopreserved using CryoStor CS10 (StemCell Technologies, Cat. 07930) and an isopropanol freezing vessel at -80°C overnight. The cryopreserved vials were subsequently stored in liquid nitrogen tanks for long-term storage. Working Cell Banks (WCB) of iPSCs were cryopreserved at passage 9-14 and then Distribution Cells Banks (DCB) were created between passages 18-22.

#### ***Mycoplasma testing***

The absence of mycoplasma contamination in the iPSC lines were confirmed monthly using the MycoAlert Detection Kit, a selective biochemical test (LONZA, Cat. LT07-1188).

1    ***EBNA-related gene analysis***

2    250ng of genomic DNA was extracted using the KingFisher™ DUO Prime purification system (Thermo Fisher  
3    Scientific) and the MagMAX™ DNA Multi-Sample Ultra 2.0 Kit (Applied Biosystems, A36570). An embryonic  
4    stem cell line (H9) was included alongside LBC lines a negative control. DNA Amplification was conducted  
5    using TaKaRa Ex Taq® DNA Polymerase (TaKaRa Bio, RR001) and a Bio Rad 1000 Touch Thermal Cycler.  
6    Primers that recognize EBNA1 along with housekeeping gene Glyceraldehyde 3-phosphate dehydrogenase  
7    (GAPDH), which was used as a housekeeping gene, were included in this study (Table 2). PCR was run for 35  
8    cycles at 95°C for 30 seconds, 60°C for 30 seconds, and 72°C for 30 seconds.

9

10    ***TCRB and TCRG T-Cell Clonality Assay***

11    TCRB and TCRG T-Cell Clonality testing was conducted using Gene Rearrangement and Translocation assays  
12    from Invivoscribe Technologies, Inc. Genomic DNA was harvested from all iPSC lines using the MagMAX™ DNA  
13    Multi-Sample Ultra 2.0 Kit (Cat. A36570) from Applied Biosystems and it was re-suspended to a final  
14    concentration of 100-400µg per ml in dilution buffer. Three Clonal Control DNA and one Polyclonal Control  
15    DNA provided with the kit were used. PCR was carried out as per the manufacturer's protocol. PCR products  
16    were analysed using 6% TBE gel electrophoresis with gel red staining.

17

18    ***Karyotyping***

19    Human PBMC-iPSCs were incubated in Colcemid (100ng/mL; Life Technologies) for 30 minutes at 37°C and  
20    then dissociated using TrypLE for 5 minutes. They were then washed in phosphate buffered saline (PBS) and  
21    incubated at 37°C in 5mL hypotonic solution (1g KCl, 1g Na Citrate in 400mL water) for 30 minutes. The cells  
22    were centrifuged for 2.5 minutes at 1500RPM and re-suspended in fixative (methanol: acetic acid, 3:1) at  
23    room temperature for 5 minutes. This was repeated twice, and finally cells were re-suspended in 500µl of  
24    fixative solution and submitted to the Cedars-Sinai Clinical Cytogenetics Core for G-band karyotyping.  
25    Karyotyping of each iPSC line was conducted at early and late passage, between passages 6-22. Approximately  
26    20 metaphase spreads were counted per line.

27

28    ***Immunocytochemistry***

29    iPSCs were plated on Matrigel™ (Corning, Cat. 354230) -coated glass coverslips or optical-bottom 96-well  
30    plates (ThermoFisher Scientific, Cat. 165305) and subsequently fixed in 4% paraformaldehyde (10 minutes,  
31    room temperature (RT)). The blocking buffer used was 5% goat serum (Millipore, Cat. S26-100ML) and 5%  
32    donkey serum (Millipore, S30-100ML) with 0.15% Triton X-100 in PBS, except for SSEA4 and OCT4 staining, for  
33    which 5% goat serum with 0.15% Triton X-100 in PBS was used as the block. All cells were blocked for one  
34    hour at RT, then incubated with primary antibodies (Table 3) for either 3 hours at RT or overnight at 4°C. Cells  
35    were then rinsed and incubated in species-specific AF488 or AF594-conjugated secondary antibodies (1:500,  
36    diluted in the same block as the primary antibodies) for one hour at RT, followed by DAPI (0.5-1µg/ml; Sigma)

to counterstain nuclei (10 minutes, RT). Cells were imaged using Nikon/Leica microscopes or Image Express. The iPSCs exhibited an embryonic stem cell like morphology, and expressed a range of pluripotency markers (OCT3/4, NANOG, SOX2, TRA-1-60, TRA-1-81, SSEA4) (Figure 1B, Supplementary Figures 2-24B).

### ***Alkaline phosphatase staining***

Alkaline phosphatase staining was performed using the Alkaline Phosphatase Staining Kit II (Stemgent, Cat. 00-0055) according to the manufacturer's instructions.

### ***PluriTest***

PluriTest was used to assess the pluripotency of undifferentiated iPSCs (Figure 1C, Supplementary Figures 2-24C). Cell pellets were sent to Life Technologies Corporation for the PluriTest Service. Total RNA was isolated using the PureLink™ RNA Mini Kit (Thermo Fisher Scientific) and quantified using NanoDrop™. 100ng total RNA was used to prepare the GeneChip® for the PluriTest™. In this assay, 36,000 transcripts and variants against a >450 sample reference set are assessed for gene expression analysis. A non-iPSC sample was used in this experiment to serve as a control for non-pluripotency. The transcriptome of all samples were analysed and processed in the PluriTest™ algorithm to generate a pluripotency and novelty score. These two scores determine the pluripotency signature of the cell line which is represented in the pluripotency plot. The threshold for pluripotency was >20, and the threshold for novelty was <1.6.

### ***hPSC Scorecard Data Analysis***

Applied Biosystems TaqMan®hPSC Scorecard™ Panel (Thermo Fisher Scientific) was used as an additional technique to assess pluripotency and tri-lineage differentiation potential of iPSC lines using real-time qPCR assays (Figure 1E, Supplementary Figures 2-24E). Total RNA from undifferentiated and EB differentiated iPSC lines was isolated using MagMAX™ mirVana™ Total RNA Isolation Kit (A27828), and 1µg of RNA was used to make cDNA using the High Capacity cDNA Reverse Transcription Kit (4368813), both from Applied Biosystems. TaqMan qRT-PCR was carried out using the hPSC Scorecard 384w Fast plate (Life technologies, A15870) and QuantStudio 12k Flex, following manufacturer protocol. We analysed the gene expression data from the TaqMan®hPSC Scorecard™ Panel using the web-based hPSC Scorecard™ Analysis Software (Thermo Fisher Scientific).

### ***Embryoid Body (EB) Formation***

iPSC lines were allowed to differentiate by EB formation. Briefly, iPSCs were lifted from 3 wells of a 6 well plate using a cell scraper and seeded in a T25 flask treated with poly-HEMA to prevent cell attachment in EB media containing: IMDM basal media (Cat. 12440061), 17% KnockOut Serum Replacement (KOSR; Cat. 10828028), 1% non-essential amino acids (Cat. 11140050), 1% Antibiotic-Antimycotic (Cat. 15240062) and 110µM β-Mercaptoethanol (Cat. 21985023), all from Thermo Fisher. EBs were allowed to form by self-

aggregation, grow and differentiate for 14 days in EB culture media replacing it twice a week. Differentiation to endoderm, mesoderm and ectoderm was assessed by TaqMan® hPSC Scorecard™ Assay (Figure 1E, Supplementary Figures 2-24E).

## **STR Analysis**

Short Tandem Repeat (STR) Analysis is conducted to confirm iPSC genetic identity. For that, a frozen vial of the parent PBMCs and a frozen vial of the reprogrammed iPSC line at late passage (18-21, depending on the cell line) are sent to IDEXX BioResearch. STR profile and interspecies contamination testing is analysed. iPSC line human authentication was conducted at IDEXX BioResearch by Cell Check™. Profiling included using a nine marker STR profile (AMEL, CSF1PO, D13S317, D16S539, D5S818, D7S820, TH01, TPOX and vWA) and interspecies contamination check for human, mouse, rat, African green monkey and Chinese hamster cells. Comparative analysis was conducted between parent PBMCs and reprogrammed iPSC lines.

## **Funding**

This work was funded by MRC Dementias Platform UK Stem Cell Network Capital Equipment MC\_EX\_MR/N50192X/1 and Partnership Award MR/N013255/1, the UK Dementia Research Institute which receives its funding from DRI Ltd, funded by the UK Medical Research Council, Alzheimer's Society, and Alzheimer's Research UK, and the European Research Council (ERC) under the European Union's Horizon 2020 research and innovation programme (Grant agreement No. 681181). Funding for the Lothian Birth Cohort 1936 (LBC1936) has been received from Research Into Ageing programme grant and the Age UK-funded Disconnected Mind project. Additional funding from the UK Medical Research Council (MRC; G0701120, G1001245, MR/M013111/1; MR/R024065/1), National Institutes of Health (R01AG054628) and the University of Edinburgh is gratefully acknowledged. This work was undertaken as part of the Cross Council and University of Edinburgh Centre for Cognitive Ageing and Cognitive Epidemiology (CCACE), funded by the Biotechnology and Biological Sciences Research Council (BBSRC) and the MRC (MR/K026992/1). A portion of the personnel support for the generation and maintenance of iPSCs was supported by Cedars-Sinai Institutional Funds. The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

## **Acknowledgements**

We thank the LBC1936 participants who took part in this study, and gratefully acknowledge the contribution of the late Professor John M. Starr, who was the Lothian Birth Cohort 1936's research medical doctor from its beginning in 2004 until 2018. Additionally, we would like to thank the Edinburgh Clinical Research Facility nursing staff, Roslin Cells, and the LBC research team members for their contributions to sample collection,



sample processing, and data processing respectively. Finally, we thank The David and Janet Polak Foundation for their support of the Cedars-Sinai iPSC Core laboratory.

#### **Competing interests**

US patent US 10,221,395 B2 has been granted describing some of the methods to reprogram to iPSCs. Apart from this issued patent filing the authors have declared that no other competing financial interests exist.

#### **References**

- Bukvic, N., Gentile, M., Susca, F., Fanelli, M., Serio, G., Buonadonna, L., Capurso, A., Guanti, G., 2001. Sex chromosome loss, micronuclei, sister chromatid exchange and aging: a study including 16 centenarians. *Mutat. Res. Toxicol. Environ. Mutagen.* 498, 159–167. [https://doi.org/10.1016/S1383-5718\(01\)00279-0](https://doi.org/10.1016/S1383-5718(01)00279-0)
- Taylor, A.M., Pattie, A., Deary, I.J., 2018. Cohort Profile Update: The Lothian Birth Cohorts of 1921 and 1936. *Int. J. Epidemiol.* 47, 1042-1042r. <https://doi.org/10.1093/ije/dyy022>
- Wardlaw, J.M., Bastin, M.E., Valdés Hernández, M.C., Maniega, S.M., Royle, N.A., Morris, Z., Clayden, J.D., Sandeman, E.M., Eadie, E., Murray, C., Starr, J.M., Deary, I.J., 2011. Brain Aging, Cognition in Youth and Old Age and Vascular Disease in the Lothian Birth Cohort 1936: Rationale, Design and Methodology of the Imaging Protocol. *Int. J. Stroke* 6, 547–559. <https://doi.org/10.1111/j.1747-4949.2011.00683.x>

1 **Table 1: Summary of lines**

iPSC line names	Abbreviation in figures	Gender	Age at collection	Ethnicity	Genotype of locus	Disease
EDi021-A		M	78.8	White Scottish	N/A	N/A
EDi022-A		M	79.22	White Scottish	N/A	N/A
EDi023-A		F	79.1	White Scottish	N/A	N/A
EDi025-A		M	78	White Scottish	N/A	N/A
EDi026-A		M	79.45	White Scottish	N/A	N/A
EDi027-A		F	79.65	White Scottish	N/A	N/A
EDi028-A		M	79.1	White Scottish	N/A	N/A
EDi029-A		M	80.13	White Scottish	N/A	N/A
EDi030-A		F	78.98	White Scottish	N/A	N/A
EDi031-A		F	78	White Scottish	N/A	N/A
EDi032-A		F	79.29	White Scottish	N/A	N/A
EDi033-A		F	78.67	White Scottish	N/A	N/A
EDi034-A		F	78.68	White Scottish	N/A	N/A
EDi035-A		F	78.79	White Scottish	N/A	N/A
EDi036-A		F	79.22	White Scottish	N/A	N/A
EDi037-A		M	79.1	White Scottish	N/A	N/A

EDi038-A		M	79.19	White Scottish	N/A	N/A
EDi039-A		M	78	White Scottish	N/A	N/A
EDi040-A		M	79.67	White Scottish	N/A	N/A
EDi041-A		F	80.13	White Scottish	N/A	N/A
EDi042-A		F	79.42	White Scottish	N/A	N/A
EDi043-A		M	80.26	White Scottish	N/A	N/A
EDi044-A		F	79.85	White Scottish	N/A	N/A
EDi045-A		M	80.32	White Scottish	N/A	N/A

1 **Table 2: Characterization and validation**

2

Classification	Test	Result	Data
Morphology	Photography of phase contrast.	Normal. Colonies of small rounded cells with large nuclei.	Figure 1F; Supplementary Figures 2-24F.
Phenotype	Qualitative analysis: Immunofluorescence, Alkaline Phosphatase Staining.	OCT3/4+, NANOG+, SOX2+, TRA-1-60+, TRA-1-81+, SSEA4+, Alkaline Phosphatase+.	Figure 1A,B; Supplementary Figures 2-24A,B.
	Quantitative analysis: Pluritest.	Pluripotency score $\geq 20$ and novelty score $\leq 1.6$ .	Figure 1C; Supplementary Figures 2-24C.
Genotype	Karyotype (G-banding).	Normal XX and XY corresponding to gender (Table 1). Resolution 400 bands.	Figure 1D; Supplementary Figures 2-24D.
Identity	STR analysis.	9 loci tested. 100% match for lines where original PBMCs were available (22/24 lines).	Available with the authors.
		N/A	N/A
Mutation analysis (IF APPLICABLE)	Sequencing.	N/A	N/A
	Southern Blot OR WGS.	N/A	N/A
Microbiology and virology	Mycoplasma.	Negative.	sFig.27
Differentiation potential	TaqMan® hPSC Scorecard™ Assay.	Endoderm, mesoderm, ectoderm negative at day 0, positive at day 14.	Figure 1E; Supplementary Figures 2-24E.

<b>Donor screening (OPTIONAL)</b>	HIV 1 + 2 Hepatitis B, Hepatitis C.	N/A	N/A
<b>Genotype additional info (OPTIONAL)</b>	Blood group genotyping.	N/A	N/A
	HLA tissue typing.	N/A	N/A

1

2

1 **Table 3: Reagents details**


<b>Antibodies used for immunocytochemistry/flow-cytometry</b>			
	<b>Antibody</b>	<b>Dilution</b>	<b>Company Cat # and RRID</b>
Pluripotency Markers	SSEA4 (mIgG3)	1:250	Stemgent (cat. 09-0006, RRID: AB_1512169)
	TRA-1-60 (mIgM <sub>κ</sub> )	1:250	Stemgent (cat. 09-0010, RRID: AB_1512170)
	TRA-1-81 (mIgM <sub>κ</sub> )	1:250	Stemgent (cat. 09-0011, RRID: AB_1512171)
	OCT4 (Rabbit, IgG)	1:250	Stemgent (cat. 09-0023, RRID: AB_2167689)
	NANOG (Rabbit, IgG)	1:250	Stemgent (cat. 09-0020, RRID: AB_2298294)
	SOX2 (Rabbit, IgG)	1:250	Stemgent (cat. 09-0024, RRID: AB_2195775)
N/A	N/A	N/A	N/A
Secondary antibodies	Donkey anti-Mouse IgG AF488 Donkey anti-Rabbit IgG AF594 Goat anti-Mouse IgG, IgM, IgA AF488	1:500	Life Technologies (cat. A-21202) Life Technologies (cat. A-21207) Life Technologies (cat. A-10667)
<b>Primers</b>			
	<b>Target</b>	<b>Forward/Reverse primer (5'-3')</b>	
Episomal Plasmids (qPCR)	Epstein-Barr virus nuclear antigen (EBNA)	GGTCCCGAGAATCCCCATCC/ TTCATGGTCGCTGTCAGACAG	
N/A	N/A	N/A	
House-Keeping Genes (qPCR)	Glyceraldehyde 3-phosphate dehydrogenase (GAPDH)	GTGGACCTGACCTGCCGTCT/ GGAGGAGTGGGTGTCGCTGT	
N/A	N/A	N/A	
N/A	N/A	N/A	

1




2 **Figure 1: Characterization for iPSC line EDi021-A**

3


A. AP

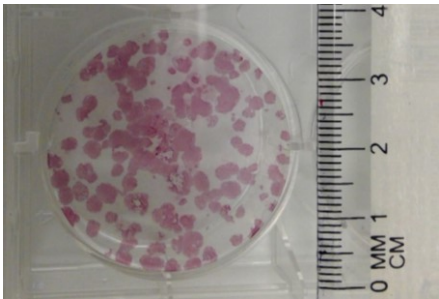


B. Immunocytochemistry

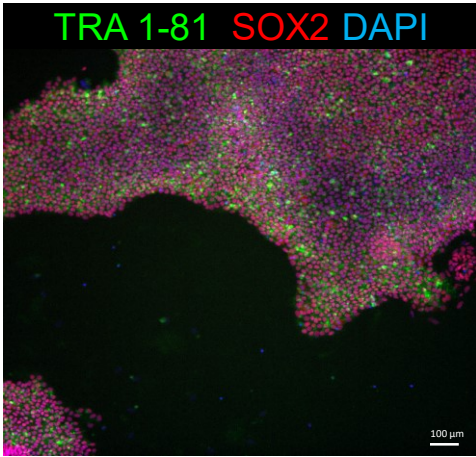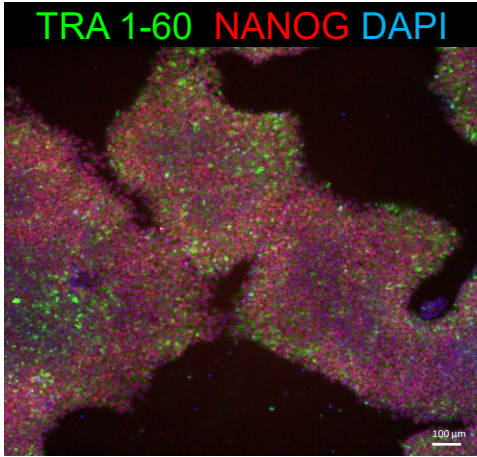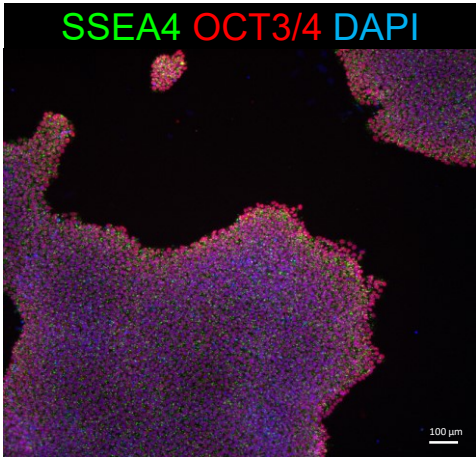

C. Pluritest

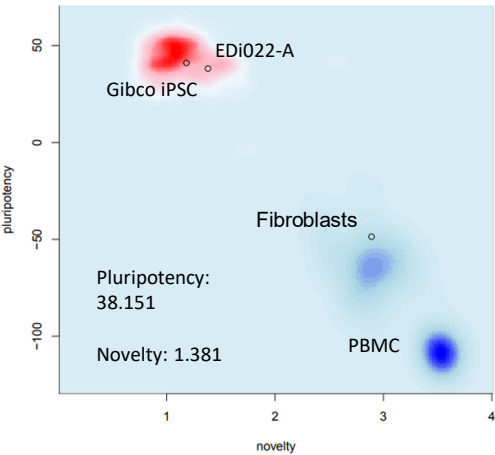


D. G-Band karyotype

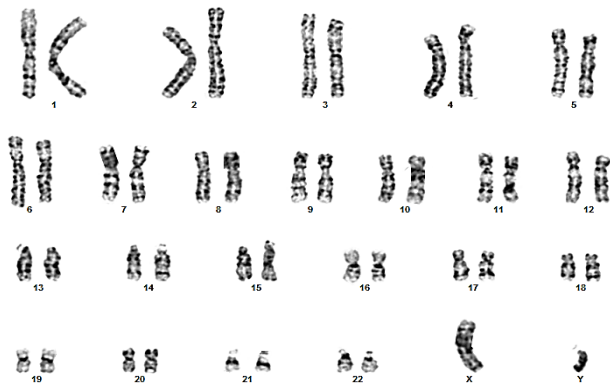

E. hPSC Scorecard


iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.66	-0.14	-0.05	-0.15	-4.84	1.99	1.83	0.65

F. Morphology  
7 days post-thaw

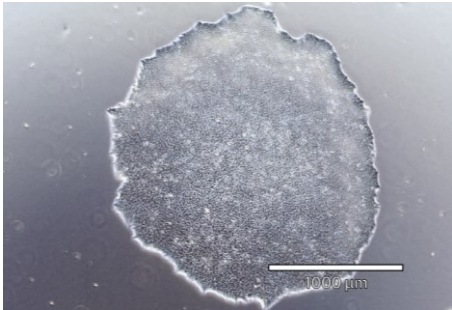




## Selection of iPSC line donors based on cognitive ageing

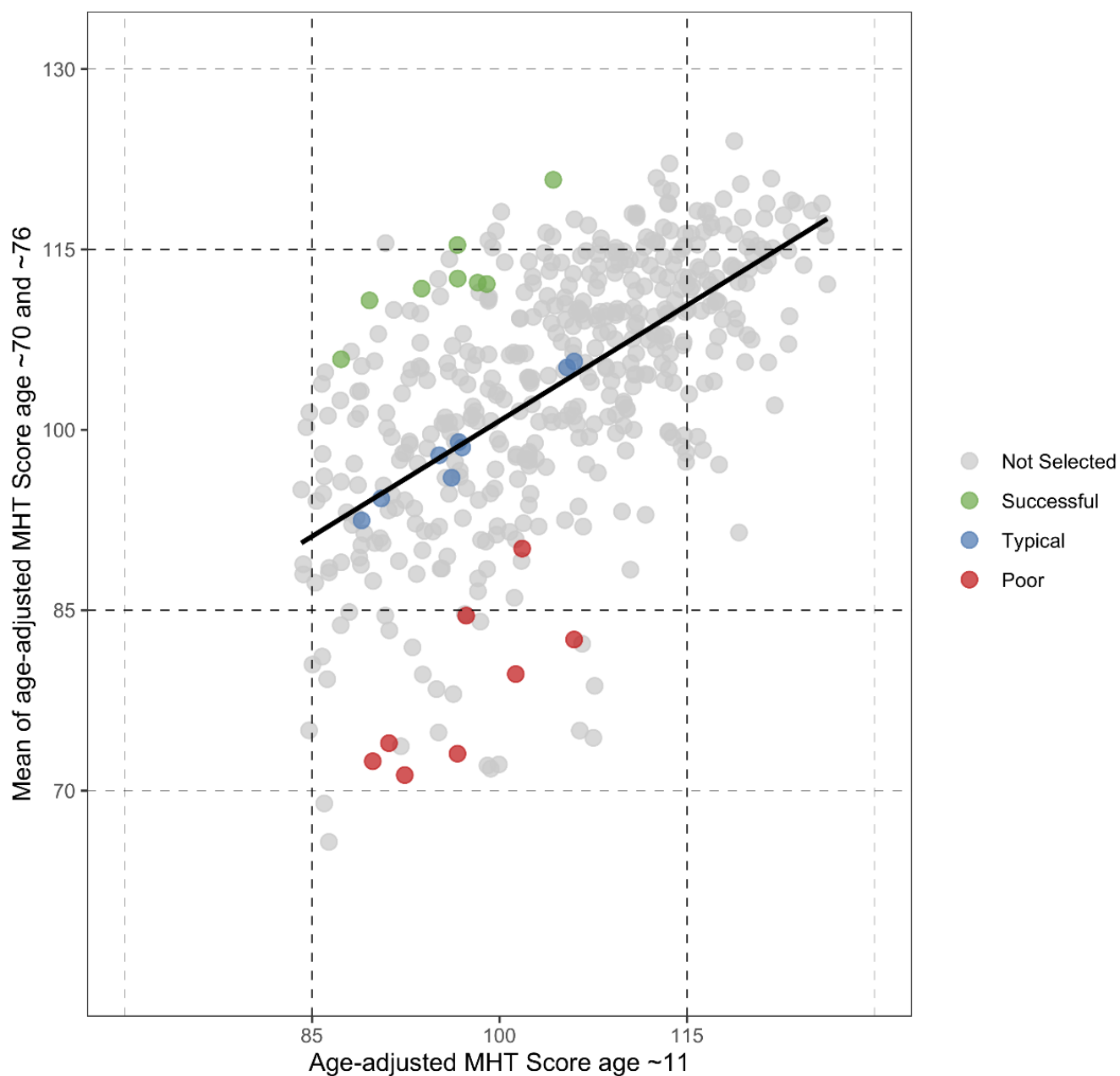




**sFig1:** Selection of iPSC candidates based on cognitive ageing profiles in the Lothian Birth Cohort 1936. Dashed lines are added to show  $\pm 1$ SD from the mean for the age 11 Moray House Test (MHT) score and  $\pm 1$  and 2SDs from the mean for later-life MHT score.

## sFig.2: Characterization for iPSC line EDi022-A


### A. AP

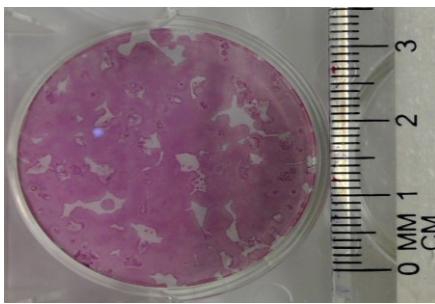

### B. Immunocytochemistry

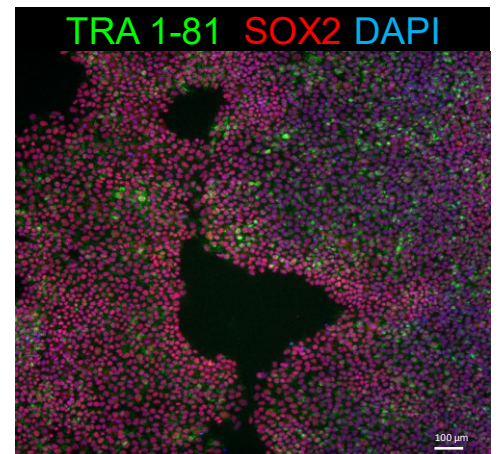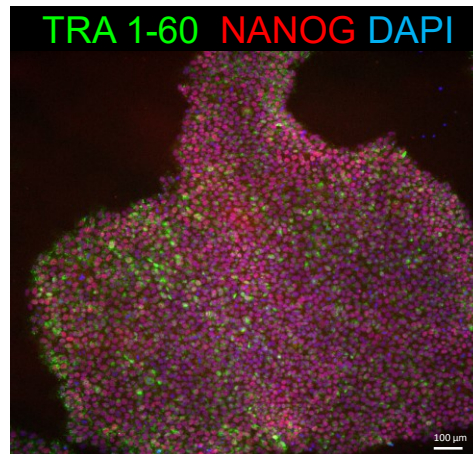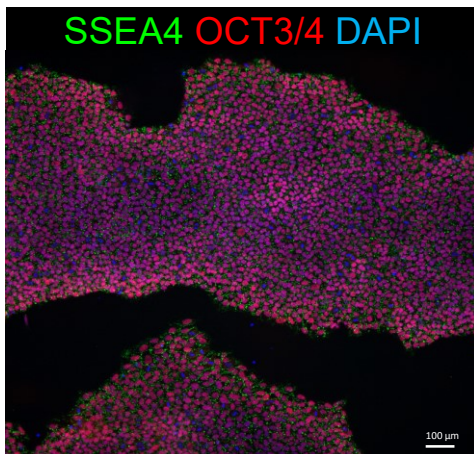

### C. PluriTest

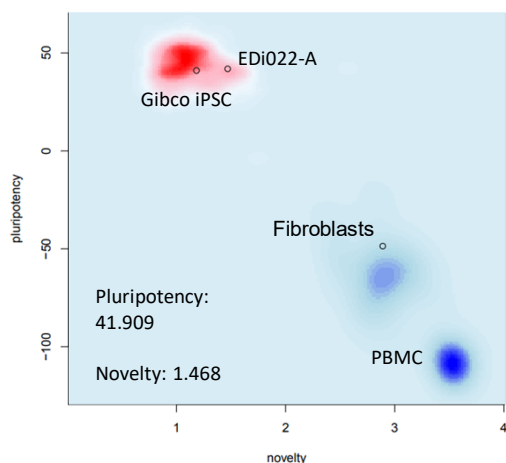

### D. G-Band karyotype




### E. hPSC Scorecard

iPSCs				Embryoid bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	-
0.29	-0.10	0.62	-0.53	-3.06	1.41	0.76	-0.32


### F. Morphology

3 days post-thaw



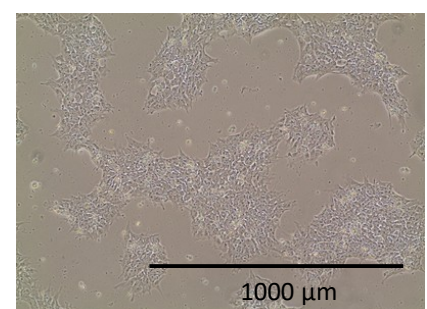


# sFig.3: Characterization for iPSC line EDi023-A


## A. AP

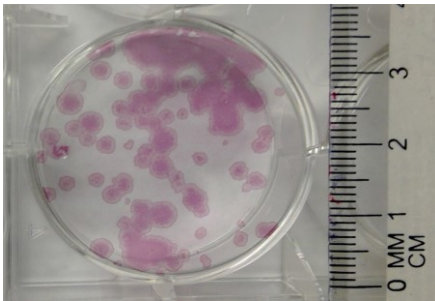

## B. Immunocytochemistry

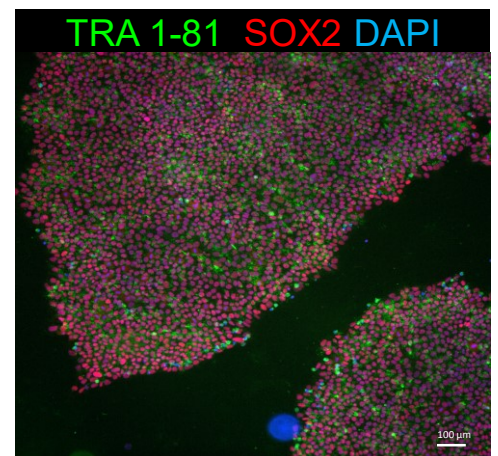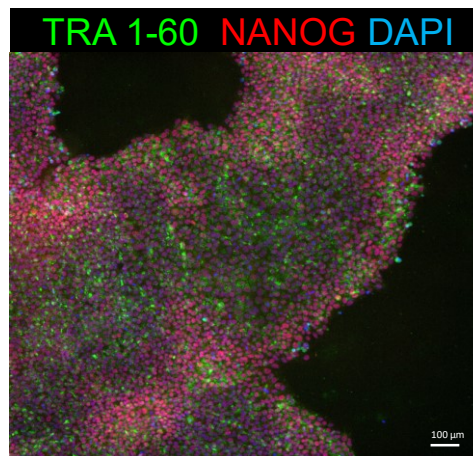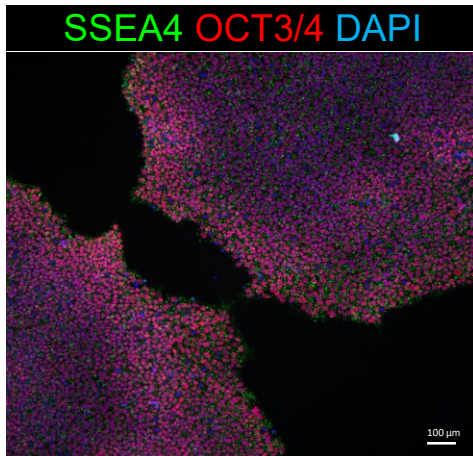

## C. PluriTest

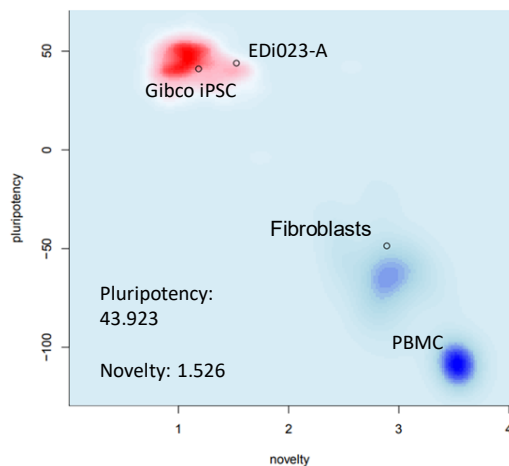


## D. G-Band karyotype



## E. hPSC Scorecard

iPSCs				Embryoid bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.28	-0.30	0.34	-0.67	-3.49	2.38	3.38	1.12


## F. Morphology 3 days post-thaw

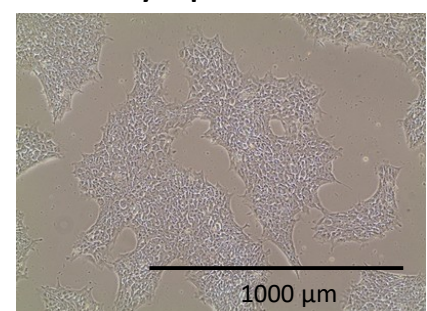






# sFig.4: Characterization for iPSC line EDi025-A


## A. AP

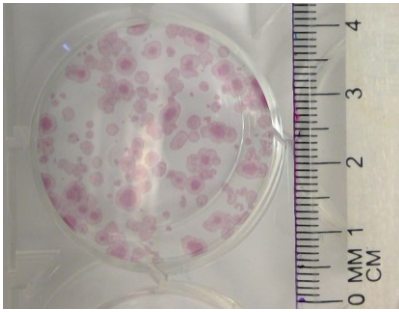


## B. Immunocytochemistry

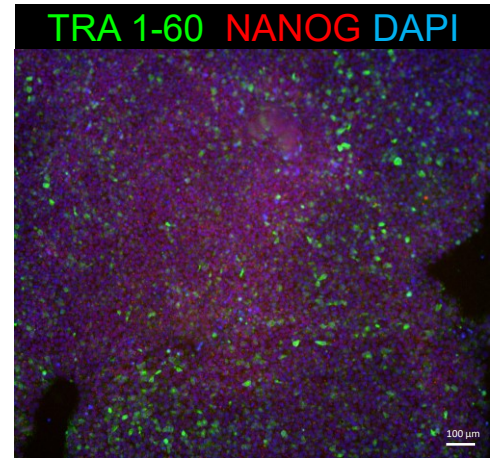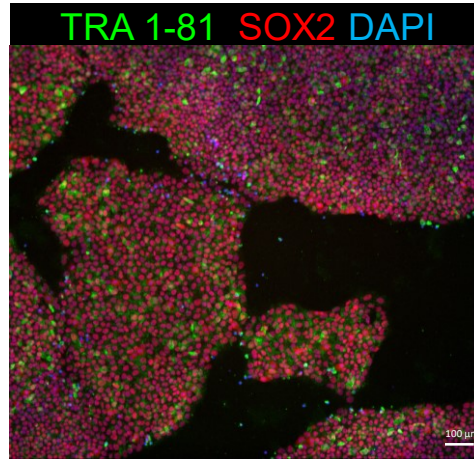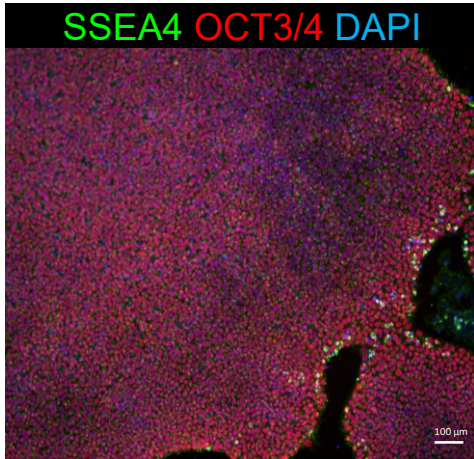

## C. PluriTest

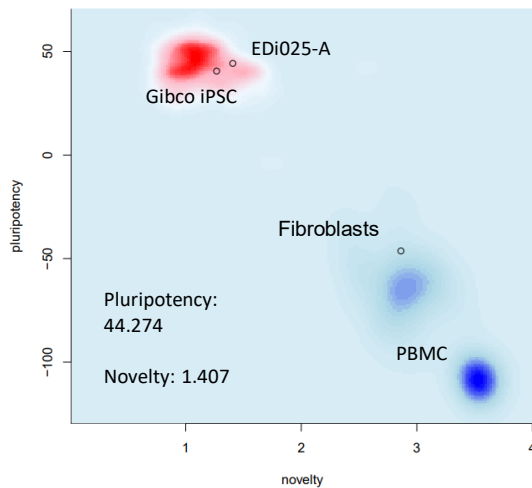

## D. G-Band karyotype

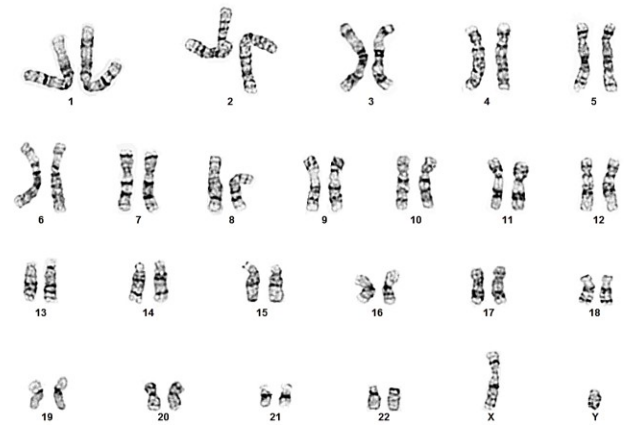


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.01	-0.33	-0.39	-1.15	-7.29	1.49	2.19	0.57


## F. Morphology

2 days post-thaw

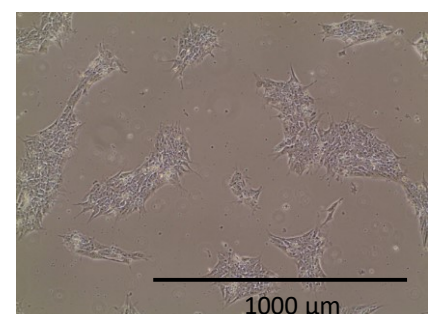


**sFig.5: Characterization for iPSC line EDi026-A**


**A. AP**

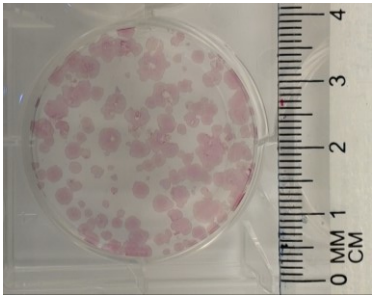


**B. Immunocytochemistry**

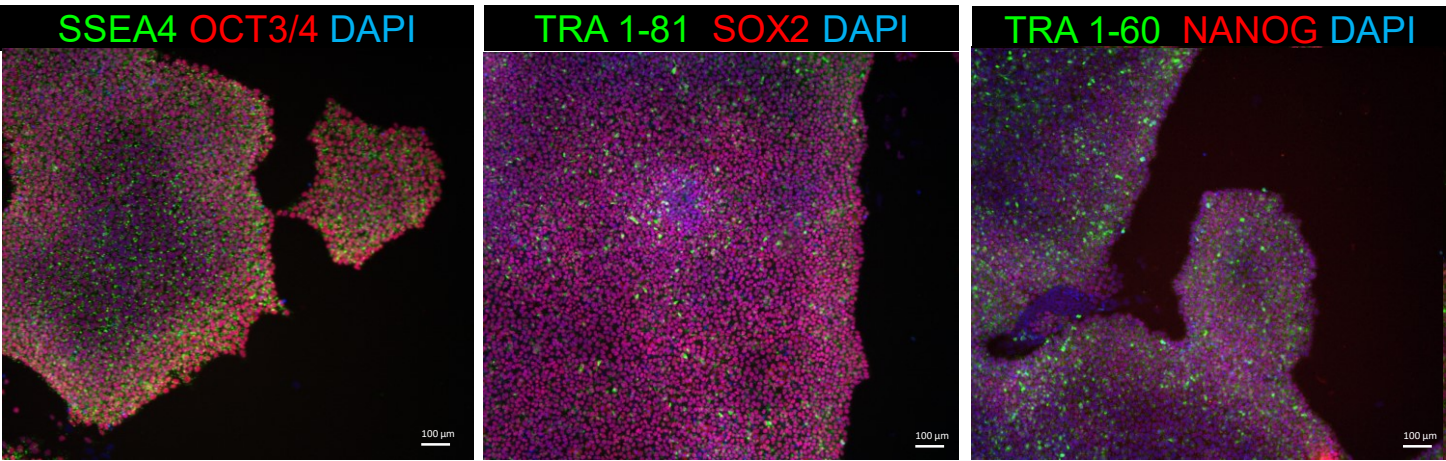

**C. Pluritest**

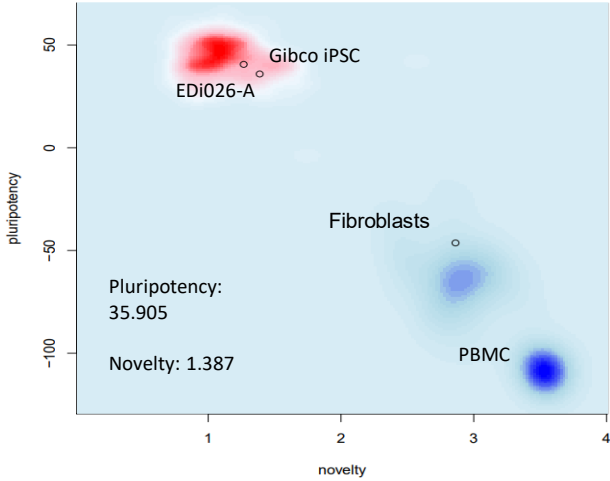


**D. G-Band karyotype**

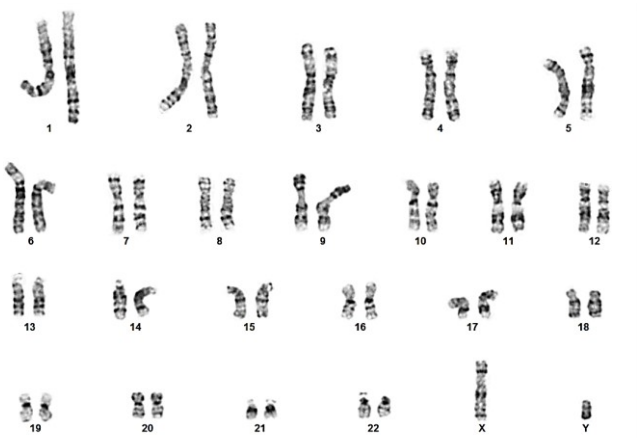

**E. hPSC Scorecard**

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.41	0.50	0.42	-0.83	-6.72	1.74	1.19	0.69


**F. Morphology**  
2 days post-thaw

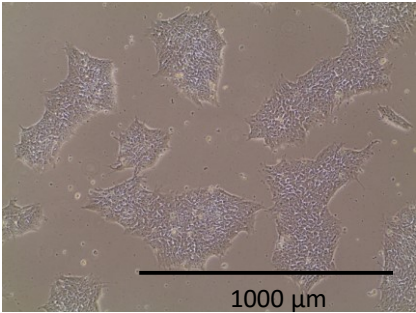




sFig.6: Characterization for iPSC line EDi027-A


A. AP

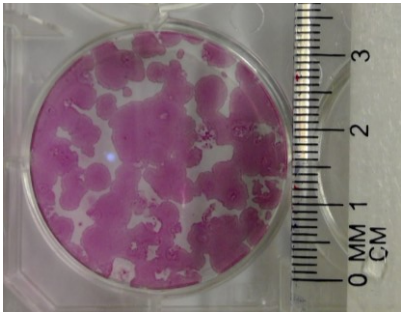

B. Immunocytochemistry

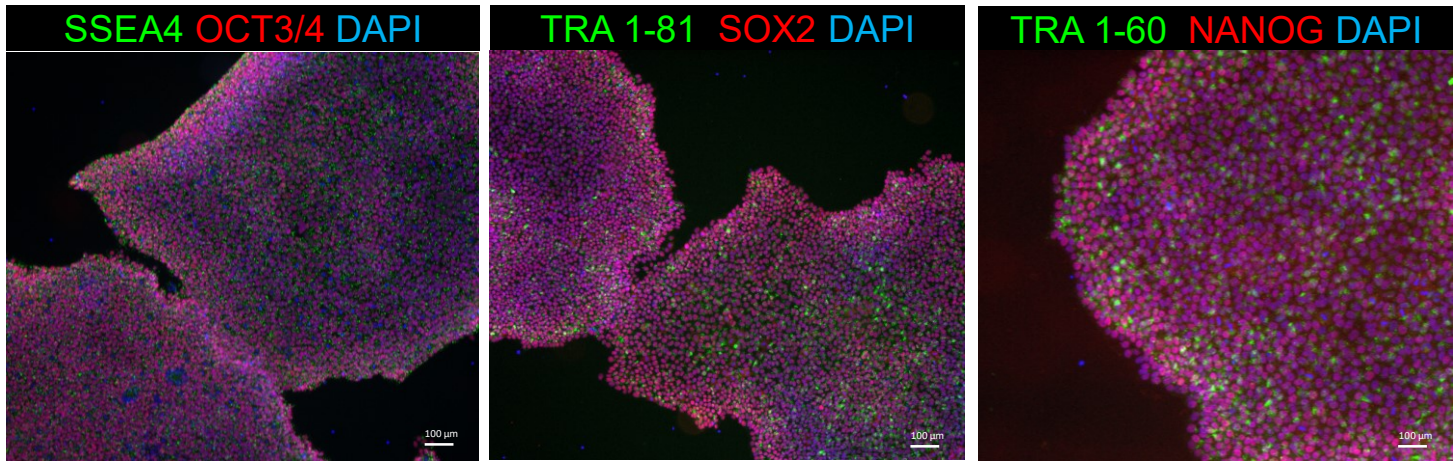

C. Pluritest

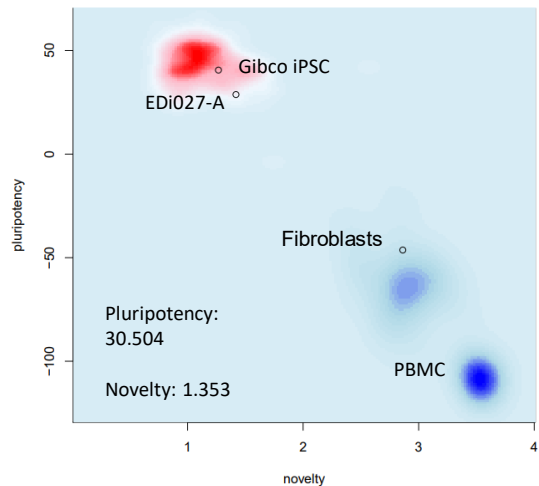


D. G-Band karyotype



E. hPSC Scorecard


iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	○	-	-	-	+	+	+
0.18	0.61	0.10	-0.75	-6.07	-1.72	-3.29	-1.78

F. Morphology  
11 days post-thaw



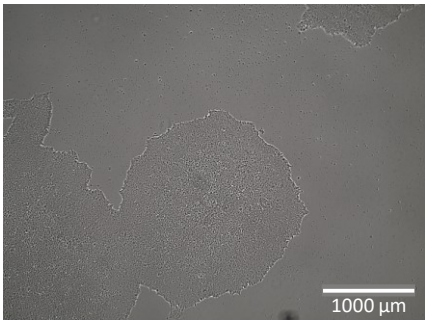


sFig.7: Characterization for iPSC line EDi028-A


A. AP

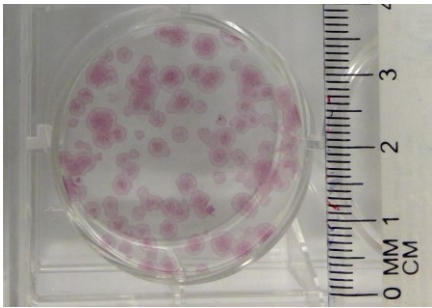


B. Immunocytochemistry

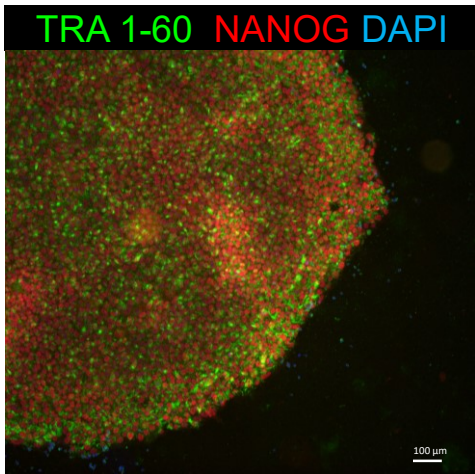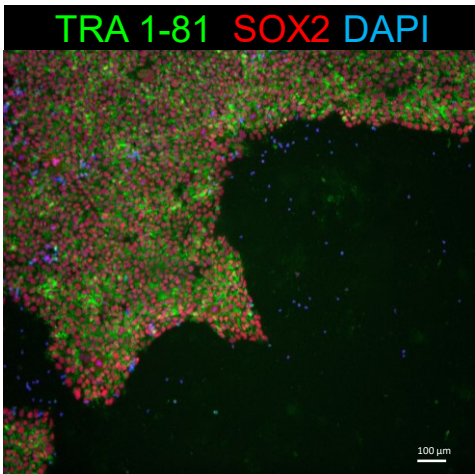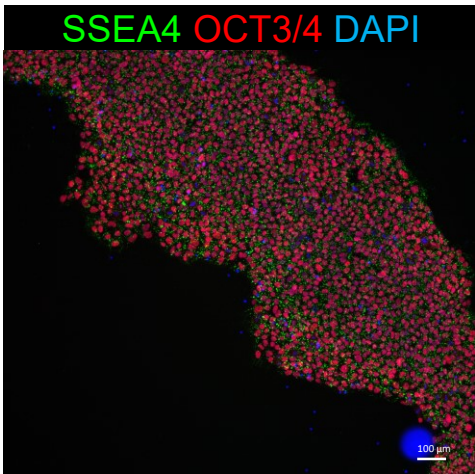

C. Pluritest

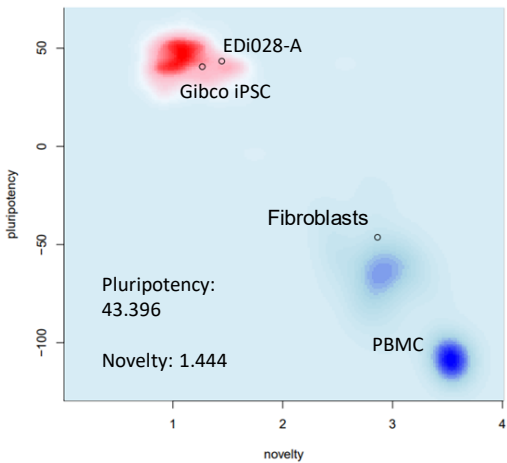


D. G -Band karyotype

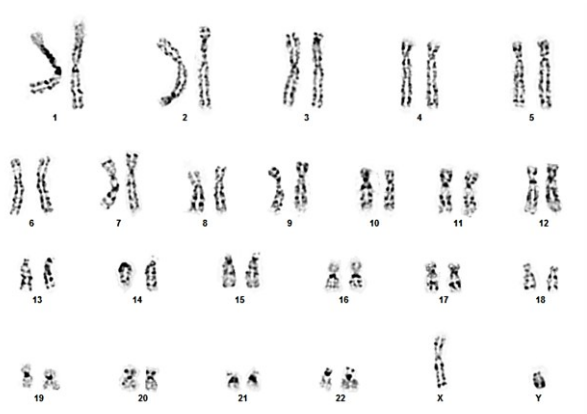

E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.17	-0.17	0.69	-1.00	-6.84	2.14	2.10	0.69


F. Morphology  
2 days post-thaw

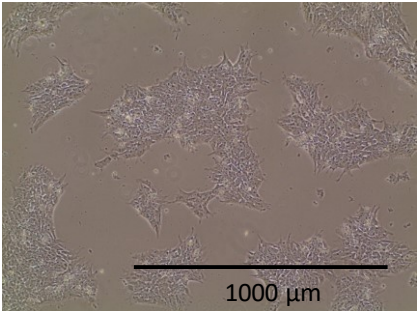






# sFig.8: Characterization for iPSC line EDi029-A


## A. AP

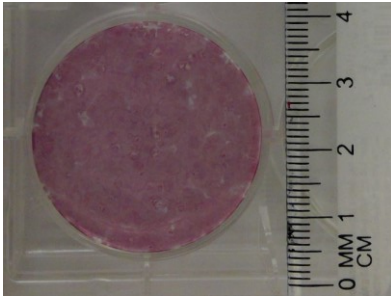


## B. Immunocytochemistry

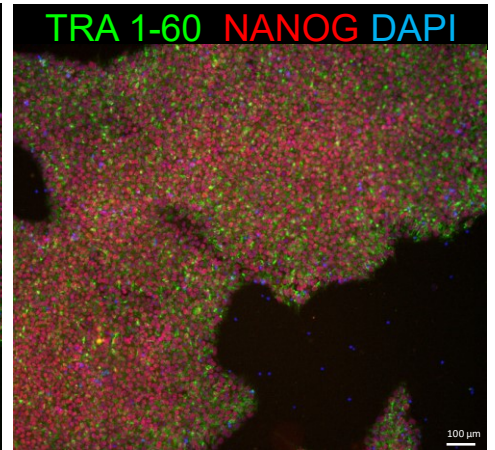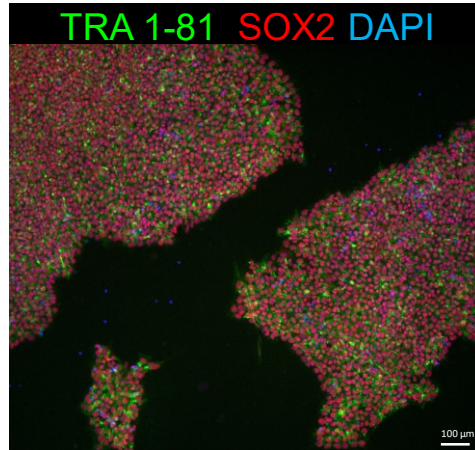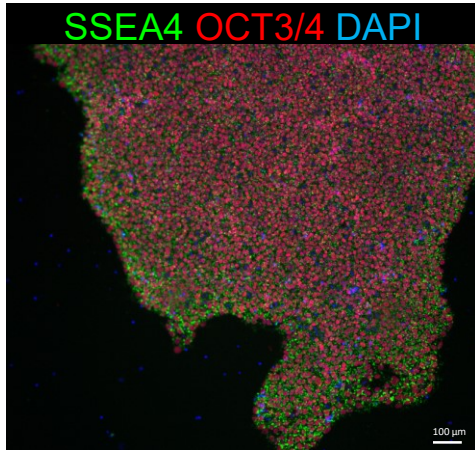

## C. Pluritest

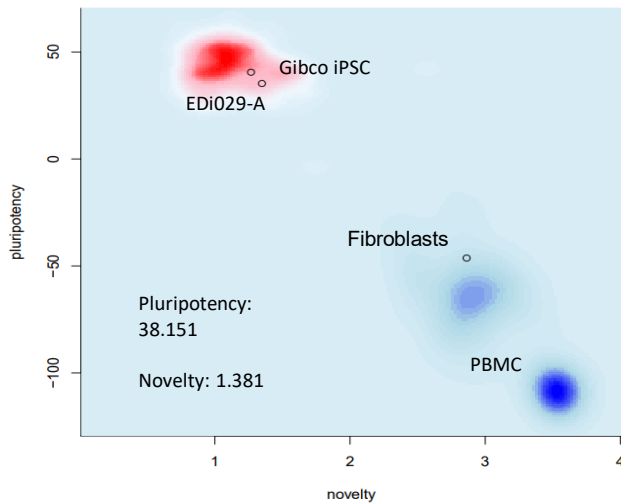

## D. G-Band karyotype

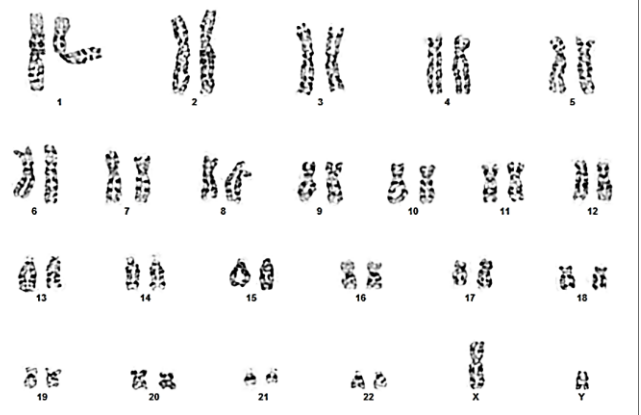


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.75	0.33	0.37	-0.81	-5.16	2.13	3.91	1.96

## F. Morphology


4 days post-thaw

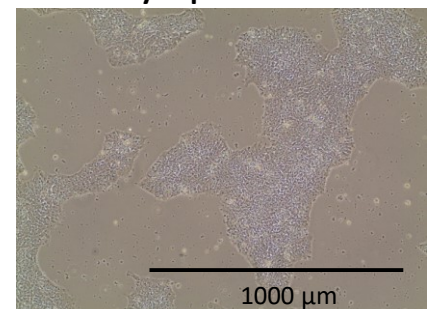






# sFig.9: Characterization for iPSC line EDi030-A


## A. AP

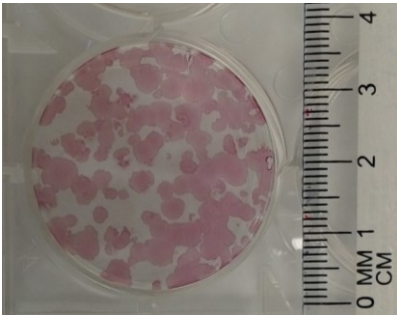


## B. Immunocytochemistry

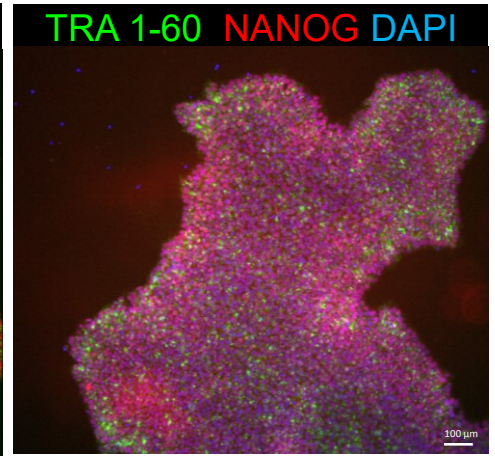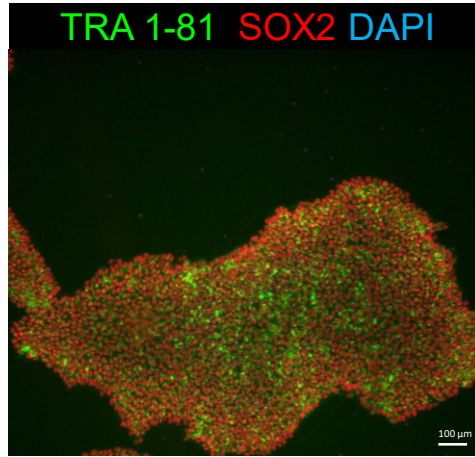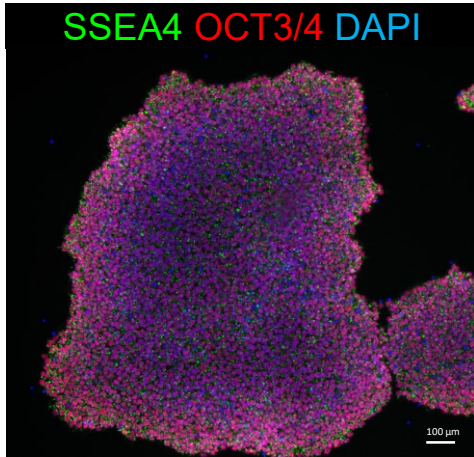

## C. Pluritest

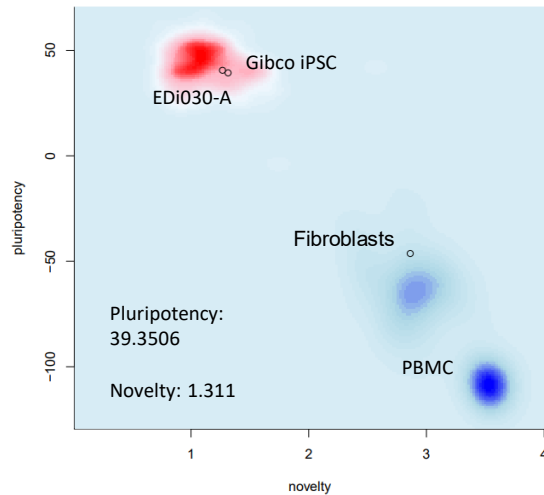

## D. G-Band karyotype

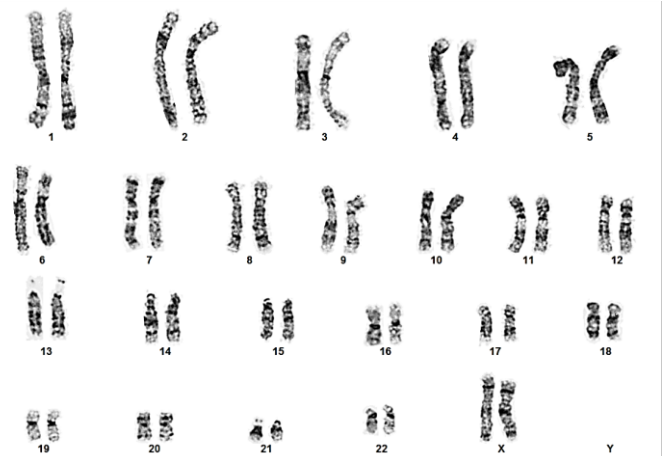


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.30	0.07	0.54	-1.48	-5.34	1.72	5.14	1.95


## F. Morphology

2 days post-thaw

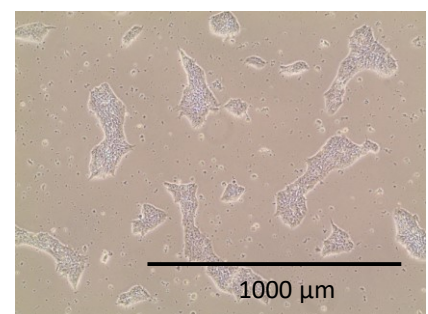


sFig.10: Characterization for iPSC line EDi031-A


A. AP

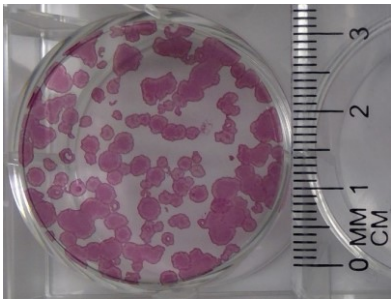


B. Immunocytochemistry

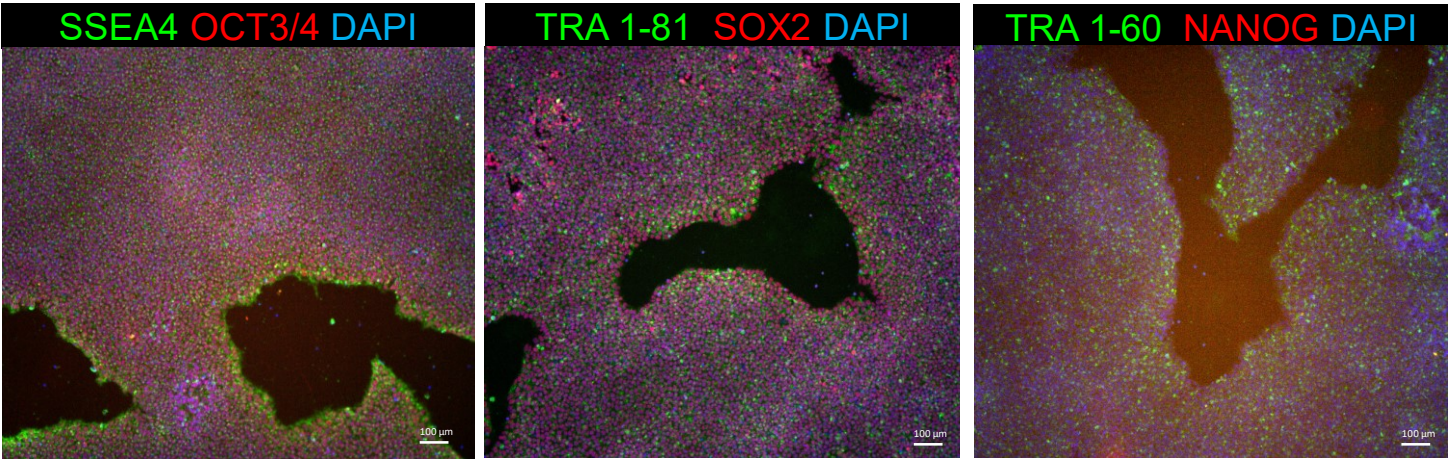

C. Pluritest

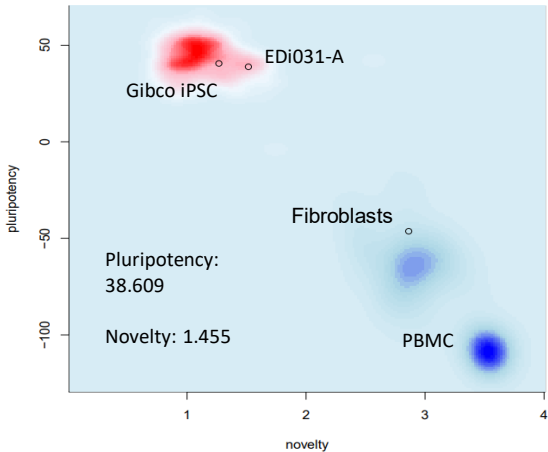

D. G-Band karyotype

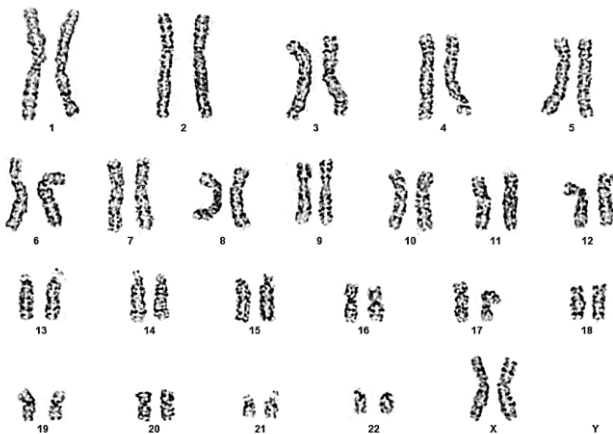


E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.19	-0.46	-0.24	-1.30	-6.33	1.77	2.65	1.65

F. Morphology


2 days post-thaw

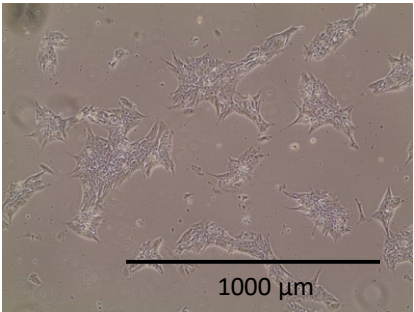






# sFig.11: Characterization for iPSC line EDi032-A


## A. AP

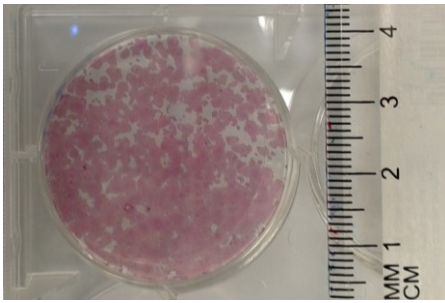


## B. Immunocytochemistry

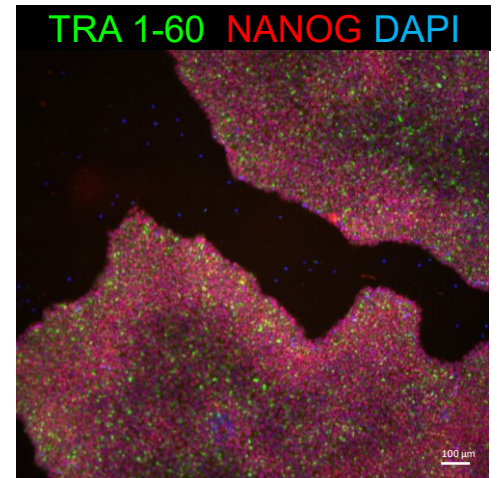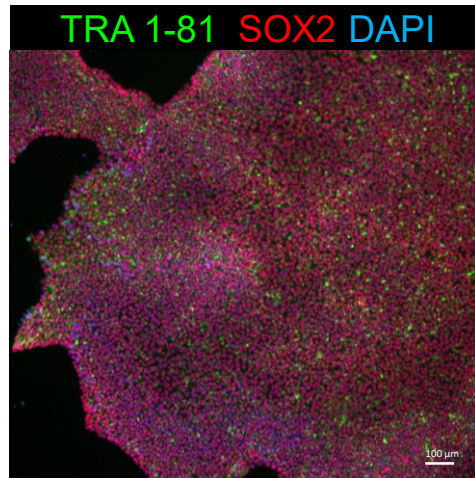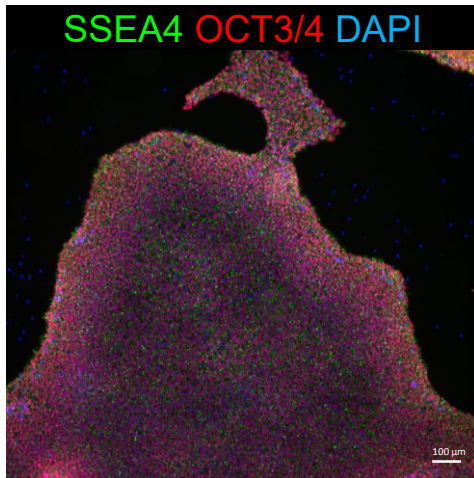

## C. Pluritest

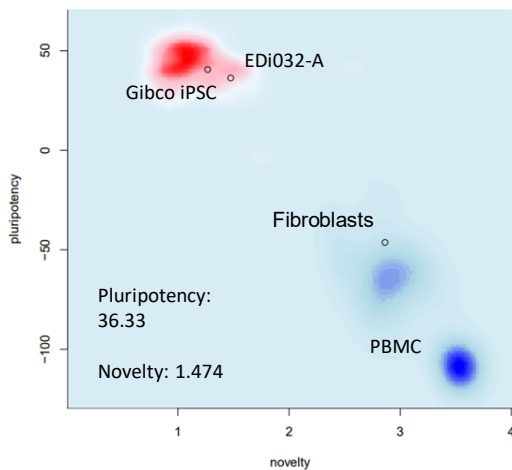

## D. G -Band karyotype

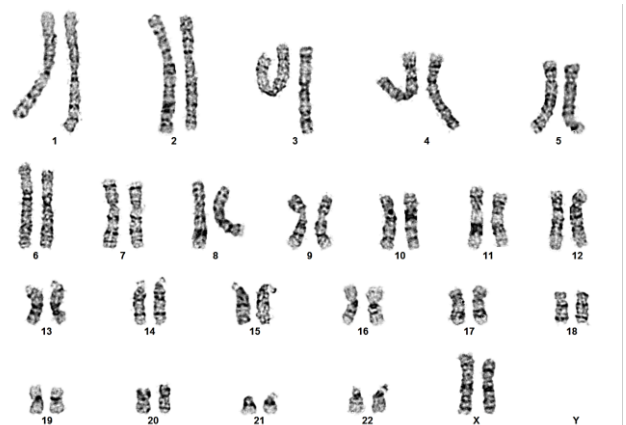


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.24	-0.06	0.08	-1.36	-6.85	1.85	2.34	0.78


## F. Morphology

4 days post-thaw



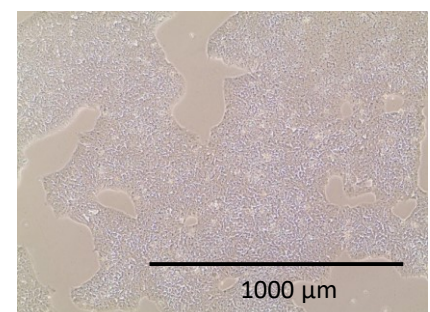


**sFig.12: Characterization for iPSC line EDi033-A**


**A. AP**

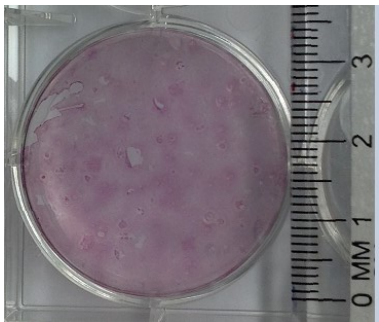

**B. Immunocytochemistry**

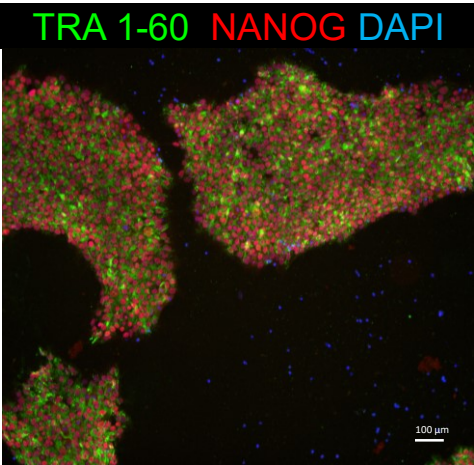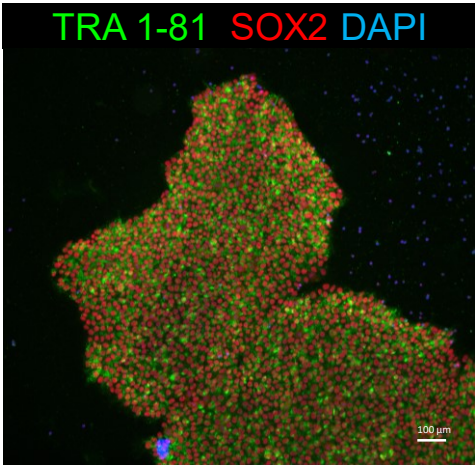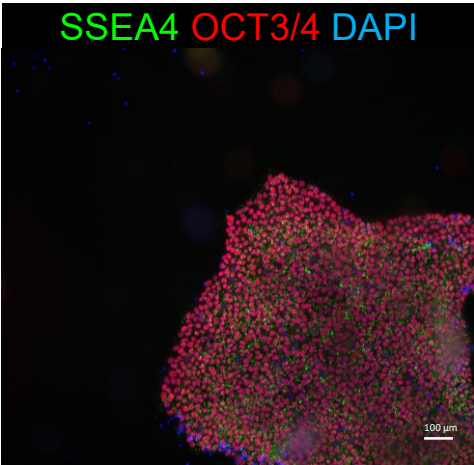

**C. Pluritest**

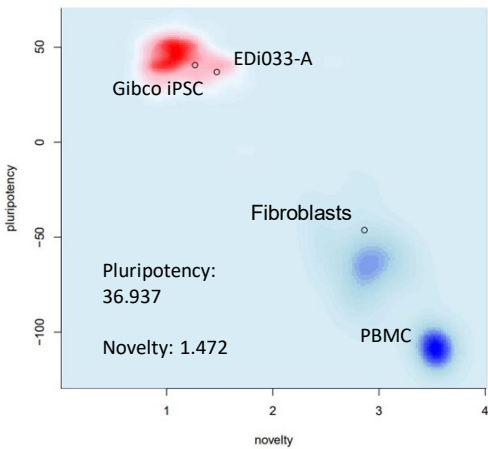

**D. G -Band karyotype**




**E. hPSC Scorecard**

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.23	-0.02	-0.18	-1.03	-3.58	1.77	3.93	1.85

**F. Morphology**


7 days post-thaw

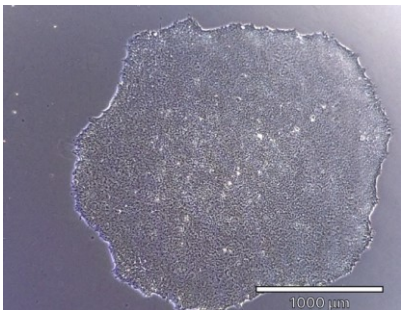






# sFig.13: Characterization for iPSC line EDi034-A


## A. AP

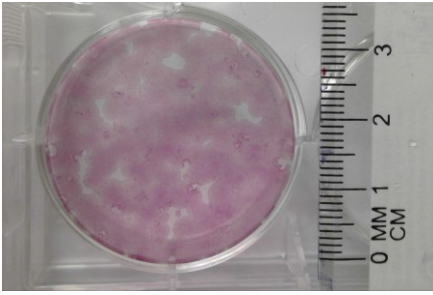


## B. Immunocytochemistry

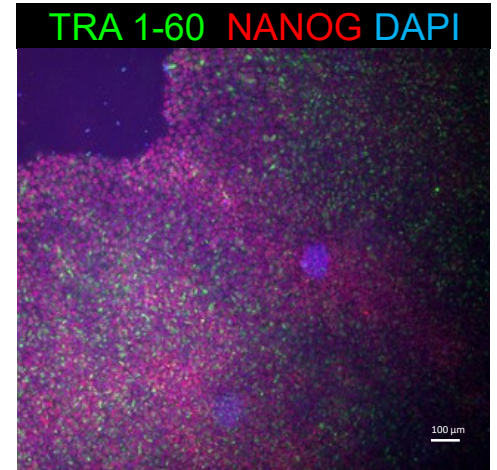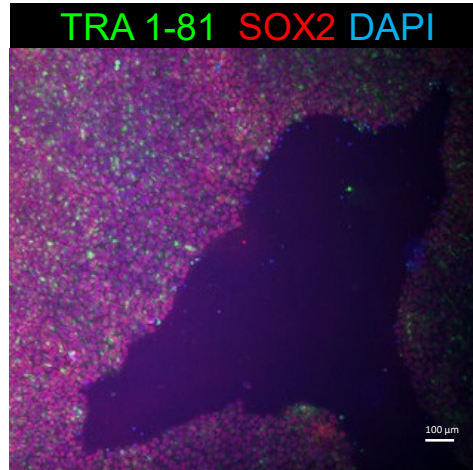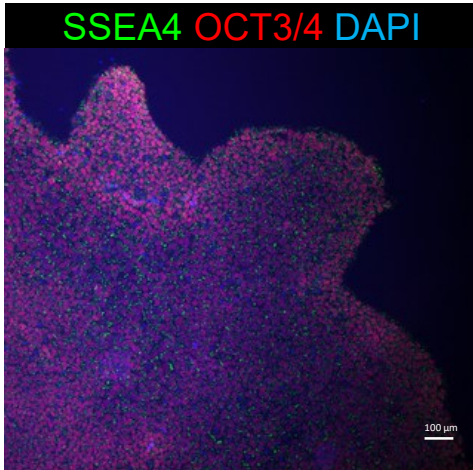

## C. Pluritest

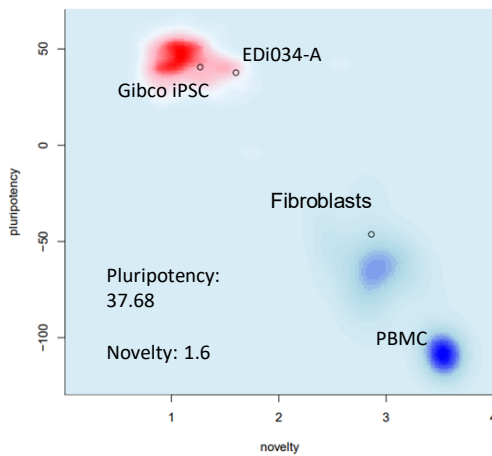


## D. G -Band karyotype

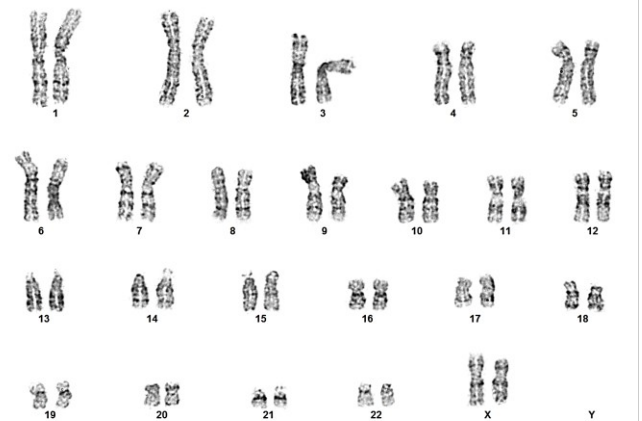

## E. hPSC Scorecard


iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.77	-0.25	-0.22	-1.57	-4.24	2.56	3.43	1.67

## F. Morphology 2 days post-thaw



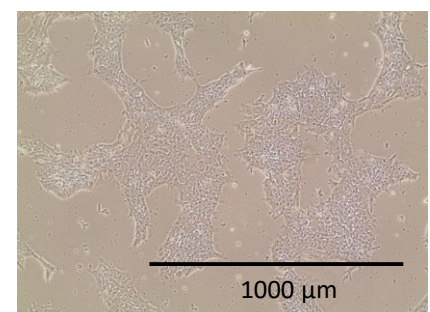


# sFig.14: Characterization for iPSC line EDi035-A


## A. AP

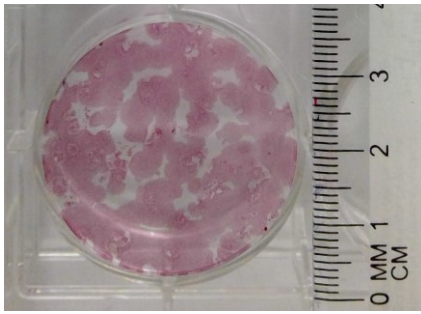


## B. Immunocytochemistry

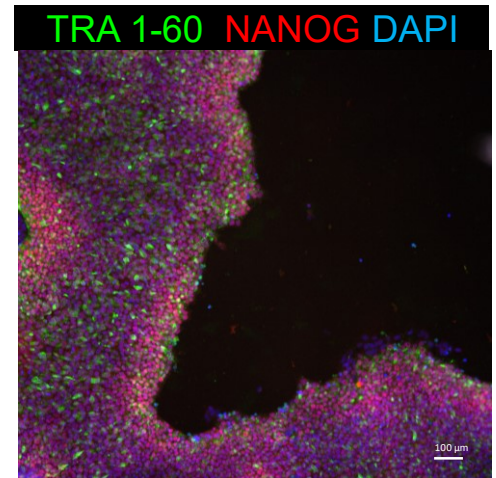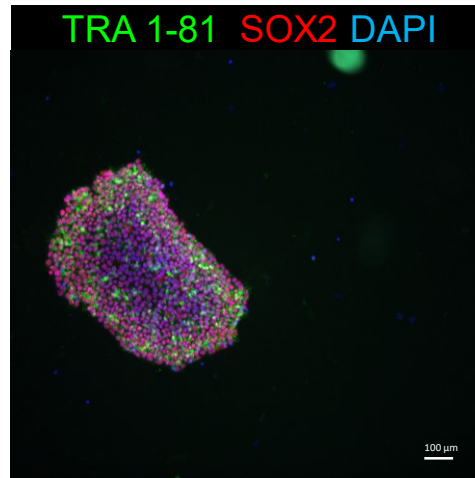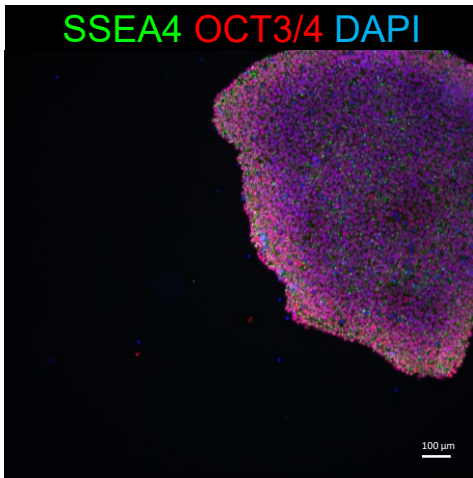

## C. Pluritest

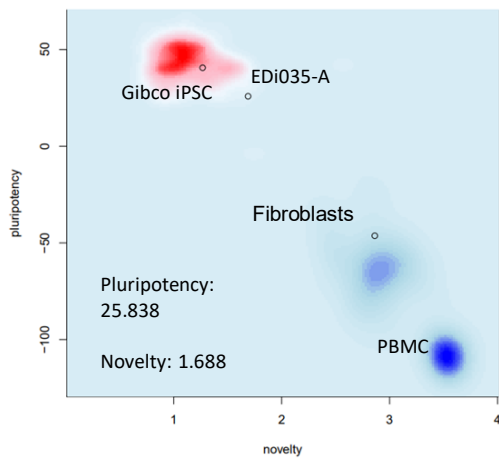


## D. G -Band karyotype

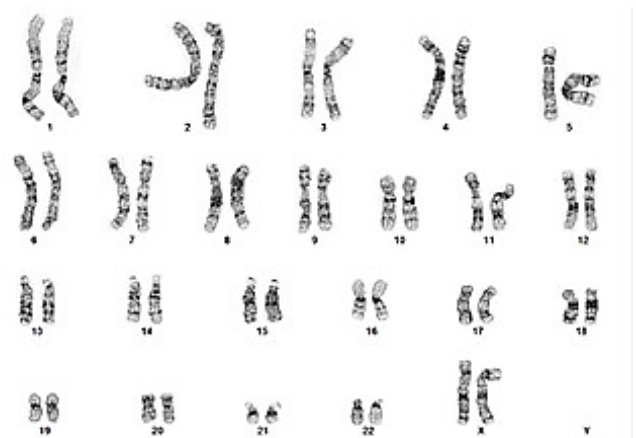

## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	-	0
0.39	0.19	0.02	-1.22	-6.69	1.81	0.64	0.29


## F. Morphology 7 days post-thaw

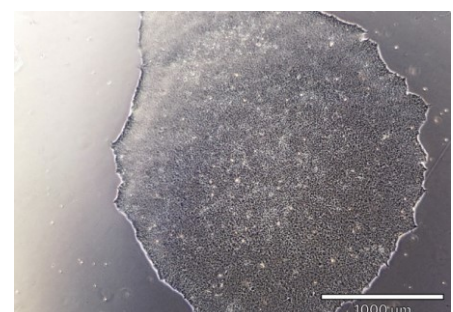






**sFig.15: Characterization for iPSC line EDi036-A**


**A. AP**

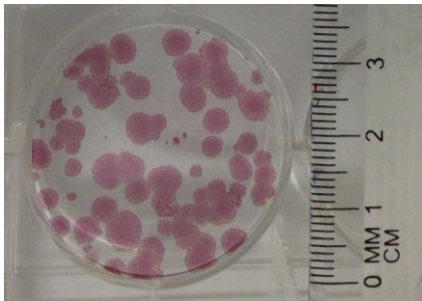


**B. Immunocytochemistry**

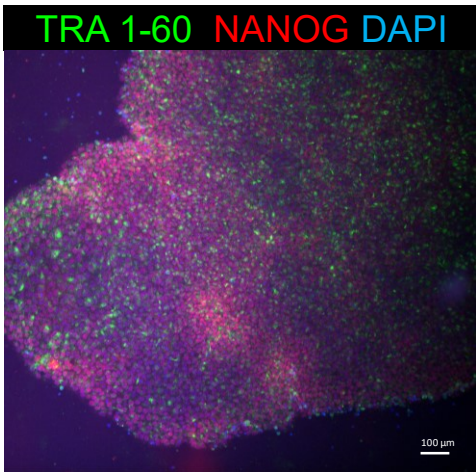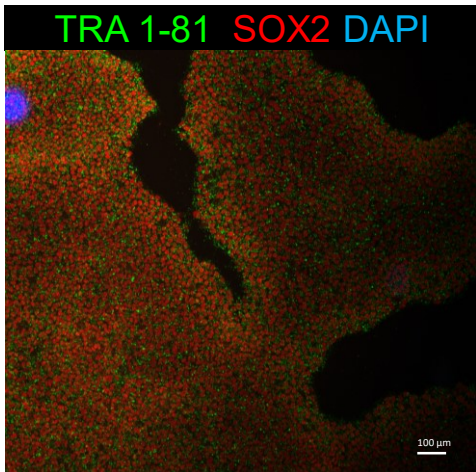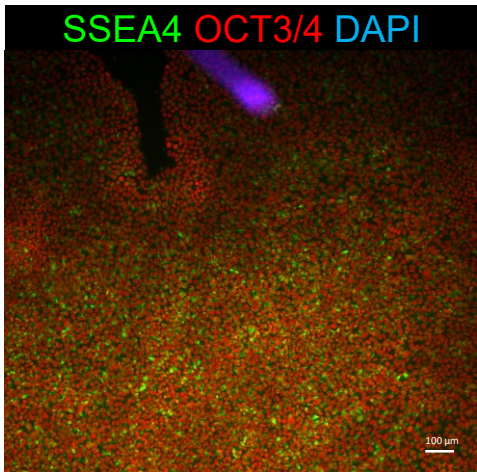

**C. Pluritest**

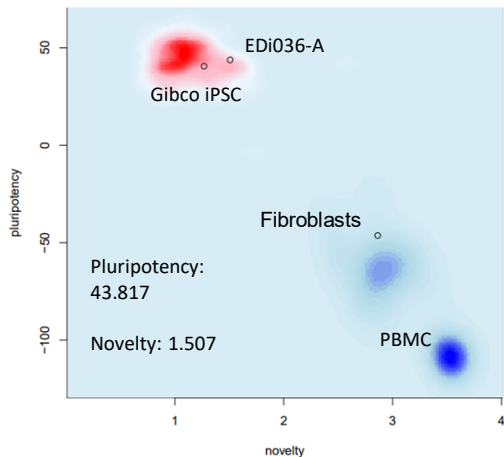

**D. G -Band karyotype**

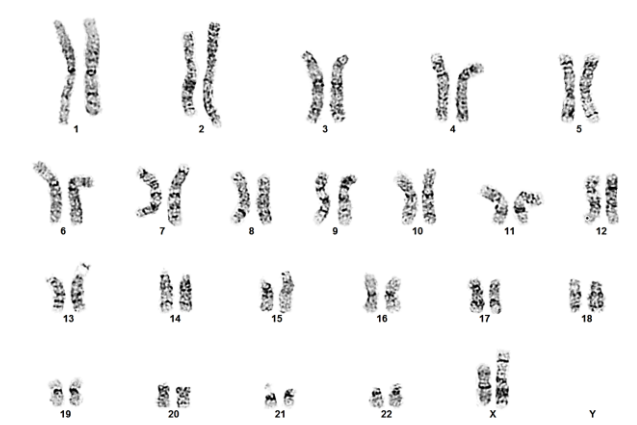


**E. hPSC Scorecard**

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.50	0.37	0.20	-1.03	-6.35	1.90	1.61	0.66


**F. Morphology**

2 days post-thaw



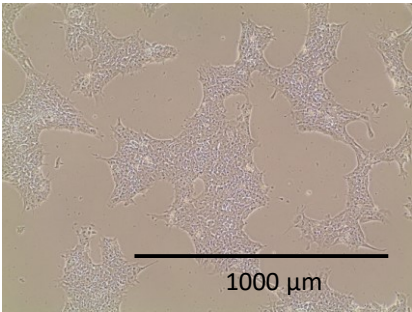


**sFig.16: Characterization for iPSC line EDi037-A**


**A. AP**

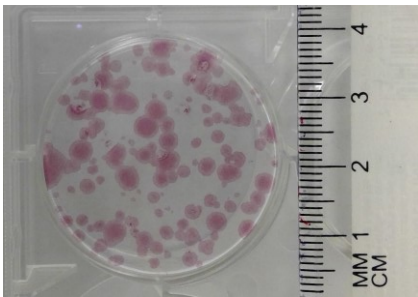

**B. Immunocytochemistry**

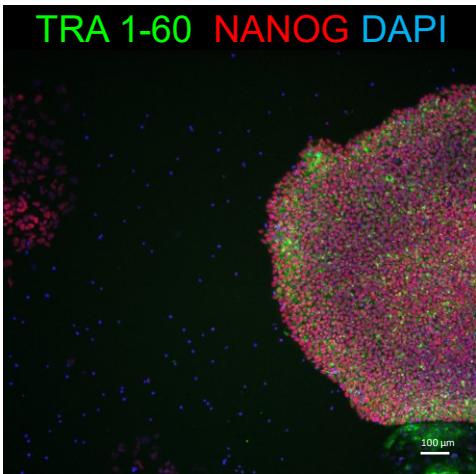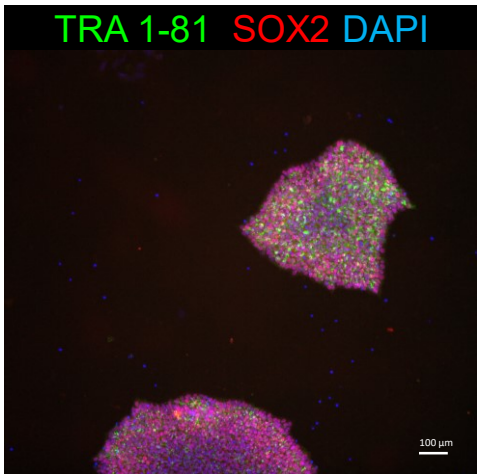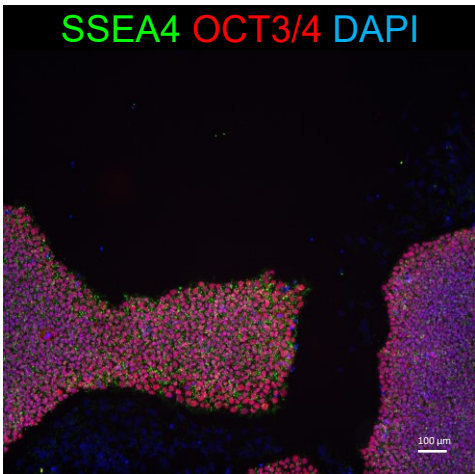

**C. Pluritest**

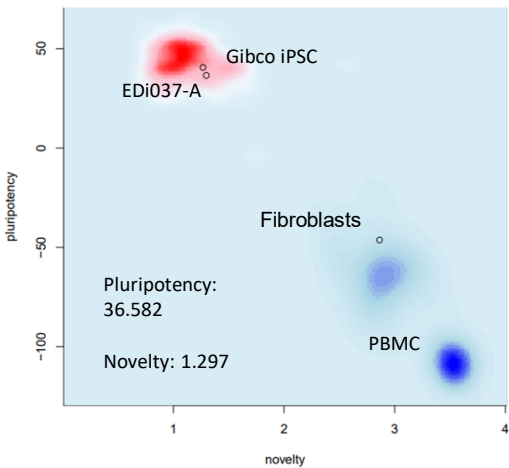

**D. G -Band karyotype**




**E. hPSC Scorecard**

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.05	0.20	0.05	-0.75	-7.38	1.65	1.80	0.87

**F. Morphology**


2 days post-thaw

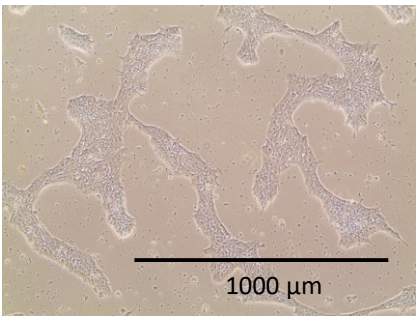






# sFig.17: Characterization for iPSC line EDi038-A


## A. AP

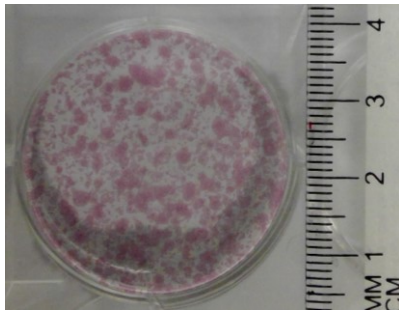


## B. Immunocytochemistry

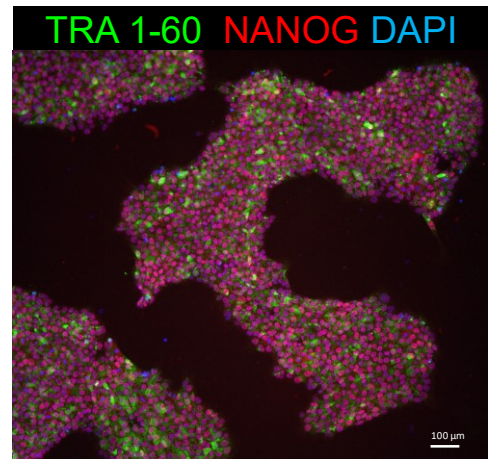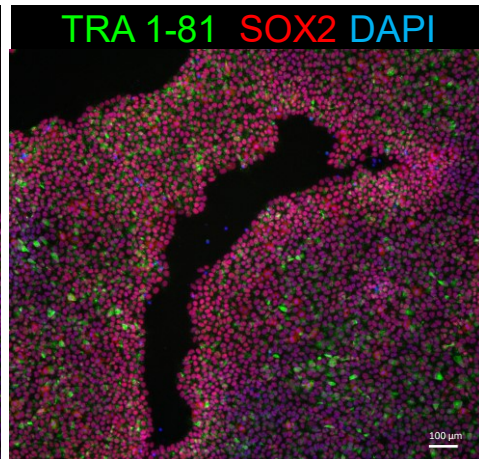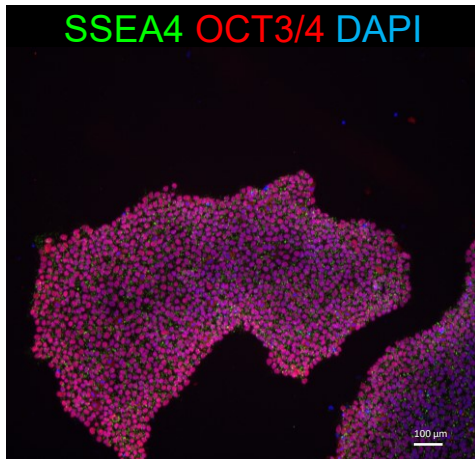

## C. Pluritest

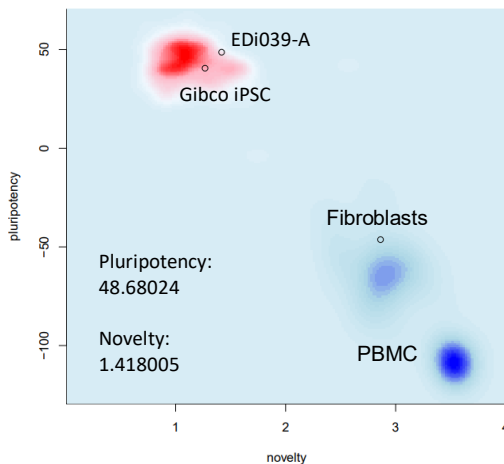

## D. G-Band karyotype

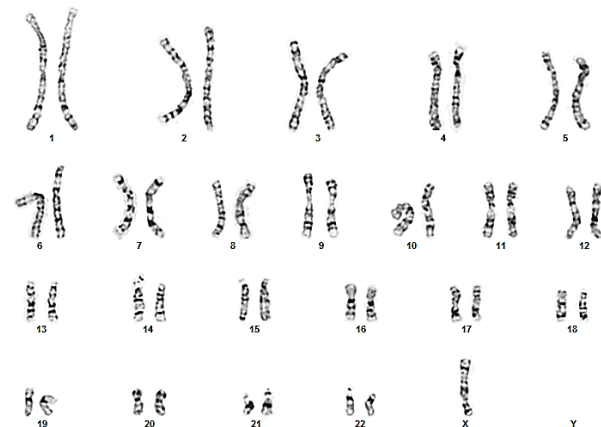


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.61	-0.09	0.48	-1.10	-6.22	1.41	6.00	1.73


## F. Morphology

11 days post-thaw

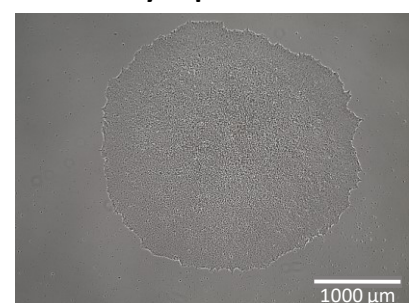


# sFig.18: Characterization for iPSC line EDi039-A


## A. AP

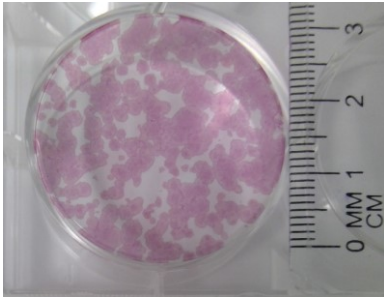


## B. Immunocytochemistry

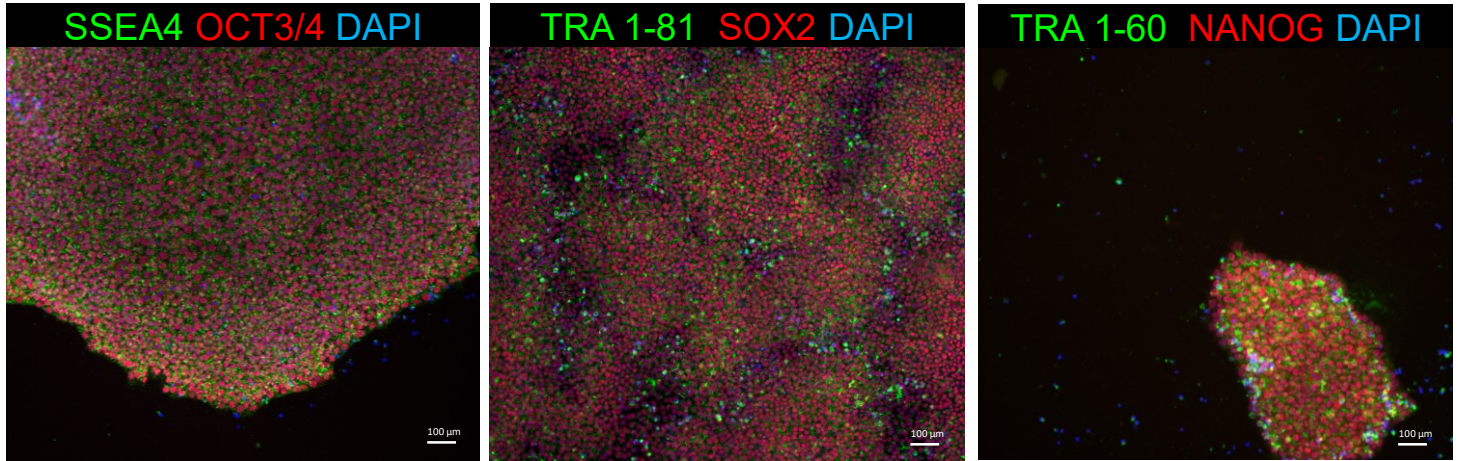

## C. PluriTest

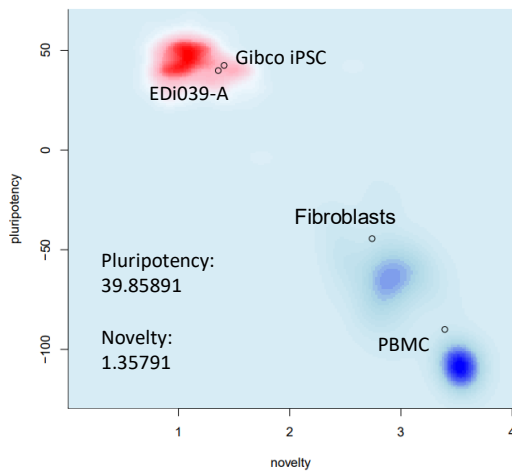


## D. G-Band karyotype

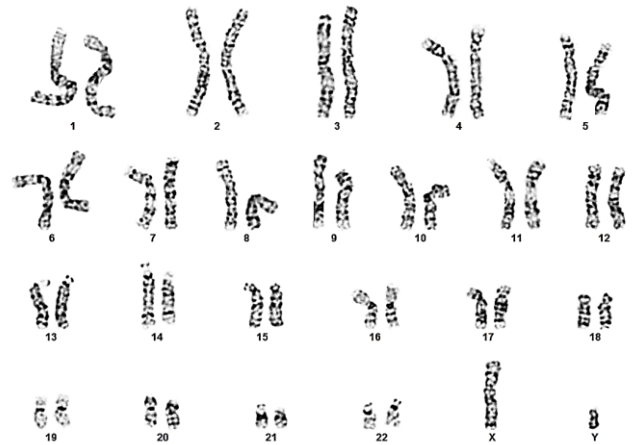

## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.19	0.14	-0.15	-1.05	-6.74	2.08	1.07	0.64


## F. Morphology 11 days post-thaw

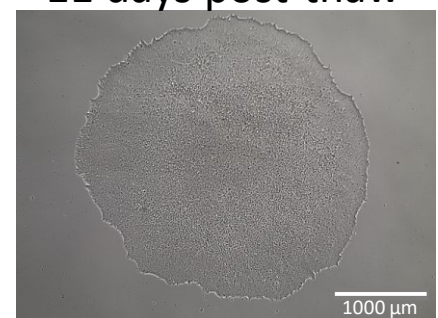






# sFig.19: Characterization for iPSC line EDi040-A


## A. AP

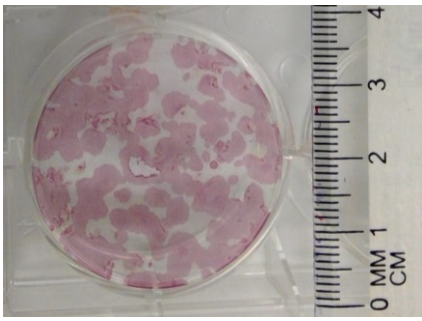


## B. Immunocytochemistry

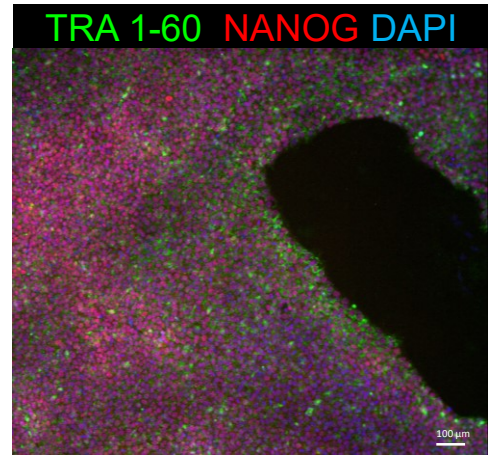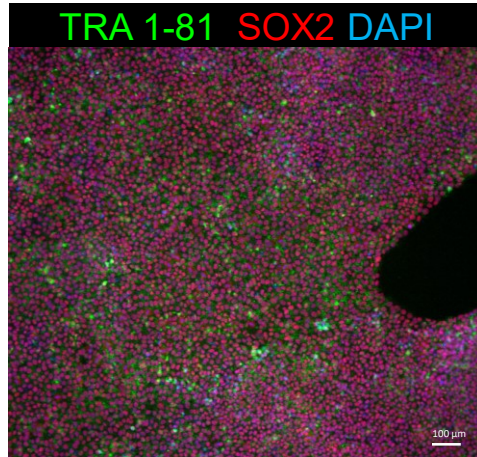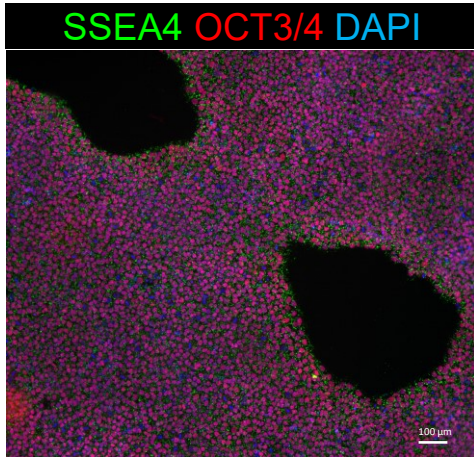

## C. Pluritest

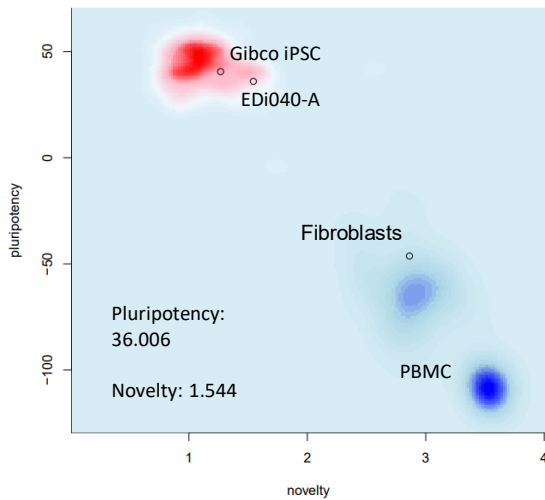

## D. G-Band karyotype

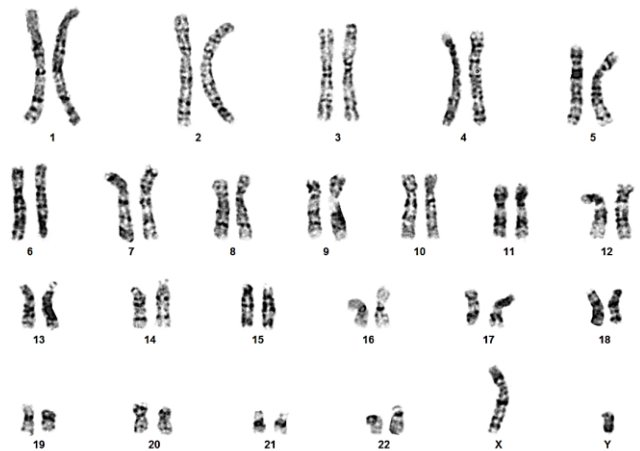


## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
0.01	-0.45	0.24	-0.94	-3.18	1.57	0.98	-0.56


## F. Morphology

2 days post-thaw

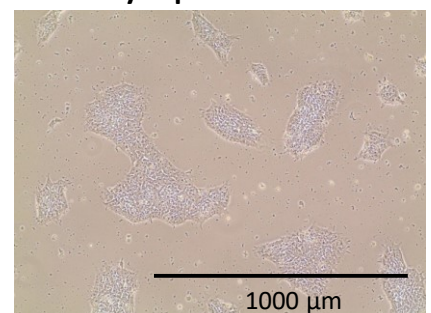


sFig.20: Characterization for iPSC line EDi041-A


A. AP

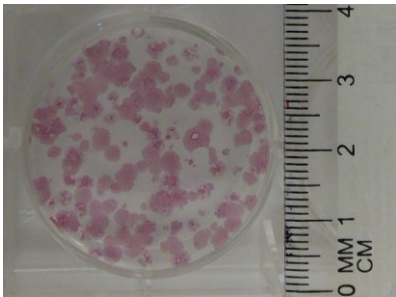


B. Immunocytochemistry

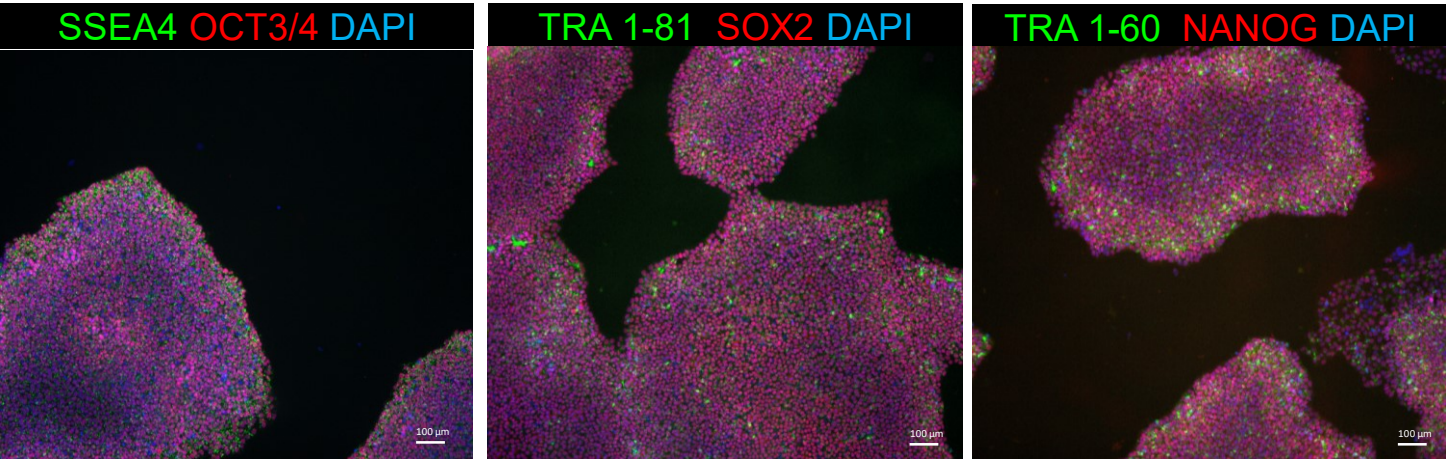

C. Pluritest

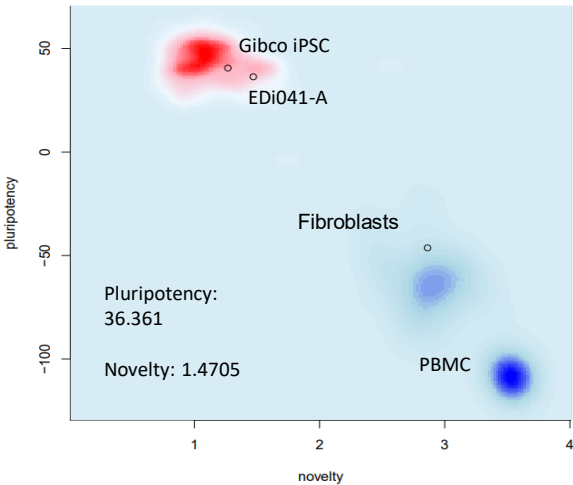

D. G-Band karyotype

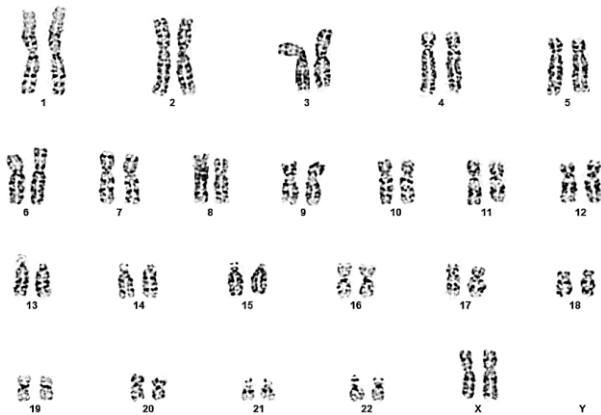


E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.18	0.05	-0.32	-1.24	-7.04	1.61	1.42	0.67

F. Morphology


3 days post-thaw

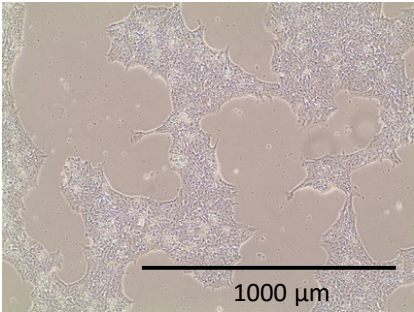






# sFig.21: Characterization for iPSC line EDi042-A


## A. AP

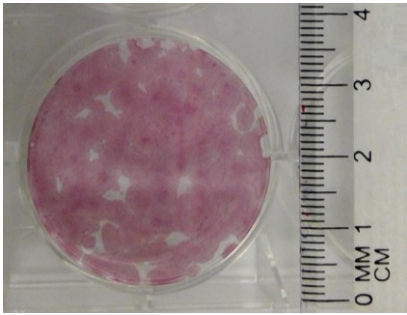

## B. Immunocytochemistry

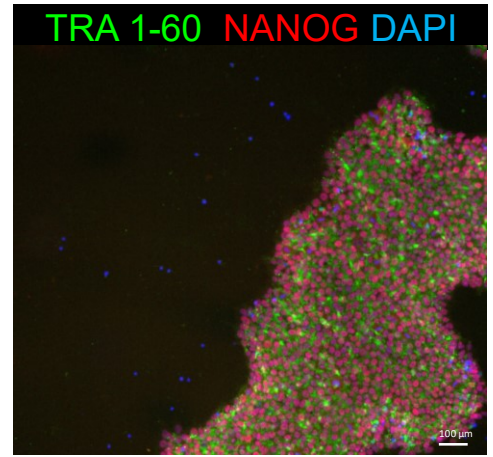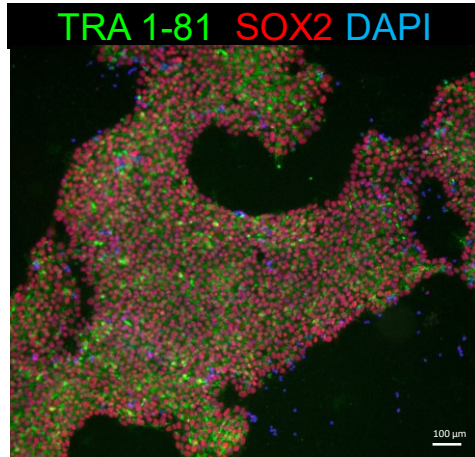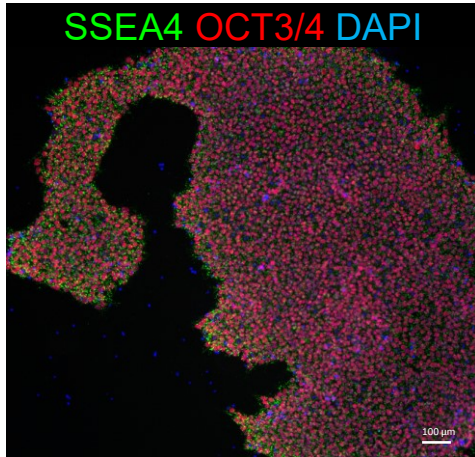

## C. Pluritest

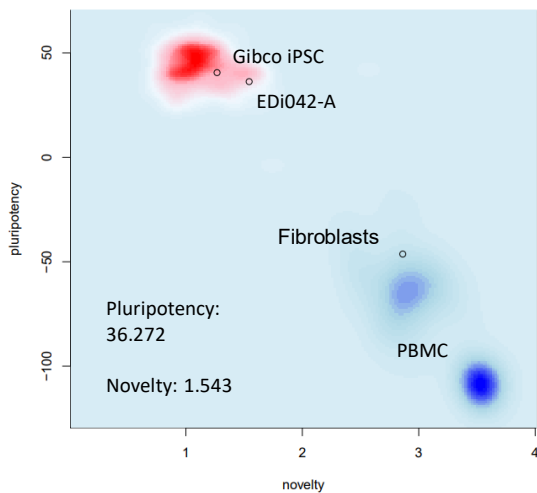

## D. G-Band karyotype




## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	-
0.05	-0.13	0.05	-1.04	-5.15	2.62	1.07	0.22


## F. Morphology

4 days post-thaw



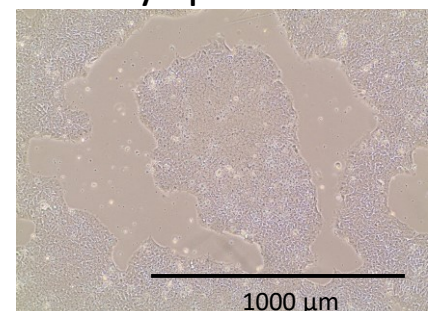


# sFig.22: Characterization for iPSC line EDi043-A


## A. AP

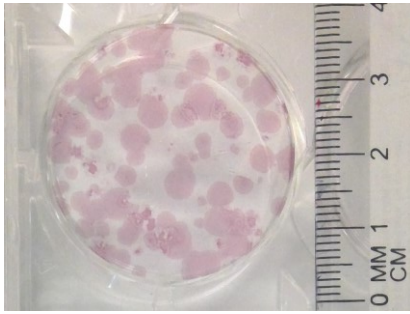


## B. Immunocytochemistry

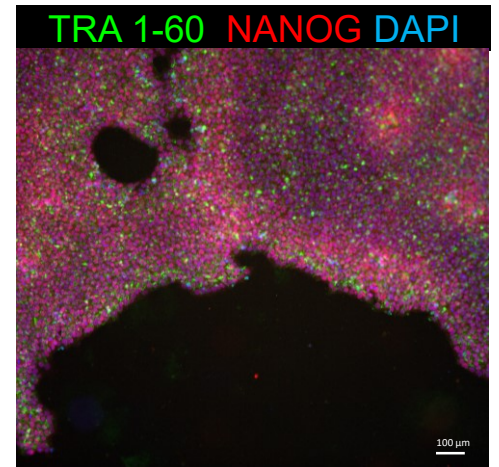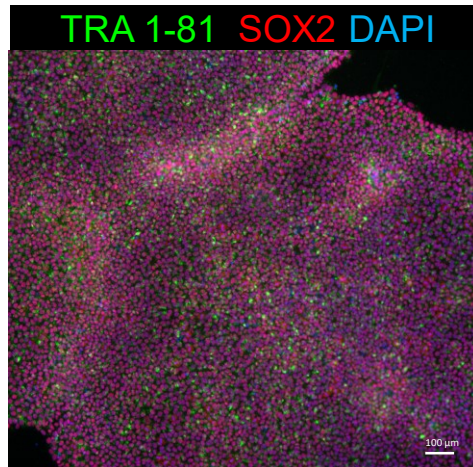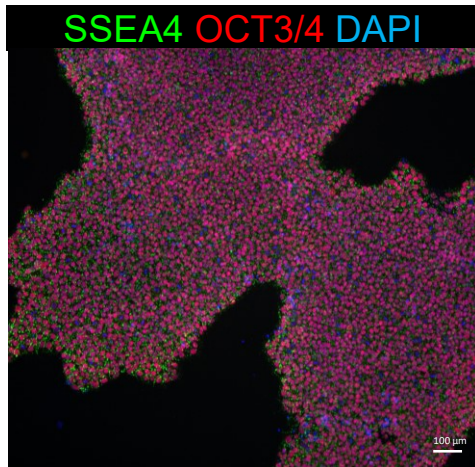

## C. Pluritest

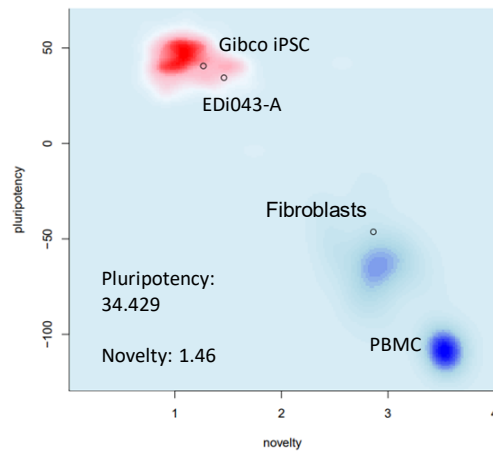


## D. G -Band karyotype

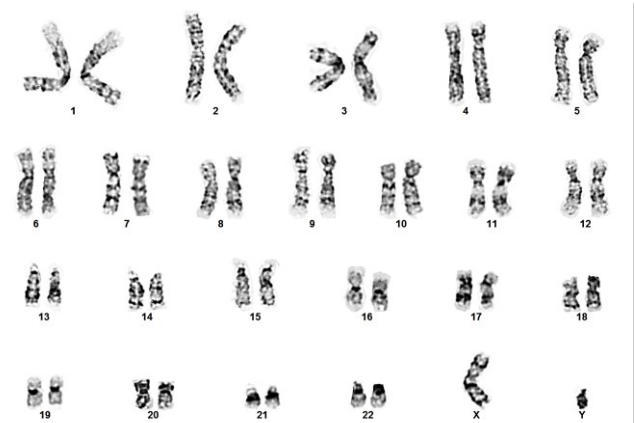

## E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.03	0.38	0.07	-1.23	-6.46	1.59	2.57	1.19


## F. Morphology 7 days post-thaw

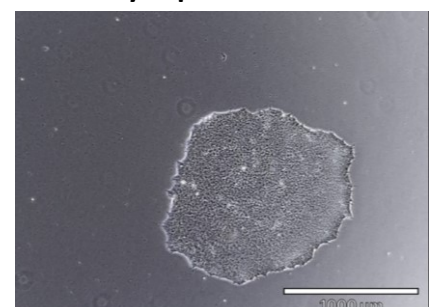




sFig.23: Characterization for iPSC line EDi044-A


A. AP

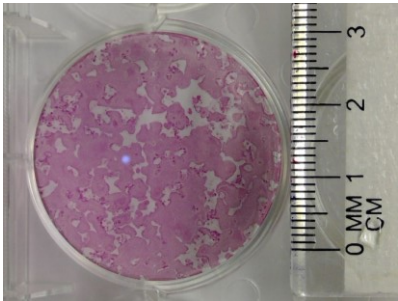


B. Immunocytochemistry

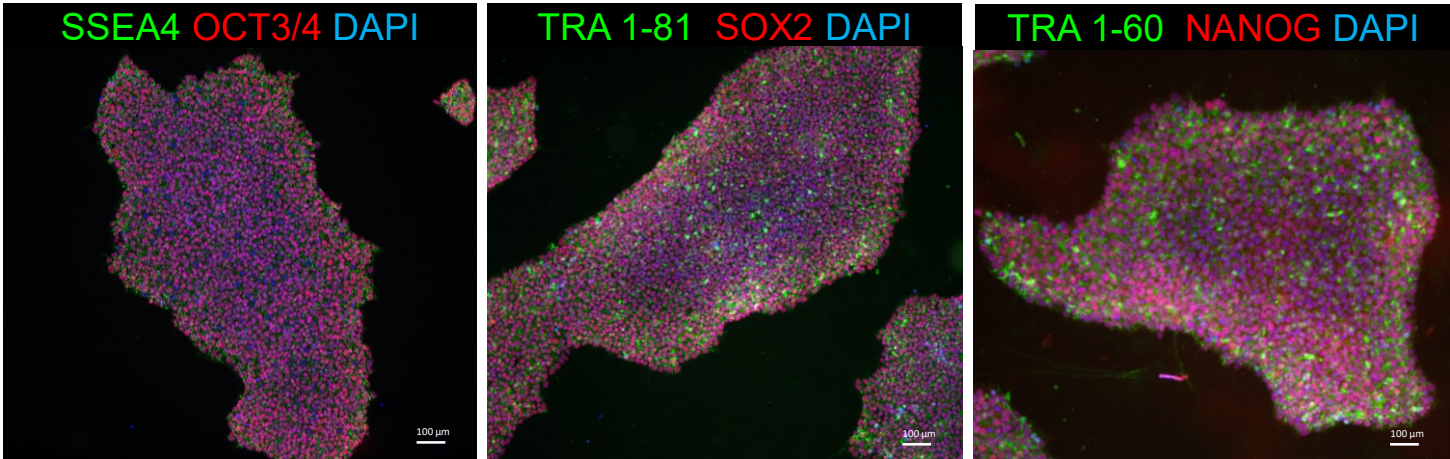

C. Pluritest

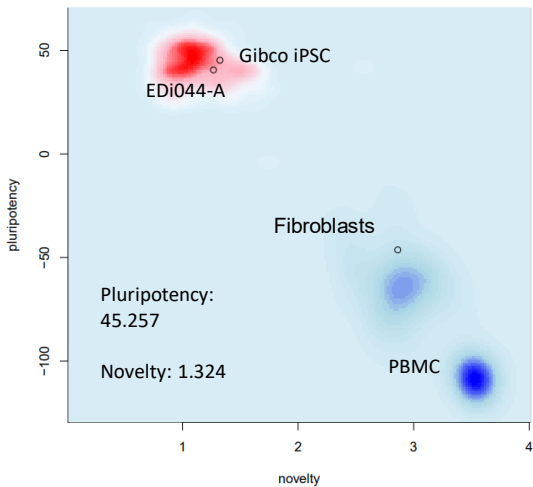

D. G-Band karyotype

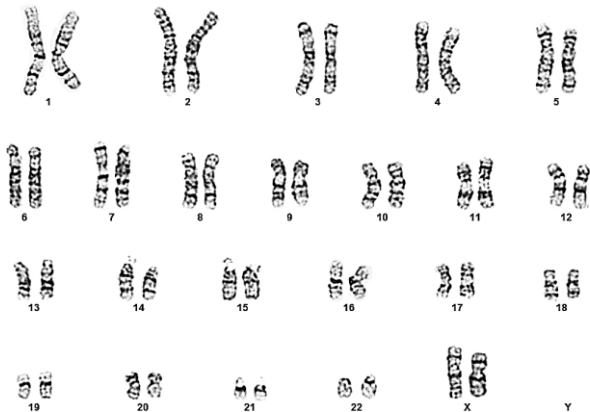


E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
○	⊖	⊖	⊖	⊖	⊕	⊕	⊕
-0.75	-0.19	0.58	-0.80	-6.29	1.65	3.42	0.94


F. Morphology

4 days post-thaw

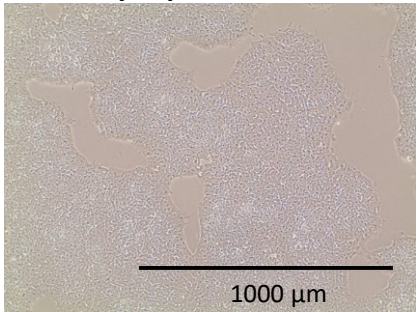


sFig.24: Characterization for iPSC line EDi045-A


A. AP

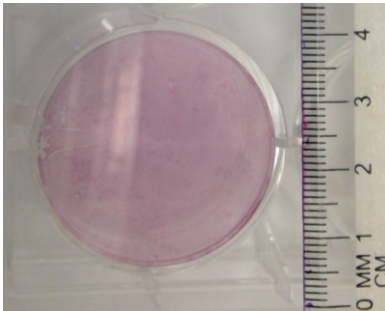


B. Immunocytochemistry

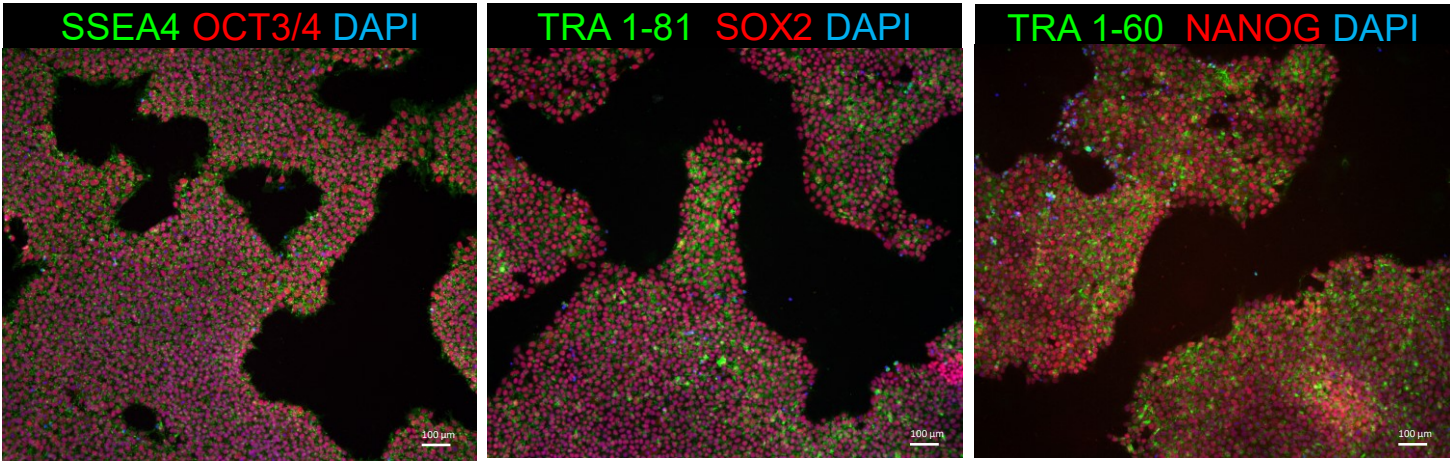

C. Pluritest

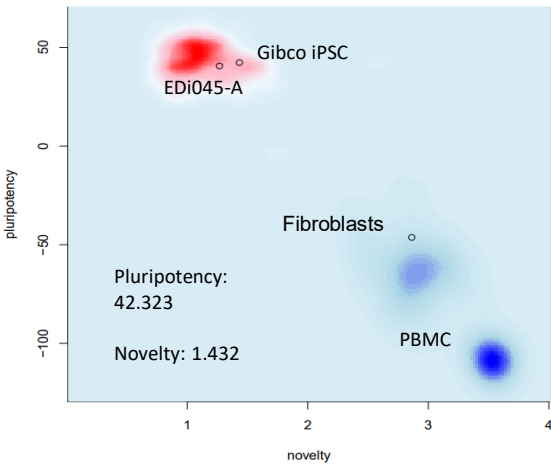

D. G-Band karyotype

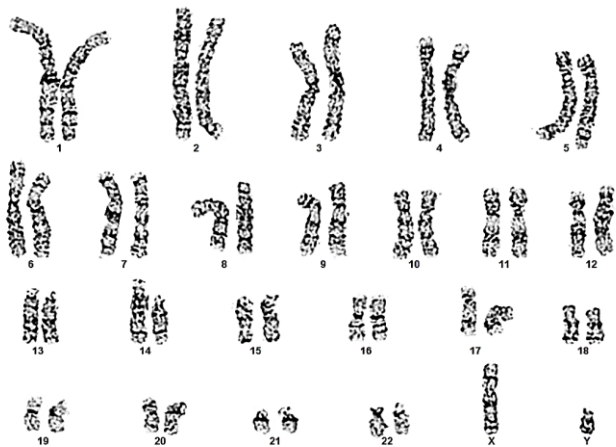


E. hPSC Scorecard

iPSCs				Embryoid Bodies			
Self-renew	Ecto	Meso	Endo	Self-renew	Ecto	Meso	Endo
+	-	-	-	-	+	+	+
-0.11	0.15	0.22	-0.81	-6.43	2.29	1.53	0.62


F. Morphology

14 days post-thaw

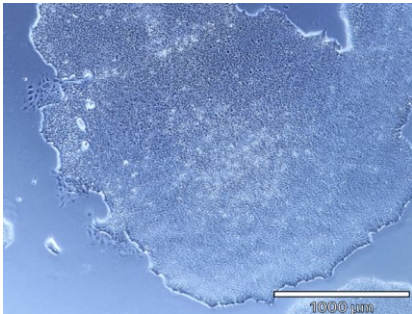









## sFig.25: EBNA assay for LBC lines

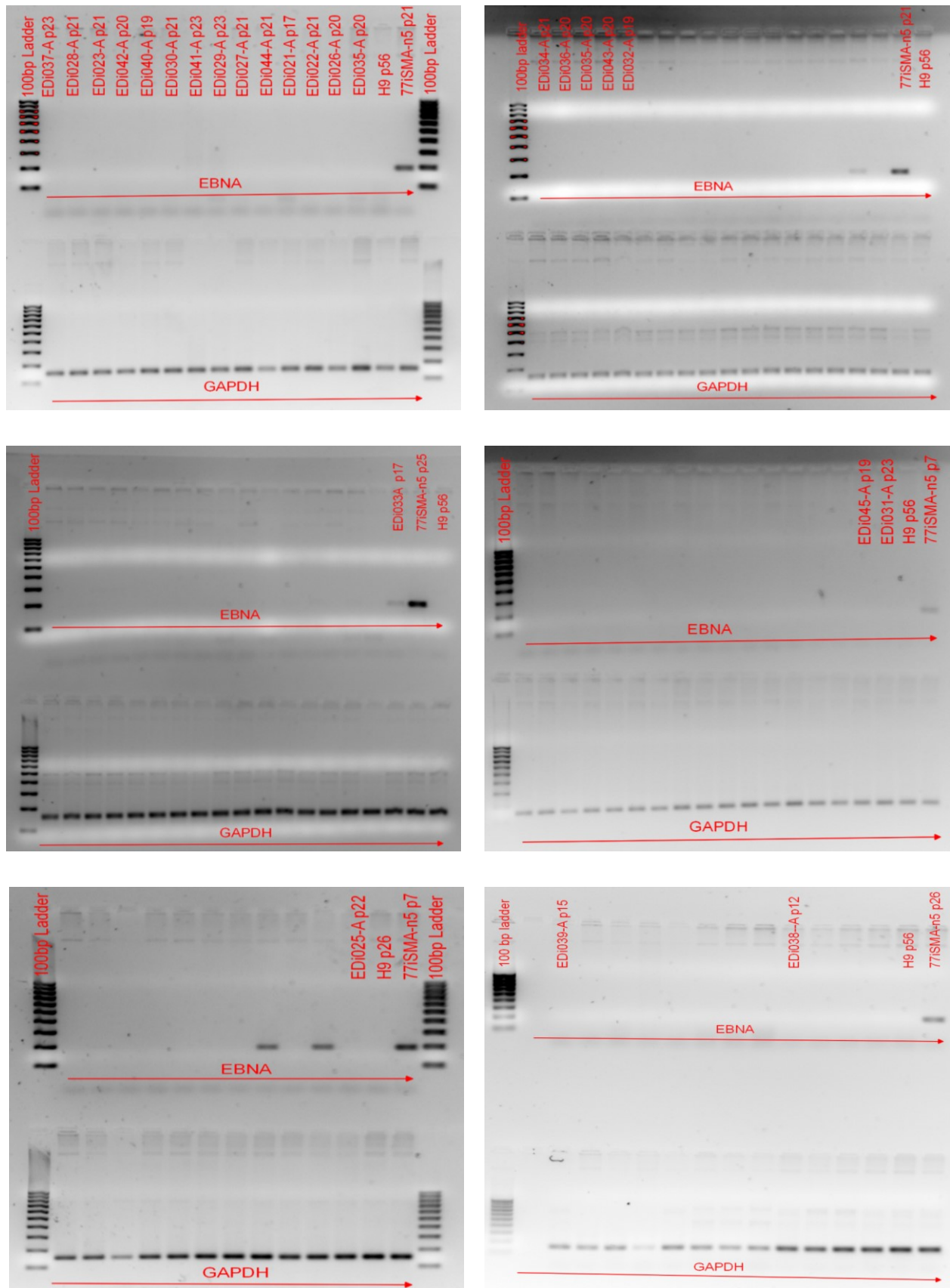

**sFig25:** 2% Agarose Gel images showing lack of EBNA persistence in all LBC iPSC lines. 77iSMA-n5 p21 = positive control iPSC for EBNA persistence. H9 p56 = Human ESC line H9, negative control for EBNA

## sFig.26: T-cell Clonality

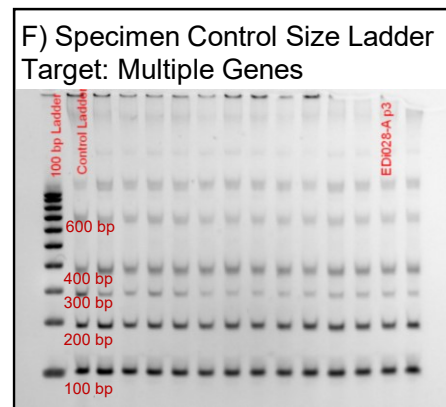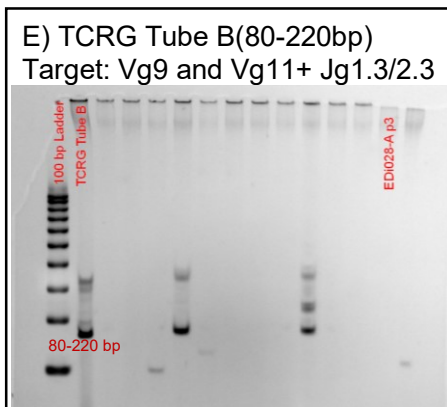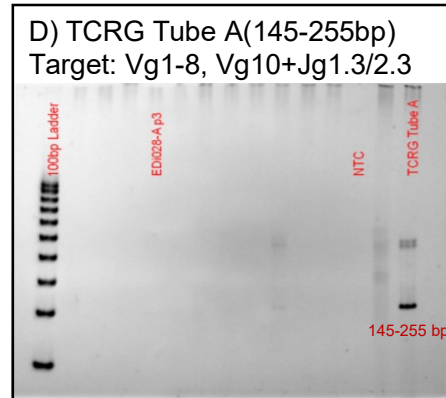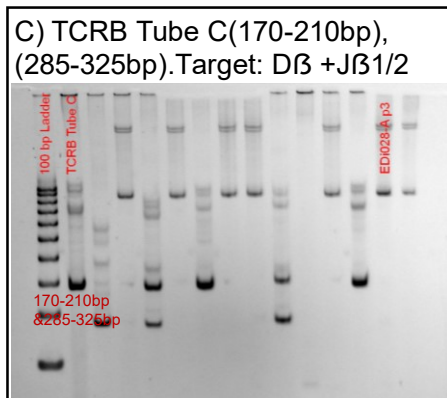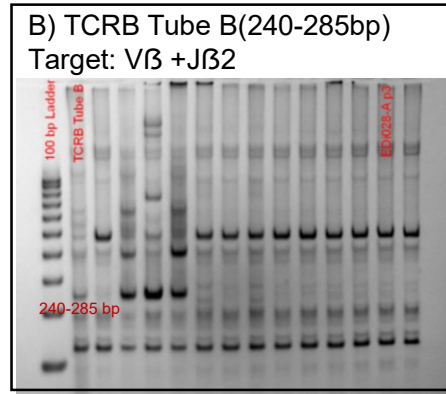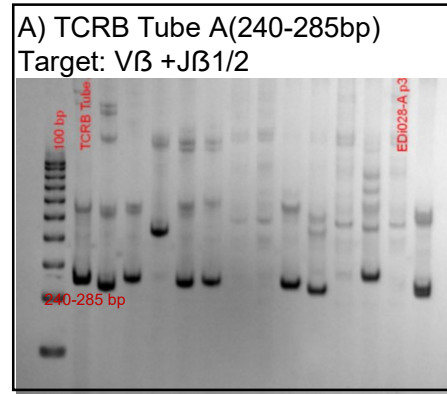


G)

Line	Passage	TCRB	TCRG
EDi021-A	p17	+	+
EDi022-A	p4	-	-
EDi023-A	p3	-	-
EDi025-A	p25	-	-
EDi026-A	p20	-	+
EDi027-A	p5	-	-
EDi028-A	p3	-	-
EDi029-A	p4	-	-
EDi030-A	p6	-	-
EDi031-A	p23	-	-
EDi032-A	p4	-	-
EDi033-A	p17	-	-
EDi034-A	p21	-	-
EDi035-A	p4	+	+
EDi036-A	p6	-	-
EDi037-A	p6	-	-
EDi038-A	p12	-	-
EDi039-A	p7	+	+
EDi040-A	p4	-	-
EDi041-A	p4	-	-
EDi042-A	p6	-	-
EDi043-A	p3	-	-
EDi044-A	p4	-	-
EDi045-A	p19	-	-

**sFig26:** Data showing T-cell clonal lineage for LBC iPSCs. For each line, three targets (V $\beta$  +J $\beta$ 1/2, V $\beta$  +J $\beta$ 2 and D $\beta$  +J $\beta$ 1/2) were tested for T-Cell Receptor Beta Chain (TCRB), and two targets (Vg1-8, Vg10+Jg1.3/2.3; Vg9 and Vg11+ Jg1.3/2.3) for T-Cell Receptor Gamma Chain (TCRG). A-B-C) Representative 6% TBE Gel image for TCRB. D-E) Representative 6% TBE Gel image for TCRG. F) Representative 6% TBE Gel image for Control Genes. G) Table listing each of the LBC iPSC lines, indicating positivity (+) or negativity (-) for each of the T-cell receptors. Lines were considered positive if any band was present in any of the samples.

## sFig.27: Mycoplasma Testing

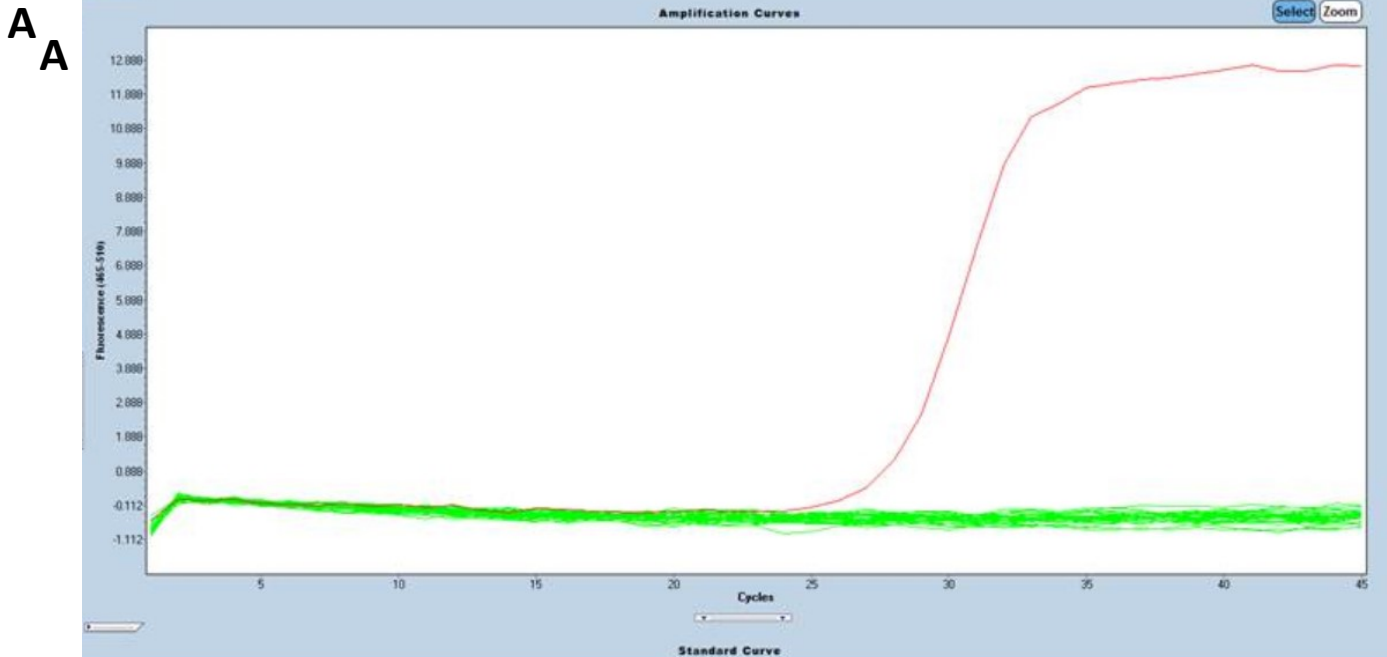

**B**



Sample	Passage	Luminescence Ratio	Mycoplasma
+ve control		48.78 (±25.06)	-
-ve control		0.07 (±0.04)	-
EDi021-A	p24	0.3	-
EDi022-A	p13	0.52	-
EDi023-A	p17	0.5	-
EDi025-A	p23	0.5	-
EDi026-A	p11	0.36	-
EDi027-A	p13	0.45	-
EDi028-A	p24	0.57	-
EDi029-A	p14	0.38	-
EDi030-A	p14	0.4	-
EDi031-A	p29	0.45	-
EDi032-A	p11	0.45	-
EDi033-A			
EDi034-A	p17	0.46	-
EDi035-A	p17	0.47	-
EDi036-A	p13	0.54	-
EDi037-A	p20	0.39	-
EDi038-A			
EDi039-A			
EDi040-A	p17	0.47	-
EDi041-A	p38	0.19	-
EDi042-A	p17	0.43	-
EDi043-A			
EDi044-A	p13	0.41	-

**sFig.27:** Data showing the results of mycoplasma testing from iPSCs stored at two sites. **A)** PCR conducted by IDEXX BioAnalytics for samples provided by Cedars-Sinai Medical Center. All 24 LBC lines shown in green, positive control in red. **B)** Table showing results of MycoAlert Lonza LT07-318 assay conducted by Dementias Platform UK for available lines. Ratio values stated are the ratio of the luminescence signal of the kit substrate to that of the kit reagent. A ratio of 0-0.999 is negative for mycoplasma, 1–1.3 is borderline (requiring retest), and >1.3 is positive for mycoplasma. All samples tested at both sites were negative for mycoplasma.



[Click here to access/download](#)

**STR analysis**  
**LBC STR Analysis.docx**



### Declaration of interests

☐ The authors declare that they have no known competing financial interests or personal relationships that could have appeared to influence the work reported in this paper.

☒ The authors declare the following financial interests/personal relationships which may be considered as potential competing interests:

US patent US 10,221,395 B2 has been granted describing some of the methods to reprogram to iPSCs. Apart from this issued patent filing the authors have declared that no other competing financial interests exist.