



THE UNIVERSITY *of* EDINBURGH

Edinburgh Research Explorer

Derivation of the clinical grade human embryonic stem cell line RCe013-A (RC-9)

Citation for published version:

De Sousa, P, Tye, BJ, Bruce, K, Dand, P, Russell, G, Collins, DM, Greenshields, A, McDonald, K, Bradburn, H, Canham, MA, Kunath, T, Downie, JM, Bateman, M & Courtney, A 2016, 'Derivation of the clinical grade human embryonic stem cell line RCe013-A (RC-9)', *Stem Cell Research*.
<https://doi.org/10.1016/j.scr.2016.04.020>

Digital Object Identifier (DOI):

[10.1016/j.scr.2016.04.020](https://doi.org/10.1016/j.scr.2016.04.020)

Link:

[Link to publication record in Edinburgh Research Explorer](#)

Document Version:

Peer reviewed version

Published In:

Stem Cell Research

Publisher Rights Statement:

Under a Creative Commons license

General rights

Copyright for the publications made accessible via the Edinburgh Research Explorer is retained by the author(s) and / or other copyright owners and it is a condition of accessing these publications that users recognise and abide by the legal requirements associated with these rights.

Take down policy

The University of Edinburgh has made every reasonable effort to ensure that Edinburgh Research Explorer content complies with UK legislation. If you believe that the public display of this file breaches copyright please contact openaccess@ed.ac.uk providing details, and we will remove access to the work immediately and investigate your claim.



Accepted Manuscript

Derivation of the clinical grade human embryonic stem cell line RCe013-A (RC-9)

P.A. De Sousa, B.J. Tye, K. Bruce, P. Dand, G. Russell, D.M. Collins, A. Greenshields, K. McDonald, H. Bradburn, M.A. Canham, T. Kunath, J.M. Downie, M. Bateman, A. Courtney

PII: S1873-5061(16)30032-0
DOI: doi: [10.1016/j.scr.2016.04.020](https://doi.org/10.1016/j.scr.2016.04.020)
Reference: SCR 761

To appear in: *Stem Cell Research*

Received date: 22 April 2016
Accepted date: 28 April 2016



Please cite this article as: De Sousa, P.A., Tye, B.J., Bruce, K., Dand, P., Russell, G., Collins, D.M., Greenshields, A., McDonald, K., Bradburn, H., Canham, M.A., Kunath, T., Downie, J.M., Bateman, M., Courtney, A., Derivation of the clinical grade human embryonic stem cell line RCe013-A (RC-9), *Stem Cell Research* (2016), doi: [10.1016/j.scr.2016.04.020](https://doi.org/10.1016/j.scr.2016.04.020)

This is a PDF file of an unedited manuscript that has been accepted for publication. As a service to our customers we are providing this early version of the manuscript. The manuscript will undergo copyediting, typesetting, and review of the resulting proof before it is published in its final form. Please note that during the production process errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.

Title: Derivation of the clinical grade human embryonic stem cell line RCe013-A (RC-9)

Authors:

P.A. De Sousa^{a,b,c}, B.J. Tye^a, K. Bruce^a, P. Dand^a, G. Russell^a, D.M. Collins^a, A. Greenshields^a, K. McDonald^a, H. Bradburn^a, M.A. Canham^c, T. Kunath^c, J.M. Downie^a, M. Bateman^a, A. Courtney^a

Affiliations:

a. Roslin Cells Limited, Nine Edinburgh Bio-Quarter, 9 Little France Road, Edinburgh, EH16 4UX, UK

b. Centre for Clinical Brain Sciences, University of Edinburgh, UK

c. MRC Centre for Regenerative Medicine, University of Edinburgh, UK

Abstract:

The human embryonic stem cell line RCe013-A (RC-9) was derived under quality assured compliance with UK regulation, European Union Directives and International guidance for tissue procurement, processing and storage according to Good Manufacturing Practice (GMP) standards. The cell line was derived from a failed to fertilise oocyte voluntarily donated as unsuitable and surplus to fertility requirements following informed consent. RCe013-A (RC-9) shows normal pluripotency marker expression and differentiation to the three germ layers in vitro and in vivo. It has a normal 46XY male karyotype and microsatellite PCR identity, HLA and blood group typing data are available.

Resource table:

| | |
|-----------------------------|---|
| Name of stem cell construct | RCe013-A |
| Alternative name | RC-9, RC9 |
| Institution | Roslin Cells Ltd. |
| Person who created resource | B. Tye, K. Bruce, P. Dand, G. Russell, D.M. Collins, A. Greenshields, J. Marshall, K. McDonald, H. Bradburn, M.A. Canham |
| Contact person and email | Paul.desousa@roslincells.com; Paul.desousa@ed.ac.uk Janet.downie@roslincells.com Aidan.courtney@roslincells.com Malcolm.bateman@roslinfoundation.com Tilo.kunath@ed.ac.uk |

| | |
|---|--|
| Date archived/stock date | 10 September 2014 (Seed lot, Passage 28) |
| Type of resource | Biological reagent: cell line |
| Sub-type | hESC, clinical grade |
| Origin | Zygote (Oocyte/1PN) |
| Key transcription factors | Oct4 (confirmed by flow cytometry), |
| Authentication | See Quality Control Certificate of Analysis (Fig. 1) |
| Link to related literature (direct URL links and full references) | N/A |
| Information in public databases | http://hpscreg.eu/cell-line/RCe013-A |
| Ethics | Informed consent obtained. Scotland A Research Ethics committee approval obtained (07/MRE00/56). Conducted under the UK Human Fertilisation and Embryology Authority licence no R0136 to centre 0202 and UK Human Tissue Authority (HTA) licensing number 22631. |

Resource details

RCe013-A (RC-9) was received as a failed to fertilise oocyte/1PN (pro-nuclear) embryo that was surplus to requirement or unsuitable for clinical use due to late development. Human embryonic stem cell (hESC) isolation, expansion and qualification was performed in a facilities whose specification, operation and monitoring complied with GMP standards enabling; i) a fully traceable procurement procedure with informed ethical consent which includes provision for commercial use, ii) detailed medical history and blood borne virus (BBV) screening of donors, and iii) compilation of a cell line history providing details on hESC manufacturing process and quality control testing regime.

Human ESC culture and processing was performed in a grade A tissue culture cabinet in a grade B clean room environment monitored for particulate and microbiological contamination during cell processing in accordance with Rules and Guidance for Pharmaceutical Manufacturers and

Distributors - The Orange Guide, compiled by the UK Medicines Healthcare Products Regulatory Authority (Go to: <https://www.gov.uk/guidance/good-manufacturing-practice-and-good-distribution-practice>). Accordingly, the facility was operating under a mature Quality Management System, compliant with ISO9001:2008 standards. Further hESC derivation was performed under licensure from the UK HFEA (R0136 to centre 0202) and HTA (Licensing Number 22631).

The embryo was grown to blastocyst stage and the cell line was derived by whole embryo outgrowth on mitotically inactivated human dermal fibroblast (HDF) feeder cells. HDFs were derived and manufactured according to GMP and had been approved for clinical use by the Food and Drug Administration, USA. During derivation on HDFs, hESCs were grown in a xeno-free cell therapy grade media (XF KODMEM) supplemented with xeno-free human recombinant bFGF. hESCs were subsequently expanded in a GMP grade serum-free medium (StemPro hESC Serum Free Medium,) and xeno-free matrix (CellStart). The former contained bovine serum albumin (BSA) from a Transmissible Spongiform Encephalopathy (TSE)-free country of origin. The cell line was cryopreserved in a GMP compliant cryopreservation solution (CryoStor CS10).

By flow cytometry, RCE013-A (RC-9) expressed the pluripotency makers Oct-4, Tra-1-60 and SSEA-4 (98.9%, 99.2% and 100.0%, respectively), whereas low expression of the differentiation marker SSEA-1 (3.0%) was observed (Fig. 1, Fig. 2). Differentiation to the three germ layers, endoderm, ectoderm and mesoderm, was demonstrated using embryoid body formation in vitro, and expression of the germ layer markers α -fetoprotein, β -tubulin and muscle actin was observed (Fig. 3, top panel). In vivo teratoma formation yielded typical hESC derived teratomas. Histological examination of fixed and stained sections clearly showed generation of cell types from ectoderm and mesoderm lineages. Endoderm differentiation was also present, but the quality of the structures meant these could not be fully characterised (Fig. 3, bottom panels).

A microsatellite PCR profile has been obtained for the cell line, and HLA Class I and II typing is available (Table 1). Blood group genotyping gave the blood group BB (Table 1). The cell line is free from mycoplasma contamination as determined by RT-qPCR.

Verification and authentication

The cell line was analysed for genome stability by G-banding and showed a normal 46XY male genotype (Fig. 4). SNP genotyping was carried out using the Illumina HumanCytoSNP-12 v2.1 BeadChip and revealed a 120kb region deleted on chromosome 8q24.23 as described in Canham et al (2015). In addition to containing no protein coding genes, this region has been observed many times in healthy individuals, occurring with a frequency of 3.85% as documented in the Database of Genomic Variants (MacDonald et al, 2014).

Materials and methods

Ethics

Derivation of hESC from surplus to requirement and failed to fertilise/develop embryos was approved by The Scotland A Research Ethics Committee and local ethics board at participating fertility clinics and conducted under licence no R0136 from the UK HFEA with informed donor consent. The processing and storage of hESC cells for human application was conducted under licence number 22631 from the UK Human Tissue Authority.

Cell culture

Fresh embryos were cultured EmbryoAssist (Origio (Medicult), Denmark) until Day 3 and BlastAssist (Origio) after Day 3 of development. Embryos were cultured at 36.5 - 37.5°C, 5.0 ± 0.5% CO₂, 5.0 ± 0.5% O₂ in drops under paraffin oil (Origio) and transferred to fresh medium at least every 2-3 days.

By Day 8 of development, embryos were placed in derivation conditions consisting of mitotically inactivated GMP grade neonatal human dermal fibroblasts (HDFs) (Forticell Biosciences, NJ, USA) on tissue culture plastic in XF KODMEM medium (Knockout-DMEM, 15% KOSR-XF, 2mM L-glutamine, 1% MEM Non essential amino acids, 2% XF Growth Factor Cocktail, 0.1mM β-mercaptoethanol (ThermoFisher Scientific, Paisley, UK) supplemented with 80 ng/ml human bFGF (ThermoFisher Scientific). When available, cell therapy system quality reagents were used. If required, assisted hatching was performed by removing the zona pellucidae mechanically using Swemed cutting tools (Vitrolife, Göteborg, Sweden).

HDF cells were cultured in DMEM (Lonza, Slough, UK), 10% Pharma grade FCS (GE Healthcare (PAA), Buckinghamshire, UK) and 2 mM L-glutamine (ThermoFisher Scientific). HDF were mitotically inactivated using gamma irradiation at 50 Gy using a Gammacell Elite 1000 machine. For use as a feeder layer, irradiated HDFs were plated at 50000 cells/cm² in XF KODMEM medium supplemented with 80 ng/ml human bFGF (ThermoFisher Scientific). Cells were cultured at 36.5 - 37.5°C, 5.0 ± 0.5% CO₂, 5.0 ± 0.5% O₂ and 50% medium exchanged 6 days a week.

The established cell line was expanded and banked using CellStart matrix and Stempro hESC Serum Free Medium (cell therapy system quality reagents, ThermoFisher Scientific). This contained BSA from a TSE-free country of origin. Passaging was performed mechanically using an EZ passage tool (ThermoFisher Scientific). hESC lines were expanded to 25-30 wells of a 6-well plate and cryopreserved in 0.5-1 ml Cryostor CS10 (Biolife Solution, Washington, USA) using an EF600-107 controlled rate freezer (Grant Instruments, Cambridge, UK) before being stored in a -150°C freezer (Panasonic Biomedical, Loughborough, UK).

Mycoplasma

In process mycoplasma detection was performed using Applied Biosystems PrepSEQ™ Mycoplasma Nucleic Acid Extraction Kit and MicroSEQ™ Mycoplasma Real-Time PCR Detection Kit (ThermoFisher Scientific (Applied Biosystems)) according to the manufacturer's instruction. European pharmacopoeia (EP) mycoplasma testing was carried out by Moredun Scientific Ltd. (Edinburgh, UK), under a quality and technical agreement.

Endotoxin

Endotoxin levels were determined using the Kinetic-QCL assay (Lonza) and an incubating plate reader (BioTek ELx808) according to the manufacturer's instructions. Briefly, an unknown sample was compared with a standard curve of known levels of control endotoxin. An assay was deemed valid if the coefficient of correlation, $r \geq 0.980$ and the CV (%) for the standard curve was $\leq 10\%$, and the reaction time of the negative control was greater than the reaction time of the lowest standard on the standard curve.

Flow cytometry

Pluripotency was determined using the Human and Mouse Pluripotent Stem Cell Analysis kit (BD, Oxford, UK). Oct 3/4 and SSEA-4 were included as pluripotency markers, and SSEA-1 as a

differentiation marker. FITC conjugated Tra-1-60 (BD) was used as an additional pluripotency marker. Fixed and permeabilised cells were analysed using a FACS Aria flow cytometer (BD) or a Guava easyCyte flow cytometer (Millipore, Watford, UK). Percentage expression of each marker was compared to isotype control or unstained cells.

Viability

Viability was determined using the Guava ViaCount assay. Briefly, the Guava Viacount reagent (Millipore) containing a nuclear and a viability dye, was mixed with a single cell suspension, incubated for 5 minutes and analysed using the Guava easyCyte flow cytometer (Millipore). Total cell count, viable cell count and percentage viable cells was obtained.

In vitro differentiation

hESCs were pre-treated for 1 h with 10 μ M ROCK inhibitor in Stempro hESC SFM (ThermoFisher Scientific) and embryoid bodies (EBs) generated in ultra low attachment plates (Corning) for 7 days before being transferred into EB medium (20% FBS (GE Healthcare (PAA)), 80% KO-DMEM 1 mM L-glutamine, 3.5 μ M β -mercaptoethanol, 1 % nonessential amino acids (all ThermoFisher Scientific)), on glass slide tissue culture chambers (Nunc, ThermoFisher Scientific) coated with 0.2% gelatin (Sigma Aldrich, Dorset, UK) at 0.1 ml/cm² for 14 days.

Immunocytochemistry

Cells were fixed in 4% paraformaldehyde (ThermoFisher Scientific (Alfa Aesar)), permeabilised using 100% ethanol (ThermoFisher Scientific) and stained with AFP (1:500; Sigma), β -tubulin III (1:1000; Sigma) and muscle-specific actin (1:50; DAKO, Glostrup, Denmark) and secondary antibody anti mouse IgG-AlexaFluor 488 (1:200; Sigma). Images were acquired using a Zeiss S100 Axiovert fluorescence microscope or Nikon eC1 confocal microscope

In vivo differentiation

The developmental potential of hESC lines to form teratomas consisting of tissues representative of all three germ layers was evaluated following transplantation under kidney capsule in NOD scid gamma mice. After three months, the animals were culled and assessed for teratoma formation. Teratomas were fixed in 4% paraformaldehyde, embedded in paraffin wax and serial sections of 7 μ m thickness were cut according to standard procedures. For histological assessment, the tissue sections were dewaxed, rehydrated and stained with Masson staining. Tissue sections were analysed using bright field and microscopy and digital images were recorded.

SNP Genotyping and Analysis

DNA samples were assayed using the Illumina HumanCytoSNP-12 v2.1 BeadChip. Genotyping data was initially assessed using GenomeStudio genotyping module (v1.94, Illumina). Karyostudio (v1.4, Illumina) was employed to perform automatic normalisation and to identify genomic aberrations utilising default settings of the built-in cnvPartition algorithm (3.07, Illumina) to generate B-allele frequency and smoothened Log R ratio plots for detected regions. These parameters are designed to detect CNVs greater than 75 kb and CN-LOH regions larger than 1 MB with a confidence value greater than 35. All identified regions were first cross-matched to the Database of Genomic Variants (DGV; <http://dgv.tcag.ca>) to identify naturally-occurring structural variations in the human. CNVs that were not identified on the DGV were then checked against a list of ES cell-associated culture adaptation genomic variants published by the International Stem Cell Initiative (Amps et al, 2011). See also Canham et al, 2015 for further details.

Genomic analysis and outsourced assays

All outsourced assays were carried out under a Quality and Technical Agreement. DNA was extracted using the QIAamp DNA Mini kit (Qiagen, Manchester, UK) according to manufacturer's recommendations and provided in recommended quantities to the service providers.

Microsatellite PCR, or Short Tandem Repeat analysis, was used to determine cell line identity and was carried out by Public Health England. A profile was obtained for the following core alleles: vWA, D16S539, Amelogenin, THO1, CSF1PO, D5S818, D7S820, D13S317 and TPOX.

Human Leukocyte Antigen (HLA) tissue typing was carried out by the Scottish National Blood Transfusion Service.

Blood group genotyping was carried out by the Molecular Diagnostics laboratory at NHSBT.

Karyotype analysis was carried out by The Doctors Laboratory (London, UK) or the Western General Cytogenetics Laboratory (Edinburgh, UK). Live cells at 60-70% confluency were shipped in warm containers, fixed and analysed by standard G-banding analysis. For clinical grade lines, 30 spreads were analysed.

Viral screening for cytomegalovirus (CMV), Human T-cell lymphotropic virus (HTLV), Human immunodeficiency virus (HIV)-1, Hepatitis C virus (HCV), Hepatitis B virus (HBV), Epstein-Barr virus (EBV) and syphilis was carried out by The Doctors Laboratory.

European pharmacopoeia (EP) sterility testing was carried out by Moredun Scientific Ltd. using the culture method.

Figures and tables



Quality Control Test Certificate

Seed Lot Test Results Sample Point 2

| | | | |
|---------------------------|------------------------|---------|---|
| Certificate Number | QCC-15-011 | Version | 1 |
| Sample ID: | RC-9 P28 CS P022SL-001 | | |
| Sample Grade | CLINICAL | | |
| Sample Point: | SP2 (MS/PDP/10) | | |
| Project Number: | P022 | | |
| Date of Cryopreservation: | 10 Sep 14 | | |

| Assay | Test Method | Test Reference | Date of Assay | Specification | Result |
|--------------------------------|----------------------------|----------------------------------|---------------------|---|----------------------------|
| Pluripotency / Differentiation | Flow Cytometry | SOP/QCP/25 v5 and SOP/QCP/31 v 6 | 10 Sep 14 | >65% | SSEA-4 – 99.8% |
| | | | | >65% | Oct 3/4 – 97.4% |
| | | | | >50% | Tra-1-60 – 99.7% |
| | | | | <15% | SSEA-1 – 0.0% |
| Viability | Flow Cytometry | SOP/QCP/69 v 5 | 10 Sep 14 | >60% | 67.7% |
| Embryoid Body Formation | Culture and Immunostaining | SOP/QCP/58 v2 and SOP/QCP/7 v5 | 11 Sep to 07 Oct 14 | Markers for all 3 germ layers - ectoderm, endoderm and mesoderm detected. | All 3 germ layers detected |
| Karyology* | G-banding | TDL Genetics | 11 Sep 14 | Normal Karyotype (46XY or 46XX) | 46XY |

*Outsourced to an approved Third Party.

Document History:

| VERSION | AMENDMENT | CURRENT VERSION |
|---------|---------------|-----------------|
| NA | New Document. | 1 |

Certificate Prepared by (QC): Asd

Date: 09 Mar 16

Certificate Reviewed by (QC): Asd

Date: 09 Mar 16



Quality Control Test Certificate
Seed Lot Test Results Sample Point 3

| | | | |
|----------------------------------|-------------------------------|----------------|----------|
| Certificate Number | QCC-15-019 | Version | 1 |
| Sample ID: | RC-9 P28 CS P022SL-001 | | |
| Sample Grade | CLINICAL* | | |
| Project Number: | P022 | | |
| Sample Point: | SP3 (MS/PDP/10) | | |
| Date of Cryopreservation: | 10 Sep 14 | | |

| Assay | Test Method | Test Reference | Date of Assay | Specification | Result |
|--------------------------------|---|---|------------------------|---|--------------|
| Endotoxin EP | Kinetic Chromogenic LAL Ph.Eur <2.6.14> | RCL SOP/QCP/12 v9 | 16 Oct 14 | <5EU/ml | <0.0627 |
| Mycoplasma EP** | Culture Method Ph.Eur <2.6.7> | SP-GSM.6956 v5 (MSL) | 10 Oct 14 to 07 Nov 14 | Absence of detectable Mycoplasma contamination | Pass |
| Sterility EP** | Direct Inoculation Ph.Eur <2.6.1> USP<71> | SP-GSM.6958v4 (MSL) | 15 Oct 14 to 29 Oct 14 | Absence of detectable bacterial and fungal growth | Pass |
| Pluripotency / Differentiation | Flow Cytometry | RCL SOP/QCP/25 v5 and RCL SOP/QCP/31 v6 | 01 Oct 14 | SSEA-4 >65% | 100.0% |
| | | | | Oct 3/4 >65% | 97.6% |
| | | | | Tra-1-60 >50% | 99.2% |
| | | | | SSEA-1 <15% | 3.0% |
| Viability*** | Flow Cytometry | RCL SOP/QCP/69 v5 | 01 Oct 14 | >60% | 64.97%*** |
| Microsatellite Genotyping** | PCR (Public Health England) | RCL SOP/QCP/6 v3 | 30 Oct 14 to 19 Dec 14 | Confirm STR profile matches historical profile | ID confirmed |

*For clarification the bank is Clinical grade, testing for Pluripotency/differentiation and viability was done within a Research Grade lab.

**Outsourced to an approved Third Party.

***Original test failed and investigated under URI0005 and IR/RC/0363 and deemed to be due to data analysis not he viability of the cells.

Confidential

Page 1 of 2

Fig. 1. Quality Control Certificate of Analysis for RCe013-A (RC-9) P28 seed lot.

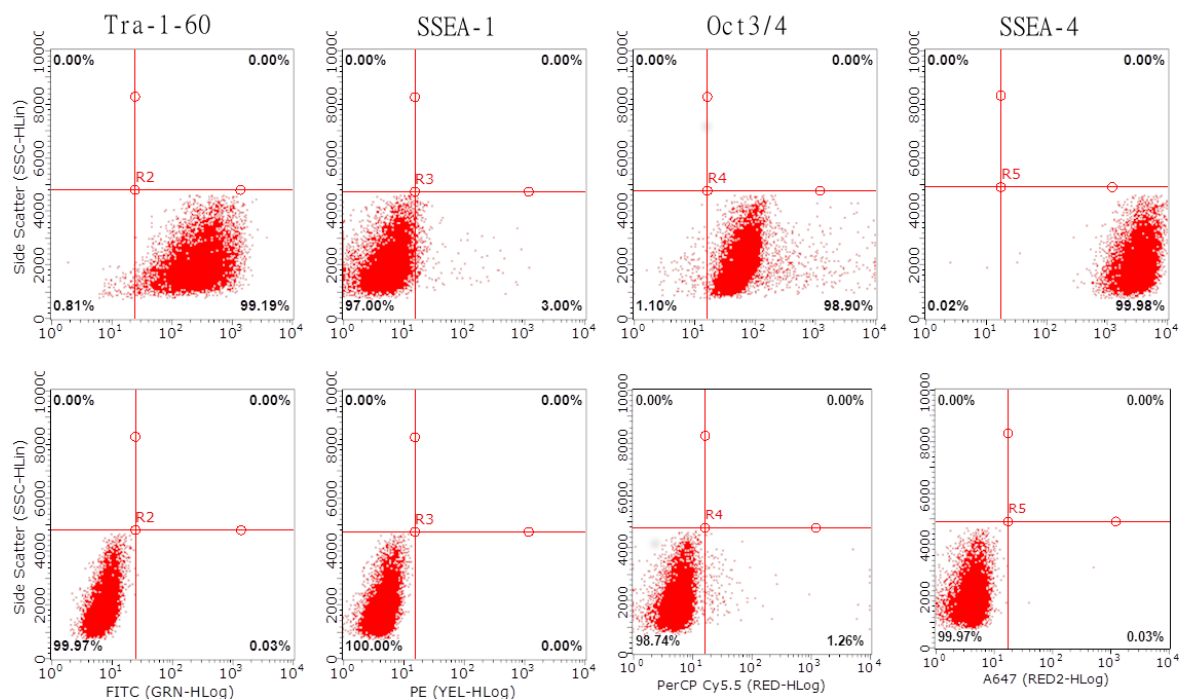
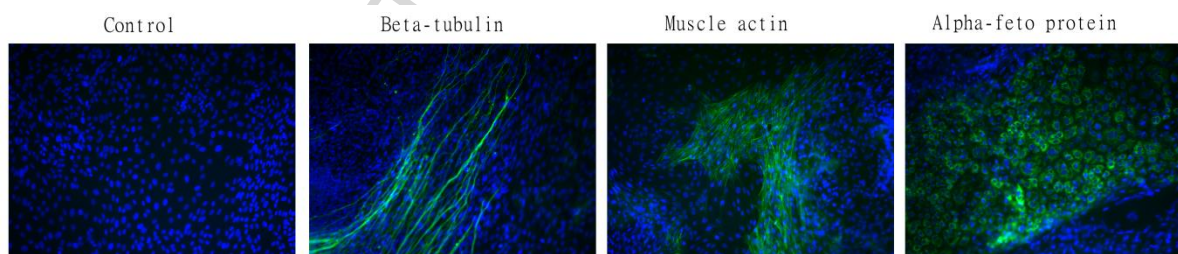


Fig. 2. RCE013-A (RC-9) was subjected to flow cytometry analysis for markers of pluripotency with specific antibody (top row) or isotype control (bottom row) as indicated above the dot plot. Percentage staining is indicated in the respective quadrants. This analysis was carried out at passage 31 during re-banking of the cell line.



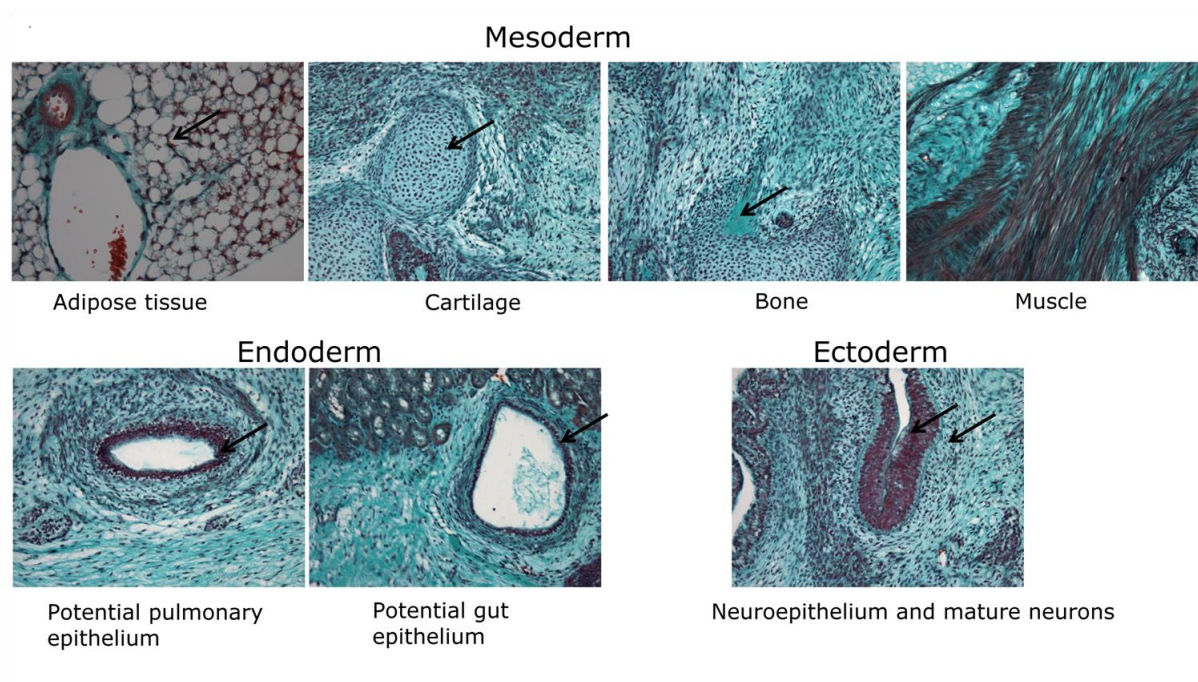


Fig. 3. RCE013-A (RC-9) can differentiate into the three germ layers in vitro and in vivo. Embryoid body mediated differentiation (top panel) resulted in expression of β -tubulin III (ectoderm), muscle actin (mesoderm), and α -fetoprotein (endoderm). Specific staining is shown in green, cell nuclei are counterstained with DAPI (blue). Histological sections of teratomas formed under the kidney capsule of mice yielded clear evidence for mesoderm and ectoderm differentiation and indicated that endoderm structures were also present. Tissue structures identified are indicated by arrows.

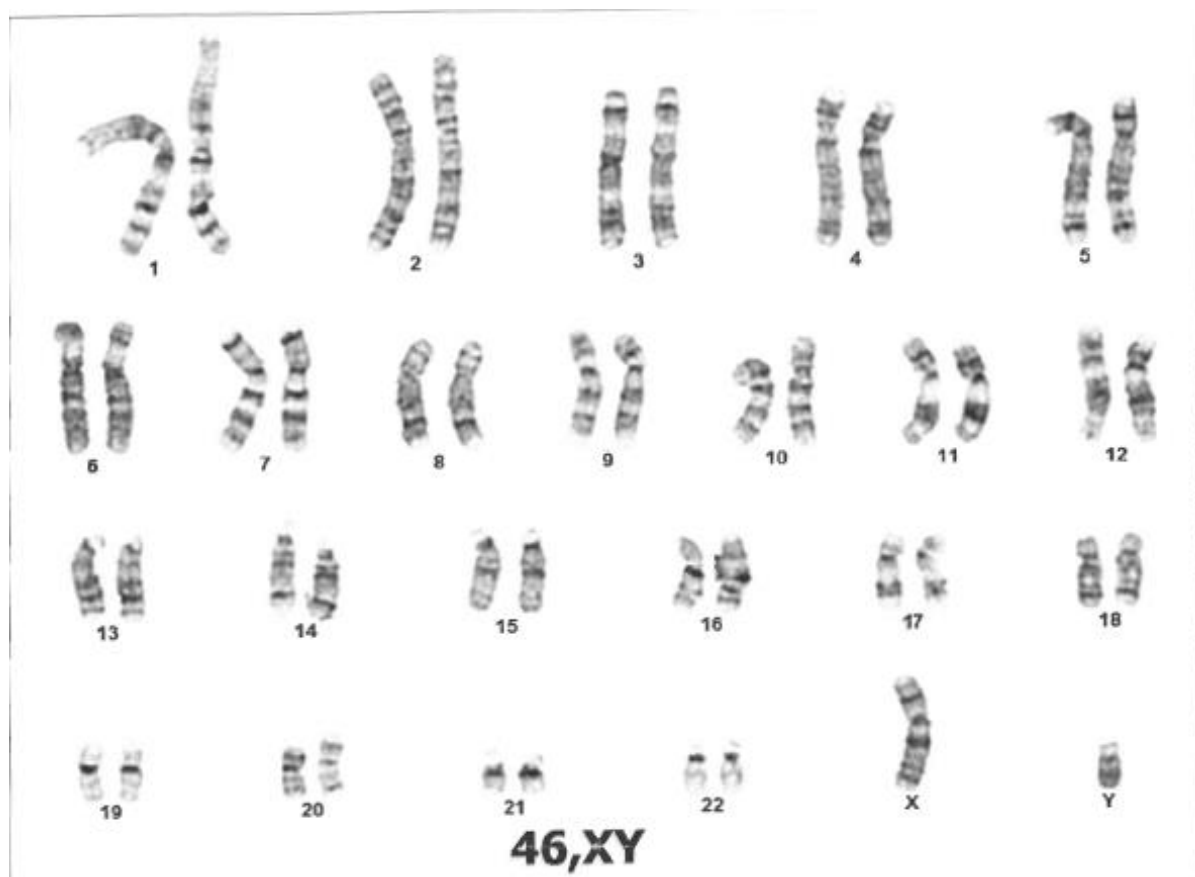


Fig. 4. RCe013-A (RC-9) at passage 28 was analysed by Giesma staining of 30 metaphase spreads and showed a normal 46XY male karyotype at banking.

Table 1. Microsatellite PCR, blood group and HLA tissue typing results for RCe013-A (RC-9).

| Microsatellite PCR results | | | | | | | |
|-----------------------------------|--------------|-----------|-----------|-----------|-----------|-----------|-----------|
| D3S1358 1 | D3S1358 2 | vWA 1 | vWA 2 | D16S539 1 | D16S539 2 | D2S1338 1 | D2S1338 2 |
| 15 | 16 | 15 | 18 | 12 | 12 | 18 | 23 |
| Amelogenin 1 | Amelogenin 2 | D8S1179 1 | D8S1179 2 | D21S11 1 | D21S11 2 | D18S51 1 | D18S51 2 |
| X | Y | 13 | 15 | 30.2 | 32.2 | 13 | 22 |
| D19S433 1 | D19S433 2 | TH01 1 | TH01 2 | FGA 1 | FGA 2 | CSF1PO 1 | CSF1PO 2 |
| 14 | 15 | 6 | 10 | 20 | 20 | 10 | 10 |
| D5S818 1 | D5S818 2 | D7S820 1 | D7S820 2 | D13S317 1 | D13S317 2 | TPOX 1 | TPOX 2 |
| 10 | 13 | 9 | 12 | 12 | 13 | 8 | 8 |
| Blood group genotyping | | | | | | | |
| RhD | RhC | Rhc | RhE | Rhe | Fy a | Fy b | Fy GATA |
| pos | pos | pos | neg | pos | neg | pos | neg |
| Jka | Jkb | K | k | M | N | S | s |
| neg | pos | neg | pos | pos | pos | pos | pos |

| | | | | |
|--------------------------|------|---|------|-----|
| Kp a | Kp b | Do a | Do b | ABO |
| neg | pos | neg | pos | BB |
| HLA tissue typing | | | | |
| HLA Class I Type | | HLA-A*01, A*02; B*07, B*08; Cw*07 | | |
| HLA Class II Type | | HLA-DRB1*01, DRB1*15; DRB5*01; DQB1*05, DQB1*06 | | |

Acknowledgements

Research culminating in the derivation of this line was funded by a grant (PM07321) from Scottish Enterprise Economic Development Agency to PDS, MB, and AC. Research culminating in SNP genotype analysis of the cell line was funded by a Medical Research Council RMRC project grant (MR/K017276/1) and The Cure Parkinson's Trust to TK and MAC.

We thank the transgenic service core facility at the Scottish Centre for Regenerative Medicine (Edinburgh) for the teratoma work.

References

- Canham MA et al. 2015. The Molecular Karyotype of 25 Clinical-Grade Human Embryonic Stem Cell Lines. *Sci Rep*. Nov 26; 5, 17258.
- MacDonald JR et al. 2014. The Database of Genomic Variants: a curated collection of structural variation in the human genome. *Nucleic Acids Res*. 42, D986–92.
- Amps K et al. 2011. Screening ethnically diverse hESC identifies a chromosome 20 minimal amplicon conferring a growth advantage. *Nat Biotech* 29: 1131-44